


Micromass Q-Tof *micro* Mass Spectrometer

Operator's Guide

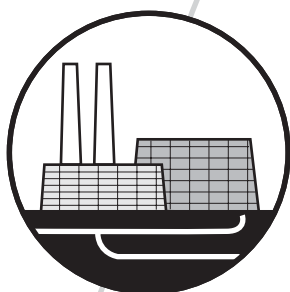

Waters

34 Maple Street
Milford, MA 01757

71500058502 Revision A

NOTICE

The information in this document is subject to change without notice and should not be construed as a commitment by Waters Corporation. Waters Corporation assumes no responsibility for any errors that may appear in this document. This document is believed to be complete and accurate at the time of publication. In no event shall Waters Corporation be liable for incidental or consequential damages in connection with, or arising from, the use of this document.

© 2002 WATERS CORPORATION. ALL RIGHTS RESERVED. THIS DOCUMENT OR PARTS THEREOF MAY NOT BE REPRODUCED IN ANY FORM WITHOUT THE WRITTEN PERMISSION OF THE PUBLISHER.

Micromass and Waters are registered trademarks, and MassLynx is a trademark of Waters Corporation.

All other trademarks or registered trademarks are the sole property of their respective owners.



Safety Information

General

The Micromass Q-ToF *micro* is designed solely for use as a mass spectrometer; any attempt to use it for any other purpose is liable to damage the instrument and will invalidate its warranty.

The Micromass Q-ToF *micro* mass spectrometer conforms to European standard EN61010-1:2001, Safety Requirements for electrical equipment for measurement, control and laboratory use - Part 1: General requirements.

The instrument has been designed and tested in accordance with recognized safety standards. If the instrument is used in a manner not specified by the manufacturer, the protection provided by the instrument may be impaired.

Whenever the safety protection of the instrument has been compromised, disconnect the instrument from all power sources and secure the instrument against unintended operation.

The instrument must be installed in such a manner that the user can easily access and isolate the power source.

Safety Symbols

Warnings in this User's Guide, or on the instrument, must be observed during all phases of service, repair, installation and operation of the instrument. Failure to comply with these precautions violates the safety standards of the design and intended use of the instrument.

Micromass UK Limited assumes no liability for the user's failure to comply with these requirements.

The following safety symbols may be used in the User's Guide, or on the instrument. A **Caution** is an instruction that draws the user's attention to the risk of injury or death; a **Attention** is an instruction that draws attention to the risk of damage to the instrument.



Caution: *This is a general warning symbol, indicating that there is a potential health or safety hazard; the user should refer to this User's Guide for instructions.*



Caution: *This symbol indicates that hazardous voltages may be present*



Caution: *This symbol indicates that hot surfaces may be present.*



Caution: *This symbol indicates that there is danger from corrosive substances.*



Caution: *This symbol indicates that there is danger from toxic substances.*



Caution: *This symbol indicates that there is danger from flammable substances.*



Caution: *This symbol indicates that there is danger from laser radiation.*



Attention: *This is a general caution symbol, indicating that care must be taken to avoid the possibility of damaging the instrument, or affecting its operation.*

Q-ToF micro Mass Spectrometer Information

Intended Use

The Micromass Q-ToF *micro* Mass Spectrometer can be used as a research tool to deliver authenticated exact mass in both MS and MS-MS mode. It is not for use in diagnostic procedures.

Biological Hazard

When you analyze physiological fluids, take all necessary precautions and treat all specimens as potentially infectious. Precautions are outlined in “CDC Guidelines on Specimen Handling,” *CDC – NIH Manual*, 1984.

Calibration

Follow acceptable methods of calibration with pure standards to calibrate methods. Use a minimum of five standards to generate a standard curve. The concentration range should cover the entire range of quality-control samples, typical specimens, and atypical specimens.

Quality Control

Routinely run three quality-control samples. Quality-control samples should represent subnormal, normal, and above-normal levels of a compound. Ensure that quality-control sample results are within an acceptable range, and evaluate precision from day to day and run to run. Data collected when quality-control samples are out of range may not be valid. Do not report this data until you ensure that system performance is acceptable.

Table of Contents

Preface	1
Chapter 1	
Overview	5
1.1 Instrument Description	6
1.1.1 Ionisation Techniques	6
1.1.2 Atmospheric Pressure Chemical Ionisation	6
1.1.3 Electrospray	6
1.1.4 Nanoflow Electrospray	7
1.2 Ion Optics	8
1.3 Internal Layout	9
1.3.1 Mechanical Components	10
1.3.2 Electrical Components	10
1.4 The Vacuum System	12
1.4.1 Fine Pumping	12
1.4.2 Rotary Pumping	13
1.4.3 Pressure Measurement	13
1.4.4 Vacuum Protection	13
1.5 Front Panel Connections	14
1.5.1 Desolvation Gas and Probe Nebuliser Gas	14
1.5.2 High Voltage	14
1.5.3 Heaters	14
1.5.4 Front Panel Controls and Indicators	15
1.6 Rear Panel Connections	16
1.6.1 Water	16
1.6.2 Nitrogen Gas In	16
1.6.3 Exhausts	17
1.6.4 Supply Inlet	17

1.6.5	Electronics	17
1.6.6	Rotary Pump	17
1.6.7	Event Out	17
1.6.8	Contact Closure In	18
1.6.9	Analog Channels	18
1.7	MassLynx Data System	18
1.7.1	Software	19

Chapter 2

Routine Procedures	21
2.1 Start Up Following a Complete Shutdown.....	21
2.1.1 Preparation	21
2.1.2 Pumping	22
2.1.3 MCP Detector Conditioning	23
2.1.4 Instrument Warm-up	24
2.1.5 Using the Instrument	24
2.2 Start Up Following Overnight Shutdown.....	24
2.2.1 Preparation for Electrospray Operation	25
2.2.2 Preparation for APcI Operation	26
2.2.3 Transient Pressure Trip	27
2.2.4 Power Failure	27
2.2.5 Nitrogen Supply	28
2.3 Automatic Startup and Shutdown	28
2.3.1 The Shutdown Editor	28
2.3.2 The Shutdown Editor Toolbar	29
2.3.3 Shutdown Page	29
2.3.4 The Auto Control Tasks Page	33
2.3.1 Shutdown Log	36

Chapter 3

Electrospray	39
3.1 Introduction	39
3.1.1 Post-column Splitting	41
3.1.2 Megaflow	43
3.2 Operation	43
3.2.3 Megaflow Hints	49
3.2.4 Removing the Probe	50
3.3 Sample Analysis and Calibration	50
3.3.5 General Information	50
3.4 Chromatographic Interfacing	52
3.4.1 LC-MS Sensitivity Enhancement	53

Chapter 4

Nanoflow Electrospray	55
4.1 Installing the Interface	57
4.1.1 Operation of the Camera System	60
4.1.2 Glass Capillary Option	61
4.2 Nano-LC Option	63
4.2.1 Installation	63
4.2.2 Operation	65
4.3 Changing Options	66

Chapter 5

APCI	67
5.1 Preparation	68
5.1.1 Checking the Probe	68
5.2 Obtaining an Ion Beam	69
5.3 Hints for Sample Analysis	71
5.4 Removing the Probe	72

Chapter 6

Tuning	75
6.1 Tuning	75
6.2 Tune Page	77
6.2.1 Tune Page Basics	77
6.2.2 Saving and Restoring Parameter Settings	78
6.2.3 Changing the Peak Display	79
6.2.4 Source Tuning Page	81
6.2.5 Quadrupole Tuning Page (Analyzer)	83
6.2.6 Time Of Flight Page	85
6.2.7 Other Tune Page Settings	88
6.3 Obtaining an Ion Beam	89
6.3.1 Standard Tune Parameters	90
6.3.2 Checking the Resolution	90

Chapter 7

Calibration	93
7.1 Introduction	93
7.1.1 Nominal Mass Accuracy - Lteff	94
7.1.2 Calibration	95
7.1.3 Lock Mass	96
7.1.4 Deadtime Correction	96
7.2 Generation of an Instrument Calibration	97
7.2.1 Calculation of Lteff	97
7.2.2 Acquire Data for Calibration	99
7.2.3 Check for Deadtime Distortion	100
7.2.4 Calibration	102
7.2.5 LockMass Correction	106
7.2.6 Deadtime Correction	106
7.3 Exact Mass Measurement: Additional Hints	109

Chapter 8

Data Acquisition	111
8.1 Starting an Acquisition	111
8.1.1 Starting an Acquisition from the Tune Page	111
8.1.2 Multiple Samples	114
8.2 Monitoring an Acquisition	114
8.2.1 The Acquisition Status Window	114
8.2.2 Chromatogram Real-Time Update	115
8.2.3 Spectrum Real-Time Update	115
8.2.4 System Manager	115
8.2.5 Stopping an Acquisition	115
8.3 The Experiment Setup Editor	115
8.3.1 The Experiment Setup Editor Toolbar	116
8.3.2 Adding a New Function	117
8.3.3 Modifying an Existing Function	117
8.3.4 Copying an Existing Function	118
8.3.5 Removing a Function	118
8.3.6 Changing the Order of Functions	118
8.3.7 Setting a Solvent Delay	119
8.3.8 Analog Channels	120
8.3.9 Saving and Restoring a Function List	120
8.3.10 Setting up an MS Scan Function	121
8.4 Automated Data Dependent Acquisition (DDA)	125
8.4.1 Creating a Data Dependent MS to MS/MS Switching Experiment.	126
8.4.2 Survey	127
8.4.3 MS to MSMS	128
8.4.4 MSMS to MS	131
8.4.5 MSMS Template	132
8.4.6 Include Masses	133
8.4.7 Collision Energy	135

8.4.8	Exclude Masses	137
-------	----------------------	-----

Chapter 9

Parent Ion Discovery	139
-----------------------------------	------------

9.1	Creating a Parent Ion Discovery Experiment.....	139
9.1.1	Parent Survey Window	140
9.1.2	Product Ions	142
9.1.3	Neutral Loss	143
9.1.4	MSMS Template	144
9.1.5	MS to MSMS	145
9.1.6	MSMS to MS	148
9.1.7	Include	149
9.1.8	Exclude	151
9.1.9	Collision Energy	152

Chapter 10

Maintenance and Fault Finding	155
--	------------

10.1	Electronics Maintenance	156
10.1.1	Cooling Fans and Air Filters	156
10.1.2	The Vacuum System	156
10.2	The Source.....	159
10.2.1	Cleaning the Sample Cone in Situ	160
10.2.2	Removing and Cleaning the Sample Cone	163
10.2.3	Removing and Cleaning the Source Block and Extraction Cone ...	165
10.2.4	Removing and Cleaning the RF Lens Assembly	167
10.2.5	Reassembling and Checking the Source	168
10.2.6	The Corona Discharge Pin	169
10.3	The Electrospray Probe	170
10.3.1	Replacement of the Stainless Steel Sample Capillary	171
10.4	The APcI Probe	172
10.4.1	Cleaning the Probe Tip	173

10.4.2	Replacing the Probe Tip Heater	174
10.4.3	Replacing the Fused Silica Capillary	175
10.5	The Analyzer	176
10.5.1	Removing the MS1 and Collision Cell Optical Bench Assembly ...	177
10.5.2	Dismantling and Cleaning the Entrance Prefilter	178
10.5.3	Cleaning the MS1 Analyzer Assembly	178
10.5.4	Replacing the MS1 and Gas Cell Optical Bench Assembly	178
10.5.5	The MS2 Analyzer and Detector Assembly	179
10.6	Fault Finding	179
10.6.1	No Beam	179
10.6.2	Unsteady Beam	180
10.6.3	High Back Pressure	181
10.6.4	Loss of Sensitivity	182
10.6.5	Incorrect Isotope Distributions	182
10.6.6	High Noise Levels	182
10.6.7	Poor Analyzer Vacuum	183
10.7	Cleaning Materials	184
10.8	Preventive Maintenance Check List.....	184
10.8.1	Daily.....	185
10.8.2	Weekly	185
10.8.3	Monthly	185
10.8.4	Four-Monthly	185

Appendix A

Reference Information	187
A.1 Positive Ion	188
A.2 Negative Ion	190
A.3 Preparation of Calibration Solutions	192

Appendix B

Instrument Specifications	195
B.1 Room Location	195
B.2 Environment	195
B.3 Water Supplies	196
B.4 Power Requirements	197
B.5 Gases and Regulators	198
B.6 Exhaust Outlets	199
B.7 Performance Specifications	199

List of Figures

1-1	The Micromass Q-tof micro	5
1-2	Q-Tof micro Ion Optics.....	8
1-3	Mechanical Components - Internal View.....	10
1-4	Rear View of Instrument.....	11
1-5	Q-Tof micro Vacuum System.....	12
1-6	Front Panel	14
1-7	Rear Panel Connections	16
1-8	Rear Panel Outputs	18
2-1	MCP Conditioning Dialog	23
2-2	Shutdown Page	30
2-3	Shutdown on Error Configuration Dialog	31
2-4	Configure Event In Dialog.....	32
2-5	Gas Threshold Dialog	33
2-6	Auto Control Tasks Page	33
2-7	Shutdown Startup Log	36
2-1	Log Parameters Dialog	37
3-1	Electrospray Scheme.....	39
3-2	Post Column Split	42
3-3	Comparison of Normal and Megaflow Electrospray	42
3-4	ESI Source	43
3-5	Electrospray Tip	45
3-6	Optimum Probe Position.....	46
3-7	Probe Adjustment	46
4-1	The Nanoflow Interface	55
4-2	Nanoflow Electrospray Scheme.....	56
4-3	Electrospray Source Being Removed	57

4-4	Nanoflow Source Without Cover.....	58
4-5	Nanoflow Source with Cover and End Flange	59
4-6	Nanoflow Camera and Microscope System.....	60
4-7	Nanovial Tip Position	62
4-8	Nano-LC Assembly	64
5-1	Schematic Representation of APcI	67
5-2	Illustration of Probe Adjustment.....	70
5-3	APcI Probe showing twin Thumb Nuts used for Removal.....	72
6-1	The Ion Optical System of the Q-Tof micro	76
6-2	The QToF micro Tune Page.....	77
6-3	The Customize Display Menu	79
6-4	The Customise Plot Appearance Dialog.....	79
6-5	The Source Tuning Page	81
6-6	The Quadrupole Tuning Page	83
6-7	The Time Of Flight Page	85
6-8	The TDC Settings Dialog	88
6-9	The Tune Page Acquisition Dialog.....	91
6-10	Combined Spectrum of Leucine Enkephalin Showing Peak Width at Half Height	92
7-1	TDC Settings Dialog.....	95
7-2	The Tune Page Acquisition Dialog.....	97
7-3	Indicated Mass of Leucine Enkephalin prior to Lteff Correction.....	98
7-4	Indicated Mass of Leucine Enkephalin after Lteff Correction.	99
7-5	The TOF Accurate Mass dialog.....	100
7-6	The TOF Spectrum Center Dialog.....	101
7-7	Comparison of Centered Spectra With (top) and Without (bottom) Deadtime Correction	101
7-8	The Calibration Window.....	102
7-9	The Display Calibration Graphs Dialog	103
7-10	Calibration Window Showing Residual Errors.....	104

7-11	The Calibration Parameters Dialog.....	105
7-12	The Accept Calibration Dialog	105
7-13	Calibration Window Deadtime Distortion	106
7-14	Same data as Figure 7-13 but with Np multiplier of 0.7.....	107
7-15	Two Spectra One (bottom) with no Deadtime Correction and the Other (top) with.	108
8-1	Tune Page Acquisition Dialog	112
8-2	The Scan Report Dialog.....	114
8-3	The Experiment Setup Editor	116
8-4	The Solvent Delay Dialog.....	119
8-5	The Analog Data Dialog	120
8-6	The TOFMS Scan Dialog	122
8-7	The MSMS Scan Dialog.....	124
8-8	Survey Scan Dialog: Qtof MicroSurvey Page	127
8-9	Survey Scan Dialog: MS to MSMS Page	128
8-10	Survey Scan Dialog: MSMS to MS Page	131
8-11	Survey Scan Dialog: MSMS Template Page	132
8-12	Survey Scan Dialog: Include Masses Page.....	133
8-13	Survey Scan Dialog: Collision Energy Page	135
8-14	The CE Profile Dialog	136
8-15	The Modify Charge State Dialog.....	137
8-16	Survey Scan Dialog: Exclude Masses Page.....	138
9-1	Parent Scan Dialog: Parent Survey Page	140
9-2	Example of a Screen Capture.....	142
9-3	Parent Scan Dialog: Neutral Loss Page	143
9-4	Parent Scan Dialog: MSMS Template	144
9-5	Parent Scan Dialog: MSMS to MS Page	148
9-6	Parent Scan Dialog: Include Page.....	149
9-7	Parent Scan Dialog: Exclude Dialog	151
9-8	Parent Scan Dialog: Collision Energy Page.....	152

9-9	The CE Profile Dialog	153
9-10	The Modify Charge State Dialog.....	154
10-1	Rotary Pump	158
10-2	Removing The Source Enclosure	161
10-3	The Source, With Housing Removed	162
10-4	Removing the Cone	164
10-5	Source Block and Extraction Cone Removed.....	165
10-6	Removing The RF Lens Assembly	167
10-7	The RF Lens Assembly.....	168
10-8	Removing the Stainless Steel Capillary form the Source	171
10-9	The APcI Probe Tip Heater.....	174
10-10	Replacing the Fused Silica Capillary on the APcI probe.....	175
10-11	The Prefilter and Main Analyzer of the Q-Tof micro	176

Preface

The *Micromass QToF micro Mass Spectrometer* is intended for a wide variety of users whose familiarity with Mass Spectrometers, computers and software ranges from novice to expert. This guide describes the basics of how to Startup the instrument, obtain an Ion Beam, Calibrate and Acquire Data using MassLynx Software, and maintain the instrument.

Organization

This guide contains the following:

Chapter 1 gives a brief description and overview of the instrument

Chapter 2 describes the routine procedures that are required to Startup and Shutdown the Instrument

Chapter 3 describes the the Electrospray Interface.

Chapter 4 describes the Nanoflow Interface.

Chapter 5 describes the APcI Interface.

Chapter 6 describes the Tuning page and how to Obtain an Ion Beam.

Chapter 7 describes how to calibrate the instrument.

Chapter 8 describes how to Acquire data from the MassLynx sample list and how to set up a method.

Chapter 9 describes a Parent Ion Discovery experiment.

Chapter 10 describes Maintenance and Fault Finding for the instrument.

Related Documentation

Waters Licenses, Warranties, and Support: Provides software license and warranty information, describes training and extended support, and tells how Waters handles shipments, damages, claims, and returns.

Online Documentation

MassLynx Help: Describes all MassLynx windows, menus, menu selections, and dialog boxes for the base software and software options. Also included are help Files

on Inlet Control, Interfacing, Security and any application software that may have been purchased.

MassLynx ReadMe File: Describes product features and enhancements, helpful tips, installation and/or configuration considerations, and changes since the previous version.

Printed Documentation for Base Product

Instrument User's Guides: Provides an introduction to the running and maintenance of the Instrument and also basic instructions on how to acquire data and calibrate the instrument.

MassLynx User's Guide: Provides a comprehensive introduction to the MassLynx software. Describes the basics of how to use MassLynx software to acquire data develop an acquisition method, review and process results, and print a report.

MassLynx Interfacing Guide: Provides information on how to interface MassLynx with other Software applications.

MassLynx Inlet Control Guide: Provides information on how to install and run Autosamplers, LC and GC systems, UV detectors using MassLynx.

MassLynx Security User's Guide: Describes how to add security to your MassLynx system.

Printed Documentation for Software Options

QuanLynx User's Guide: Describes the procedures for installing, configuring and using QuanLynx Software.

OpenLynx User's Guide: Describes the procedures for installing, configuring and using OpenLynx Software.

FractionLynx User's Guide: Describes the procedures for installing, configuring and using FractionLynx Software.

MetaboLynx User's Guide: Describes the procedures for installing, configuring and using MetaboLynx Software.

BioLynx and ProteinLynx User's Guide: Describes the procedures for installing, configuring and using BioLynx and ProteinLynx Software.

MicrobeLynx User's Guide: Describes the procedures for installing, configuring and using MicrobeLynx Software.

NeoLynx User's Guide: Describes the procedures for installing, configuring and using NeoLynx Software.

TargetLynx User's Guide: Describes the procedures for installing, configuring and using TargetLynx Software.

GCLynx User's Guide: Describes the procedures for installing, configuring and using GCLynx Software.

Documentation on the Web

Related product information and documentation can be found on the World Wide Web. Our address is <http://www.waters.com/micromass>.

Documentation Conventions

The following conventions can be used in this guide:

Convention	Usage
Bold	Bold indicates user action such as keys to press, menu selections, and commands. For example, "Click Next to go to the next page."
<i>Italic</i>	Italic indicates information that you supply such as variables. It also indicates emphasis and document titles. For example, "Replace <i>file_name</i> with the actual name of your file."
Courier	Courier indicates examples of source code and system output. For example, "The SVRMGR> prompt appears."
Courier Bold	Courier bold indicates characters that you type or keys you press in examples of source code. For example, "At the LSNRCTL> prompt, enter set password oracle to access Oracle."
Keys	The word <i>key</i> refers to a computer key on the keypad or keyboard. <i>Screen keys</i> refer to the keys on the instrument located immediately below the screen. For example, "The A/B screen key on the 2414 Detector displays the selected channel."
...	Three periods indicate that more of the same type of item can optionally follow. For example, "You can store <i>filename1</i> , <i>filename2</i> , ... in each folder."
>	A right arrow between menu options indicates you should choose each option in sequence. For example, "Select File > Exit " means you should select File from the menu bar, then select Exit from the File menu.

Notes

Notes call out information that is helpful to the operator. For example:

Note: *Record your result before you proceed to the next step.*

Attentions

Attentions provide information about preventing damage to the system or equipment. For example:



Attention: *To avoid damaging the detector flow cell, do not touch the flow cell window.*

Cautions

Cautions provide information essential to the safety of the operator. For example:



Caution: *To avoid burns, turn off the lamp at least 30 minutes before removing it for replacement or adjustment.*



Caution: *To avoid electrical shock and injury, unplug the power cord before performing maintenance procedures.*



Caution: *To avoid chemical or electrical hazards, observe safe laboratory practices when operating the system.*



Caution: *Operating the source without the source enclosure will result in solvent vapor escape.*



Caution: *Strong acid causes burns. Carry out this procedure in a fume cupboard using protective equipment.*



Caution: *Cleaning the source requires the use of solvents and chemicals which may be flammable.*

Chapter 1

Overview

This chapter gives a brief description of the following

- The Micromass Q-Tof *micro* Section 1.1
- The Ion Optics Section 1.2
- The Internal Layout Section 1.3
- The Vacuum System Section 1.4
- Front and Rear Panel Connections Section 1.5 and Section 1.6
- The MassLynx Data System Section 1.7



Figure 1-1 The Micromass Q-tof *micro*

1.1 Instrument Description

The Q-ToF *micro* hybrid quadrupole time of flight mass spectrometer is available with electrospray ionisation (ESI) and atmospheric pressure chemical ionisation (APCI).

Q-ToF *micro* utilises a high performance, research grade quadrupole mass analyser, incorporating a prefilter assembly to protect the main analyser from contaminating deposits, and an orthogonal acceleration time of flight (TOF) mass spectrometer. A hexapole collision cell, between the two mass analysers, can be used to induce fragmentation to assist in structural investigations.

Ions emerging from the second mass analyser are detected by the microchannel plate detector and ion counting system.

A PC computer runs the MassLynx NT software system to control Q-ToF *micro*, and to acquire and process data.

1.1.1 Ionisation Techniques

Using the Micromass Z-spray atmospheric pressure ionisation (API) source, two techniques are available.

1.1.2 Atmospheric Pressure Chemical Ionisation

Atmospheric pressure chemical ionisation (APCI) generally produces protonated or deprotonated molecular ions from the sample via a proton transfer (positive ions) or proton abstraction (negative ions) mechanism. The sample is vapourised in a heated nebuliser before emerging into a plasma consisting of solvent ions formed within the atmospheric source by a corona discharge. Proton transfer or abstraction then takes place between the solvent ions and the sample. Eluent flows up to 2 millilitres/minute can be accommodated without splitting the flow.

1.1.3 Electrospray

Electrospray ionisation (ESI) takes place as a result of imparting a strong electrical charge to the eluent as it emerges from the nebuliser. An aerosol of charged droplets emerges from the nebuliser. These undergo a reduction in size by solvent evaporation until they have attained a sufficient charge density to allow sample ions to be ejected from the surface of the droplet (“ion evaporation”).

A characteristic of ESI spectra is that ions may be singly or multiply charged. Since the mass spectrometer filters ions according to their mass-to-charge ratio (m/z), compounds of high molecular weight can be determined if multiply charged ions are formed.

Eluent flows up to 1 ml/min can be accommodated although it is often preferable with electrospray ionisation to split the flow such that 5-50 μ l/min of eluent enters the mass spectrometer.

1.1.4 Nanoflow Electrospray

The optional nanoflow interface allows electrospray ionisation to be performed in the flow rate range 5 to 1000 nanolitres per minute.

For a given sample concentration, the ion currents observed in nanoflow are comparable to those seen in normal flow rate electrospray. Great sensitivity gains are therefore observed when similar scan parameters are used, due to the great reductions in sample consumption.

1.2 Ion Optics

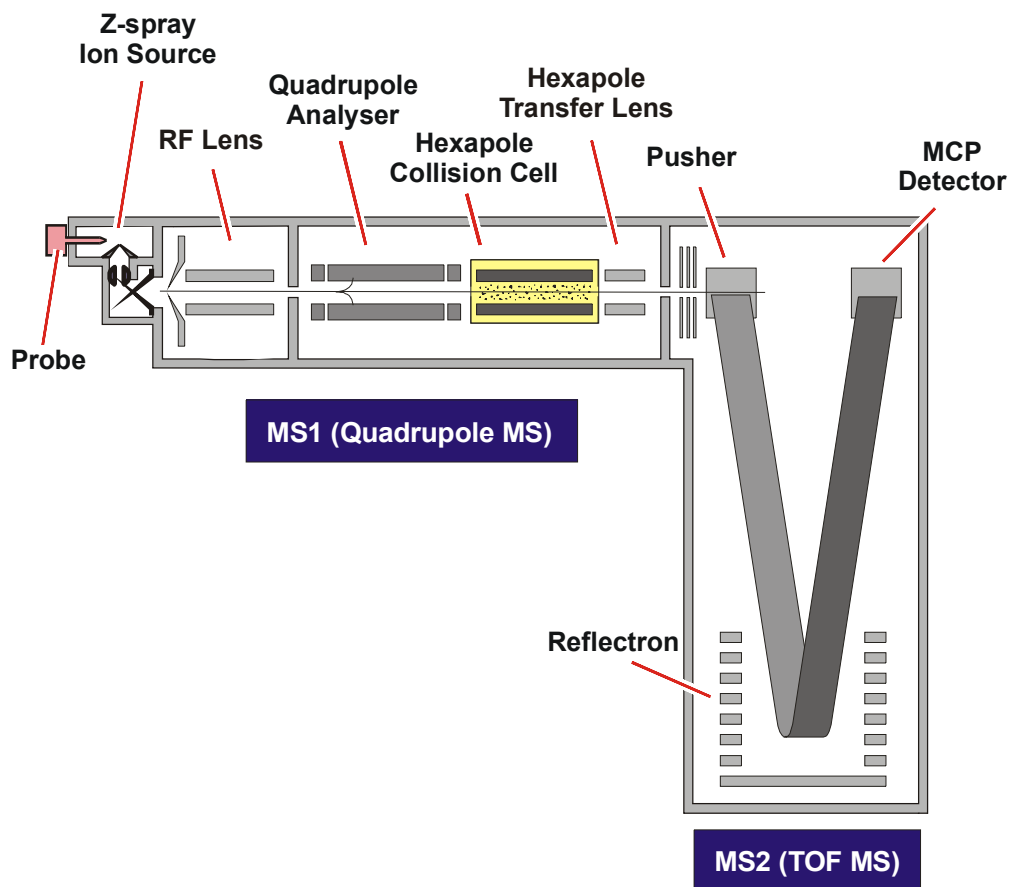


Figure 1-2 Q-ToF *micro* Ion Optics

The principal components of the ion optical system are shown in Figure 1-2.

Ions generated in the Z-spray source are transferred to the quadrupole analyser MS1 via the independently pumped RF lens. After leaving the quadrupole analyser the ions flow into the orthogonal time of flight analyser MS2. The ion beam is focused into the pusher by the acceleration, focus, steer and tube lenses. The pusher then pulses a section of the beam towards the reflectron, which then reflects ions back to the detector.

As ions travel from the pusher to the detector they are separated in mass according to their flight times, with ions of the highest mass to charge ratio (m/z) arriving later.

The pusher may be operated at repetition frequencies of up to 30 kHz, resulting in a full spectrum being recorded by the detector every 33 microseconds. Each spectrum is summed in the histogram memory of the time to digital converter until the histogrammed spectrum is transferred to the host PC.

If the user has requested an acquisition rate of 1 spectrum/second, each spectrum viewed on the host PC will be the result of summing up to 30,000 individual spectra recorded at the detector.

Unlike scanning instruments, the TOF performs parallel detection of all masses within the spectrum at very high sensitivity and acquisition rates. This characteristic is of particular advantage when the instrument is coupled to fast chromatography, since each spectrum is representative of the sample composition at that point in time, irrespective of how rapidly the sample composition is changing.

1.3 Internal Layout



Caution: *The covers should not be removed without first isolating the instrument at the electricity supply.*



Attention: *The internal layout is shown in the following diagrams for information only, and does not imply that labelled components are user-serviceable.*

1.3.1 Mechanical Components

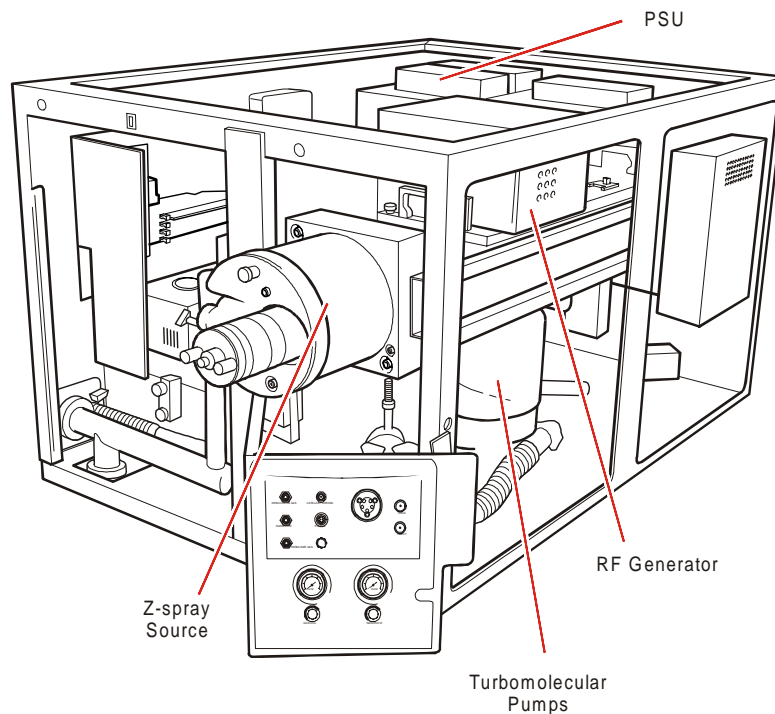


Figure 1-3 Mechanical Components - Internal View.

The main internal mechanical components of the instrument are:

- The source housing, containing the RF (hexapole) lens.
- The MS1 analyser housing, containing the quadrupole analyser, hexapole collision cell and hexapole transfer lens
- The TOF analyser housing, containing the pusher, detector and reflectron assemblies.
- One 250 litre/second turbomolecular pump, plus one split-flow turbomolecular pump.
- Two active inverted magnetron (Penning) gauges and two Pirani gauges.

1.3.2 Electrical Components

The main electronics unit is located in the lower rear section of the instrument. This contains:

- High voltage power supplies.
- These supply the probe or corona, reflectron, TOF flight tube and lens circuits.
- Low voltage power supplies.
- These supply the PCBs, high voltage supplies and turbomolecular pumps.
- Main PCBs.
- For communications, lenses and quadrupole control.

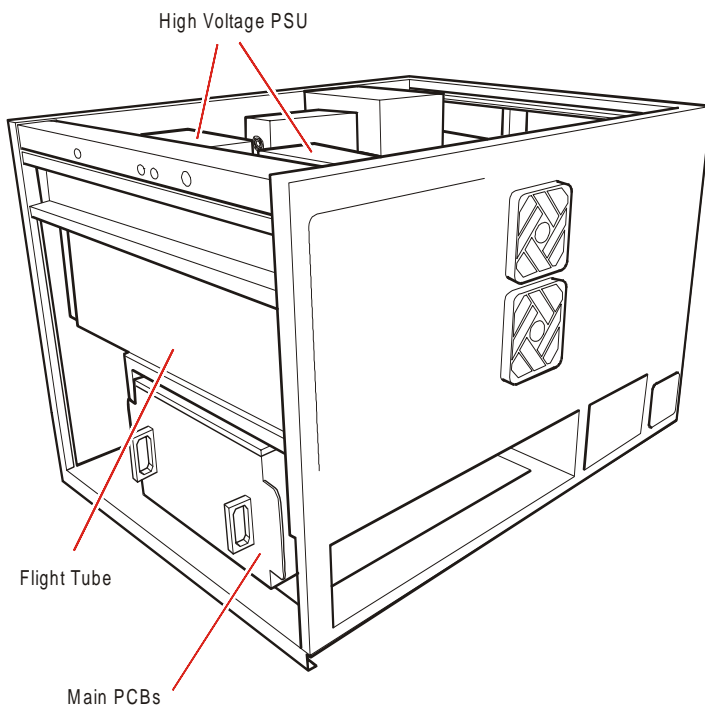


Figure 1-4 Rear View of Instrument

1.4 The Vacuum System

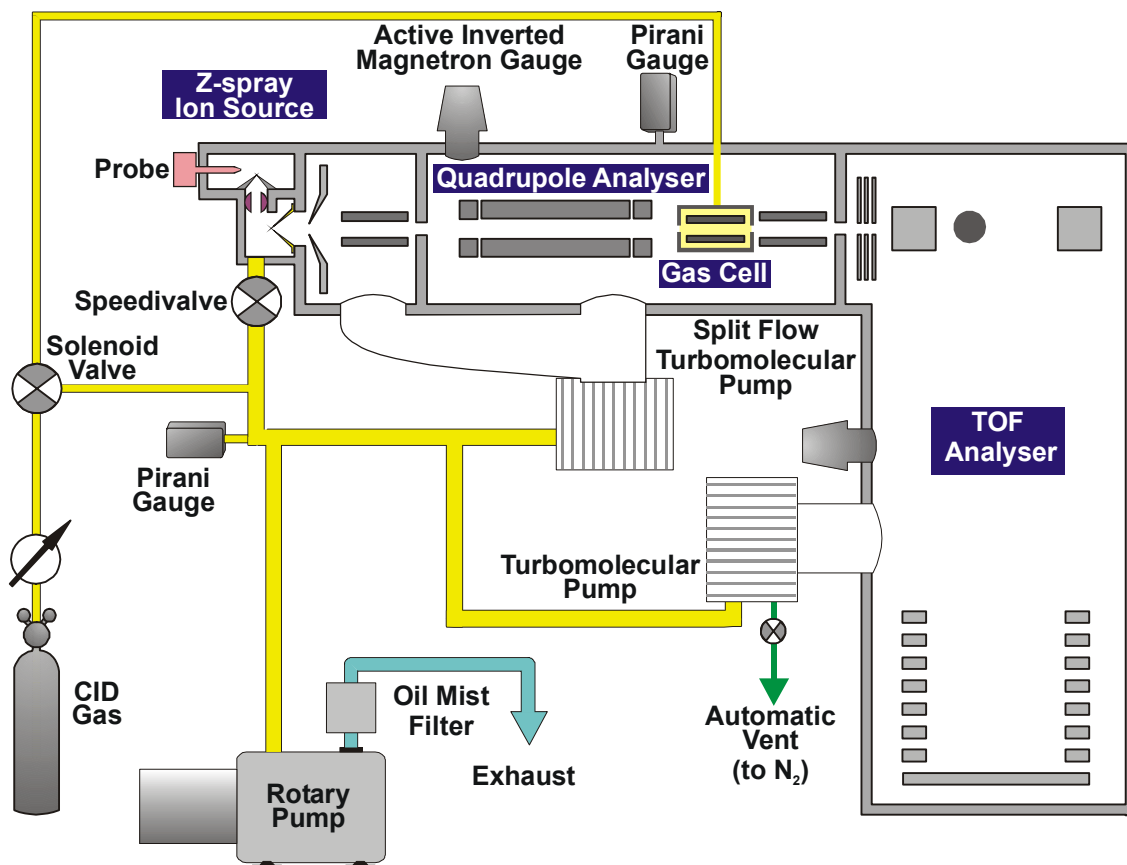


Figure 1-5 Q-ToF *micro* Vacuum System

1.4.1 Fine Pumping

Q-ToF micro is equipped with three water cooled turbomolecular pumps, providing independent fine pumping of the source hexapole, quadrupole and TOF analysers. Details of the operation and maintenance of the pumps can be found in the manufacturer's manuals provided.

1.4.2 Rotary Pumping

Source pumping and turbomolecular pump backing is by a direct drive rotary pump. The rotary pump is situated at the front of the instrument. Details of the operation and maintenance of the pump can be found in the manufacturer's manual provided.

1.4.3 Pressure Measurement

The backing pressure is monitored by an active Pirani gauge. The analyser and TOF pressures are monitored by active inverted magnetron (Penning) gauges. These gauges act as vacuum switches, switching the instrument out of Operate mode if the pressure is too high. Pressure readings may be displayed on the MassLynx NT tune page.

The analyser Penning gauge only comes on when the vacuum display window is open. At other times the gauge is off. The analyser Pirani gauge is used when the display is off, though no pressures are shown.

1.4.4 Vacuum Protection

The vacuum system is fully interlocked to provide adequate protection in the event of:

- a fault in the vacuum system.
- a failure of the power supply.
- a failure of the water supply.
- a vacuum leak.

1.5 Front Panel Connections

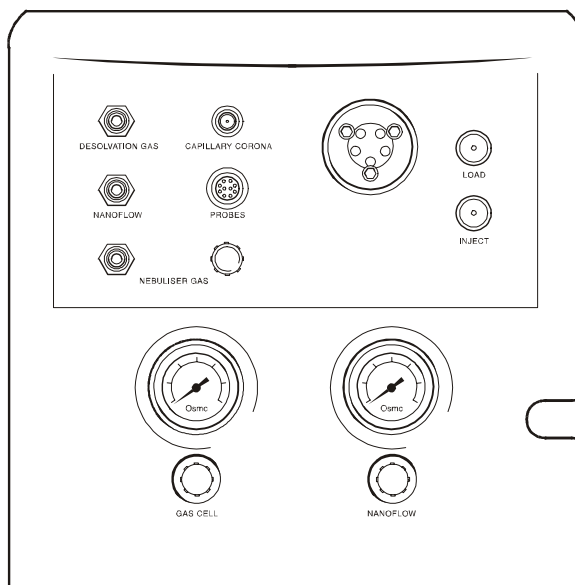


Figure 1-6 Front Panel

1.5.1 Desolvation Gas and Probe Nebuliser Gas

The PTFE gas lines for the Desolvation Gas and probe Nebuliser Gas are connected to the front of the instrument using threaded metal fittings. Cone Gas is connected internally.

1.5.2 High Voltage

The electrical connection for the ESI capillary or the APcI corona discharge pin is via the coaxial high voltage connector. This socket is labeled Capillary / Corona.

1.5.3 Heaters

The electrical connection for the APcI probe or the ESI desolvation heater is via the multi-way connector labeled Probes. This is removed from the front panel by pulling on the metal sleeve of the plug. Both the electrospray desolvation heater and the APcI probe heater use this connector.

The power for the source block heater is permanently connected. As a consequence, the source block assembly is usually very hot, and should not be touched.

1.5.4 Front Panel Controls and Indicators

Status Display

The display on the front panel of the instrument consists of two 3-colour light emitting diodes (LEDs).

The display generated by the Pump LED is dependent on the vacuum status of the instrument. The Operate LED depends on both the vacuum status and whether the operate mode has been selected from the Data System. Further information is included in Automatic Pumping and Vacuum Protection (see Routine Procedures).

Divert / Injection Valve

The divert / injection valve may be used in several ways depending on the plumbing arrangement:

- As an injection valve, with the needle port and sample loop fitted.
- As a divert valve, to switch the flow of solvent during a LC run.
- As a switching valve to switch, for example, between a LC system and a syringe pump containing calibrant.

This valve is pneumatically operated, using the same nitrogen supply as the rest of the instrument.

The two switches marked Load and Inject enable the user to control the valve when making loop injections at the instrument.

1.6 Rear Panel Connections

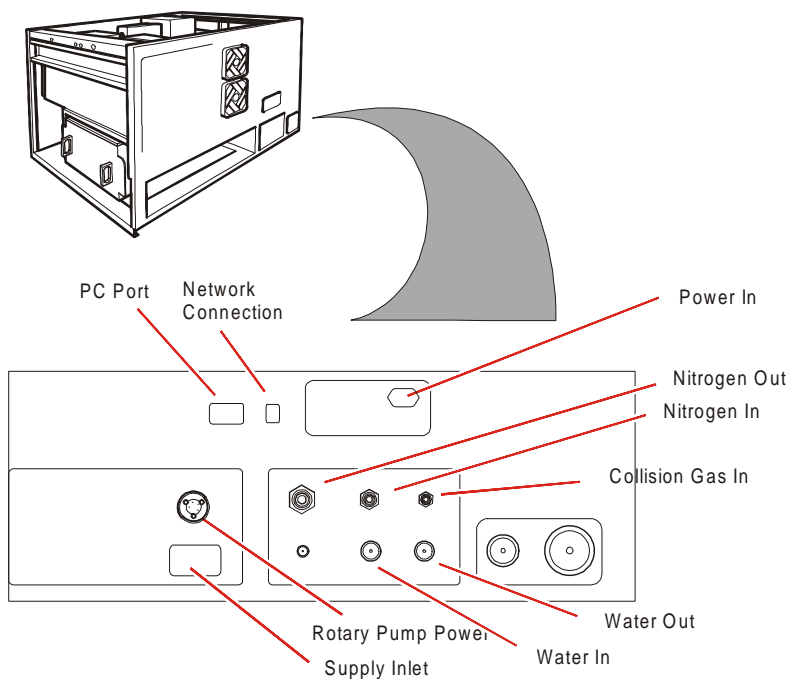


Figure 1-7 Rear Panel Connections

1.6.1 Water

Water is used to cool the turbomolecular pumps.

1.6.2 Nitrogen Gas In

The nitrogen supply (100 psi, 7 bar) should be connected to the Nitrogen Gas In push-in connector using 6mm PTFE tubing. If necessary this tubing can be connected to ¼ inch tubing using standard ¼ inch fittings.

Caution: Use only PTFE tubing or clean metal tubing to connect between the nitrogen supply and the instrument. The use of other types of plastic tubing will result in chemical contamination of the source.

1.6.3 Exhausts

The exhaust from the rotary pump should be vented to atmosphere outside the laboratory.

The gas exhaust, which also contains solvent vapours, should be vented via a separate fume hood, industrial vent or cold trap.

The gas exhaust should be connected using 10mm plastic tubing connected to the push-in fitting.

Caution: Do not connect these two exhaust lines together as, in the event of an instrument failure, rotary pump exhaust could be admitted into the source chamber producing severe contamination.

1.6.4 Supply Inlet

The mains power cord should be wired to a 230V mains outlet using a suitable plug, or to a transformer. For plugs with an integral fuse, the fuse should be rated at 13 amps (UK only).

1.6.5 Electronics

This circuit breaker switches power to the electronics. In the event of the instrument drawing more than the rated current, the circuit breaker will trip.

1.6.6 Rotary Pump

This circuit breaker switches power to the rotary and turbomolecular pumps. In the event of the pumps drawing more than the rated current, it will trip.

1.6.7 Event Out

Four outputs, Out 1 to Out 4 (Figure 1-8), are provided to allow various peripherals to be connected to the instrument. Switches S1 to S4 allow each output to be set to be either a contact closure (upper position) or a voltage output (lower position).

Out 1 and Out 2, when set to voltage output, each have an output of 5 volts. The voltage output of both Out 3 and Out 4 is 24 volts.

During a sample run an event output may be configured to close between acquisitions and is used typically to enable an external device to inject the next sample.

1.6.8 Contact Closure In

In 1 and In 2 inputs are provided ((Figure 1-8) to allow an external device to start sample acquisition once the device has performed its function (typically sample injection).

1.6.9 Analog Channels

Four analog channel (Figure 1-8) inputs are available, for acquiring simultaneous data such as a UV detector output. The input differential voltage must not exceed one volt.

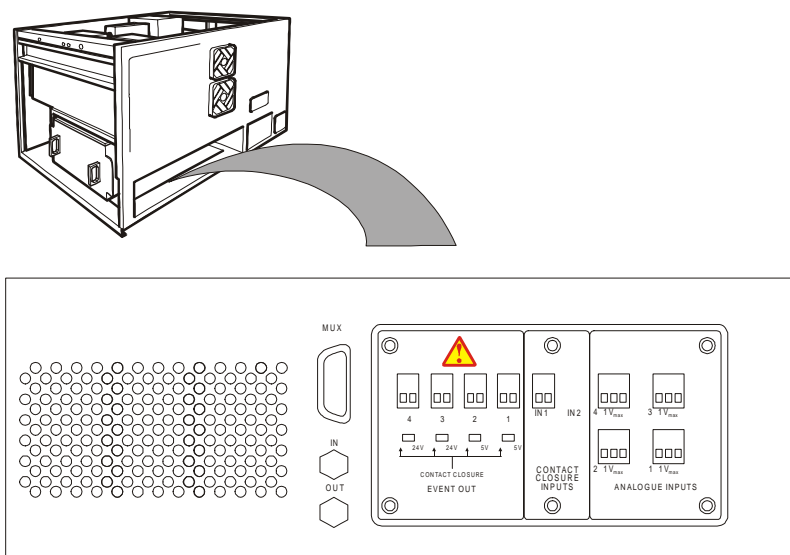


Figure 1-8 Rear Panel Outputs

1.7 MassLynx Data System

A PC computer runs the MassLynx NT software system to control Q-Tof micro, and to acquire and manipulate data from it. A high resolution colour monitor is also supplied.

Interaction with MassLynx NT is via the mouse and keyboard using menu-driven commands. Printing, file management and other routine procedures are performed using the appropriate Windows NT modules.

1.7.1 Software

The following software packages are supplied with Q-Tof micro:

- MassLynx NT.
- DataBridge, a utility to convert other format data files into MassLynx format.
- Microsoft Windows NT/2000/XP graphical environment.
- Mouse configuration.

A range of optional software modules for different applications is also available.

The MassLynx NT User's Guide describes the many facilities of the Micromass software. Documentation for the other software is also supplied.

Chapter 2

Routine Procedures

The following Routine Procedures are described in this chapter:

- Section 2.1, Start Up Following a Complete Shutdown
- Section 2.2, Start Up Following Overnight Shutdown
- Section 2.3, Automatic Startup and Shutdown

2

2.1 Start Up Following a Complete Shutdown

2.1.1 Preparation

If the instrument has been unused for a lengthy period of time, proceed as follows:

1. Check the level of oil in the rotary pump sight glass. Refill or replenish as necessary as described in the pump manufacturer's literature.
2. Connect a supply of dry, high purity nitrogen to the connector on the service panel at the rear of the instrument. Adjust the outlet pressure to 7 bar (100 psi).
3. Connect the water supply to the connections at the rear of the instrument.
4. Check that the rotary pump exhaust is connected to a suitable vent.



Attention: *Do not connect the two exhaust lines together. In the event of an instrument failure, rotary pump exhaust could be admitted into the source chamber, producing severe contamination*

5. Check that the instrument, data system and other peripheral devices (LC equipment, printer etc.) are connected to suitable mains supplies.
6. Check that the etherlink connection is made between the control PC and the instrument PC.
7. Switch on the host PC. Log on to Windows NT and wait for the system to boot up before the Q-ToF micro is switched on.
8. Switch on the mains to the mass spectrometer using the two circuit breakers situated on the service panel at the rear of the instrument.

9. Switch on the embedded PC using the switch at the rear of the instrument.
10. Log on to Micromass account (password analysis).

Note: Windows NT and MassLynx are configured to prevent unauthorized access.

11. On the host PC, select **Instrument >MS Tune** from the MassLynx Shortcut bar. This will invoke the Instrument Tune Page.

2.1.2 Pumping



Attention: To minimize wear to the lubricated components of the rotary pump, the manufacturers recommend that the pump is not started when the oil temperature is below 12°C.

1. Select **Vacuum** from the menu bar at the top of the tune page.
2. Click on **Pump**.
3. The rotary pump and the turbomolecular pumps start simultaneously.
4. The Vacuum LED on the front of the instrument shows amber as the system pumps down.

When the system has reached operating vacuum the Vacuum LED changes to a steady green.

5. If the rotary pump oil has been changed or replenished, open the gas ballast valve on the rotary pump. See the pump manufacturer's literature for details.

Rotary pumps are normally noticeably louder when running under gas ballast.

6. If opened, close the gas ballast valve when the rotary pump has run under gas ballast for 30 minutes.

2.1.3 MCP Detector Conditioning

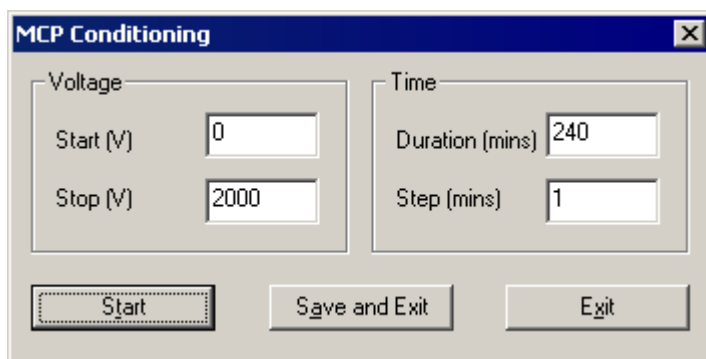


Figure 2-1 MCP Conditioning Dialog

The MCP detector must be conditioned before use, by gradually increasing the applied voltage over a long time period. This is necessary to allow escape of all absorbed water from within the microchannels.

Under normal operation the analyzer automatically vents to dry nitrogen. However, if the nitrogen supply was not connected to the instrument when last vented, or if the instrument has been left vented for more than one day, a significant amount of water vapor may have entered the analyzer. Under these circumstances it is good practice to allow the instrument to pump for 12 hours before commencing the conditioning process.

In all cases, the TOF pressure must be $2\text{e-}6$ mbar prior to commencing MCP conditioning.

MCP conditioning should be repeated after every instrument venting.

Note: *It is not necessary to recondition the detector if the instrument has been left out of the operate mode while still under vacuum.*

Note: *During routine cleaning of the source sample cone, the source isolation valve is closed in order to maintain analyzer vacuum. It is not, therefore, necessary to recondition the detector after this procedure.*

The procedure for MCP conditioning is as follows:

1. Ensure that the TOF pressure is $2\text{e-}6$ mbar.
2. Check that the MCP Detector voltage is set to zero on the tune page.
3. Switch the instrument into Operate.

4. Select **Options > MCP Conditioning** to access the MCP conditioning program (Figure 2-1).
5. Set Start to 100V, Stop to 2700V, Duration to 600 minutes and Step to 5 minutes.
6. A 'quick condition' may be performed following brief venting, after source cleaning for example.
7. Set Start to 100V, Stop to 2700V, Duration to 180 minutes and Step to 1 minute.



Attention: Failure to follow the recommended MCP conditioning procedure can severely reduce detector lifetime.

2.1.4 Instrument Warm-up

Switch the instrument into the operate mode by selecting Operate on the MassLynx tune page.

For the best mass accuracy to be obtained the instrument temperature must be stabilized for a minimum of two hours after switching into operate.

Note: *Leaving the instrument continuously in operate does not shorten the detector lifetime. It is recommended that the instrument is left in operate at all times (except of course during maintenance procedures) in order to reduce mass scale drifts due to temperature changes. Switching the instrument out of operate mode overnight is not necessary.*

2.1.5 Using the Instrument

The Q-ToF micro is now almost ready to use. To complete the start up procedure and prepare for running samples, follow the instructions in Start Up Following Overnight Shutdown in the following pages.

2.2 Start Up Following Overnight Shutdown

The instrument will have been left in the operate mode under vacuum.

It is recommended that the data system is left on overnight. However, if the data system has been switched off, switch it on as described in the preceding section.

2.2.1 Preparation for Electrospray Operation

If the corona discharge pin is fitted, proceed as follows:

1. Deselect **Operate** from the tune page to put the instrument into standby mode.
2. Disconnect the gas and electrical connections from the front panel.
3. Unscrew the probe thumb nuts and remove the probe.
4. Undo the three thumb screws and remove the probe adjustment flange and glass tube.



Caution: The ion source block can be heated to temperatures of 150°C, and will be maintained at the set temperature when the source enclosure is removed. Touching the ion block when hot may cause burns to the operator

5. Disconnect the APcI high voltage cable from the socket positioned at the bottom right corner of the source flange.
6. Remove the corona discharge pin from its mounting contact, and fit the blanking plug.
7. Replace the glass tube and adjustment flange.
8. Ensure that the source enclosure is in place.

Note: The Z-spray source enclosure consists of the glass tube and the probe adjustment flange.



Caution: Operating the source without the source enclosure will result in solvent vapor escape and the exposure of hot surfaces and high voltages.

9. With the corona discharge pin removed, the plug fitted and the source enclosure in place, proceed as follows:
10. Connect the source's gas line to **Desolvation Gas** on the front panel. Tighten the nut to ensure a good seal.
11. Check that the lead of the probe adjustment flange is plugged into the socket labelled **Probes** on the front panel.
12. Connect the electrospray probe's gas line to **Nebuliser Gas** on the front panel.
13. Connect the liquid flow of a LC system or syringe pump to the probe.
14. Insert the probe into the source and tighten the two thumb nuts to secure the probe firmly.

15. Plug the probe lead into **Capillary / Corona** on the front panel.
16. If necessary, change the ionization mode using the Ion Mode command.
17. Set **Source Block Temp** to 100°C and **Desolvation Temp** to 120°C.



Attention: The maximum operating temperature for the source heater is 150°C. Do not set **Source Block Temp** higher than 150°C.

2.2.2 Preparation for APcI Operation

If the corona discharge pin is not fitted, proceed as follows:

1. Deselect **Operate** from the tune page to put the instrument into standby mode.
2. Disconnect the gas and electrical connections from the front panel.
3. Unscrew the probe thumb nuts and remove the probe.
4. After a period of operation at high flow rates, allow the glass source enclosure to cool before removal.
5. Undo the three thumb screws and remove the probe adjustment flange and glass tube.



Attention: The ion source block can be heated to temperatures of 150°C, and will be maintained at the set temperature when the source enclosure is removed. Touching the ion block when hot may cause burns to the operator.

6. Remove the blanking plug from the discharge pin mounting contact and fit the corona discharge pin, ensuring that the tip is in line with the tip of the sample cone.
7. Connect the APcI high voltage cable between Capillary / Corona and the socket positioned at the bottom left corner of the source flange.
8. Replace the glass tube, adjustment flange and moulded cover.



Caution: Operating the source without the source enclosure will result in solvent vapor escape and the exposure of hot surfaces and high voltages.

9. With the corona discharge pin fitted and the source enclosure in place, proceed as follows:
10. Insert the APcI probe into the source and tighten up the two thumb screws.
11. If necessary, change the ionization mode using the **Ion Mode** command.

12. Set **Source Temp** to 150°C.



Attention: The maximum operating temperature for the source heater is 150°C. Do not set **Source Block Temp** higher than 150°C.

Note: Do not start the liquid flow until the gas flow and probe heater are switched on with the probe inserted.

2.2.3 Transient Pressure Trip

The transient trip is designed to protect the instrument from potentially damaging pressure surges and operates routinely whenever the pressure rises.

Should the vacuum gauge(s) detect a pressure surge above the preset trip level (normally set at 1e-5 mbar by software) the following events occur:

- The green Pump lamp becomes amber.
- If in the operate mode, the system turns off the critical source, analyzer and detector voltages, and the green Operate lamp becomes amber.
- Acquisition continues though, of course, no real data are recorded.

When the vacuum recovers:

- The amber Pump lamp becomes green.
- If previously in the operate mode, voltages are restored and Operate reverts to green.

The period during which the trip was operative will appear in a raw total ion chromatogram as a period of reduced baseline noise.

Further deterioration of the system pressures results in a “vacuum fault” condition and the system is shut down (see below).

2.2.4 Power Failure

In the event of an unexpected failure of the electrical supply the instrument is vented safely. If power is unlikely to be restored quickly, follow the shutdown procedure described later in this chapter. When power is restored follow the startup procedure.

Should the power fail and then be restored while the instrument is unattended, the system will continue to vent, and will require to be pumped down in accordance with the start-up procedure.

2.2.5 Nitrogen Supply

Replacement of nitrogen cylinders should be conducted in accordance with the operation, handling and storage instructions provided by the local gas supplier.

Toggle the API gas button to Off, to close the nitrogen inlet valve prior to disconnecting the supply.

Set the nitrogen inlet pressure to 7 bar (100 psi).



Attention: Under no circumstances should the nitrogen pressure exceed 10 bar (140 psi).

2.3 Automatic Startup and Shutdown

MassLynx comes with automatic Startup and Shutdown files. They are found in the C:\MassLynx\Shutdown directory and are called ShutDownxxx.acl and StartUpxxx.acl where xxx refers to the instrument configuration. E.g. ShutDownESI_ACE.acl for an instrument configured as an ACE system.







2.3.1 The Shutdown Editor

The shutdown editor allows the automatic startup and shutdown procedures to be modified or new procedures to be created. These can be run automatically before or after a batch if the relevant boxes are checked in the Batch Control frame of the Shutdown Page (Section 2.3.3).

Select **Edit Shutdown** or **Startup** from the MassLynx **Instrument Shortcut Bar**.

This Invokes the Shutdown Editor which has two tabbed pages, the Shutdown Page (described in Section 2.3.3 on page 29) and the Auto Control Tasks Page (described in Section 2.3.4 on page 33)

2.3.2 The Shutdown Editor Toolbar

Toolbar Button	Menu Equivalent	Purpose
	File... New	Create a new Startup or Shutdown file
	File... Open	Open an existing Startup or Shutdown file
	File... Save or Save As	Save a Startup or Shutdown file
	File... Print	Print a Startup or Shutdown file
	Control List..Run List	Run a Startup or Shutdown file
	Control List...Stop List	Stop a Startup or Shutdown file

2.3.3 Shutdown Page

Select **Edit Shutdown** or **Startup** from the MassLynx **Instrument Shortcut Bar**.

By default the Shutdown Page is displayed (Figure 2-2) with the following parameters.

Enable startup before batch	Enables/disables the running of a task file before the start of a batch of samples. Startup Enabled will appear on the right of the Status bar on the bottom of the MassLynx window when checked.
Startup file	The file name of the task file that is run before the start of a batch.
Browse startup file	Brings up a file dialog allowing a task file to be selected for the startup before batch procedure.
Enable shutdown after batch	Enables/disables the running of a task file after the end of a batch of samples. Shutdown Enabled will appear on the right of the Status bar on the bottom of the MassLynx window when checked.

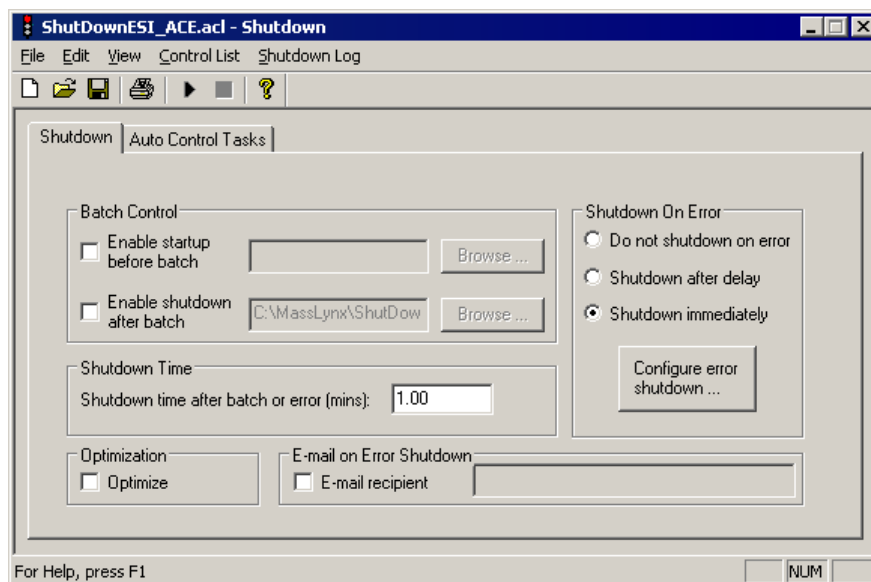


Figure 2-2 Shutdown Page

Shutdown file	The file name of the task file that is run after the end of a batch.
Browse shutdown file	Brings up a file dialog allowing a task file to be selected for the shutdown after batch procedure.
Shutdown Time	This is the delay (in minutes) between the batch finishing and the shutdown procedure initiating.
Shutdown On Error	This allows shutdown on error to be enabled/disabled and to determine whether the shutdown tasks should be initiated immediately, or after the time entered in the Shutdown Time edit box.
Configure Error Shutdown	Brings up the Shutdown on Error Configuration dialog (see Section on page 30).
Optimization	Check this box to optimize the procedure.
E-mail on Shutdown	Check box and enter the e-mail address for shutdown information to be sent.

Configure Shutdown on Error

Clicking on the Configure Error Shutdown button displays the following dialog.

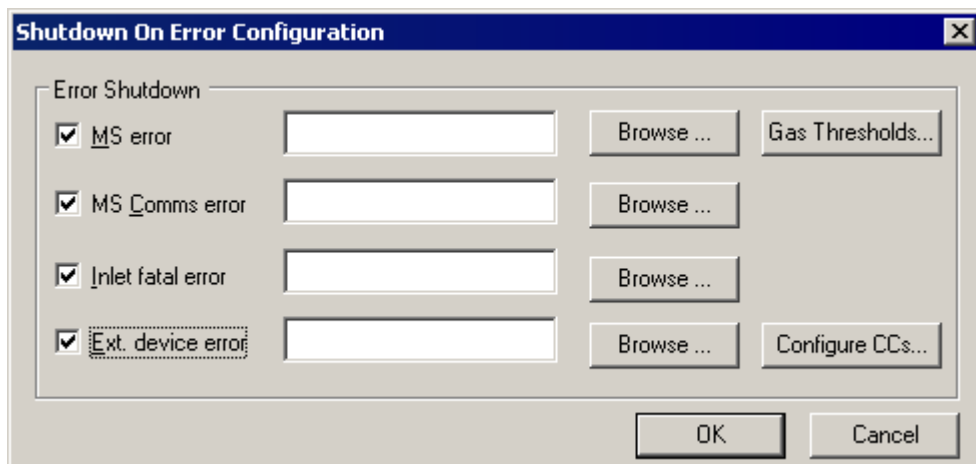


Figure 2-3 Shutdown on Error Configuration Dialog

MS error	Enables/disables running a <i>task file</i> when an <i>MS error</i> occurs.
MS error shutdown file	The <i>task file</i> to run on an <i>MS error</i> .
MS error browse file	Allows user to browse for <i>task files</i> .
MS Comms Error	Enables/Disables running task file when MS comms error occurs
MS comms error shutdown file	The <i>task file</i> to run on an <i>MS comms error</i> .
Inlet fatal error	Enables/disables running a task file when a <i>LC error</i> occurs.
Inlet error shutdown file	The <i>task file</i> to run when a <i>LC error</i> occurs.
Inlet error browse file	Allows user to browse for <i>task files</i> .
Ext. device error	Enables/disables running a <i>task file</i> when an <i>external device error</i> occurs.
Ext. device shutdown file	The <i>task file</i> to run when an <i>external device error</i> occurs.
Ext. device browse file	Allows user to browse for <i>task files</i> .
Gas Threshold	Invokes the Gas Threshold dialog (see Section on page 32)
Configure CCs	Invokes the Configure Event In dialog (see Section on page 31).

Configure Contact Closures

Clicking on the **Configure CCs** displays the following dialog.

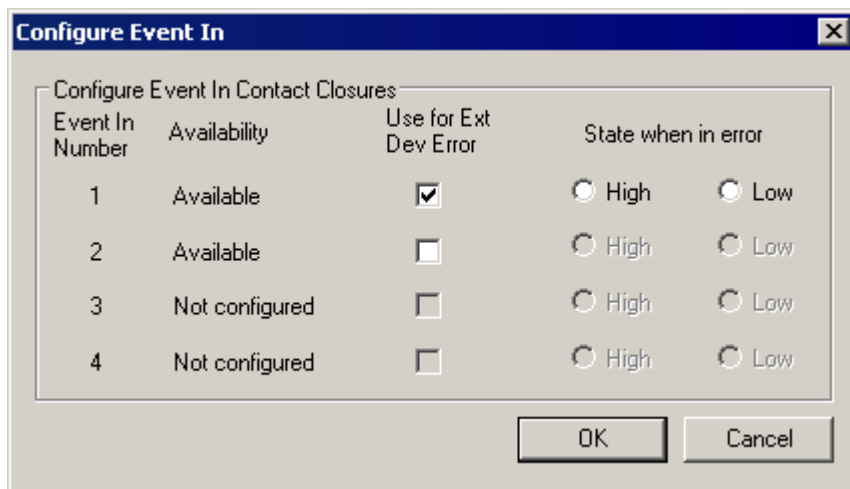


Figure 2-4 Configure Event In Dialog

This window is used to configure the Event In contact closures on the back of the MS. These can be used to detect errors in external devices enabling the External device error flag to be set. Only contact closures not used to signal the completion of an injection, or not used by MUX systems, are enabled. Each of the CCs available has the following controls.

Event In Number
Availability

Relates to the *Event In* number on the rear of the MS. This can either be '*Used By Inlet*', '*Available*', '*Used by MUX*' or '*Not Configured*'. The controls are only enabled if the CC is '*Available*'.

Use For Ext. Dev. Error

Enables/disables the use of this contact closure to signal an error in an external device.

State When in Error

Determines what state the CC will be in when an error is present in an external device.

Gas Thresholds

Clicking on the **Gas Thresholds** button displays the following dialog.

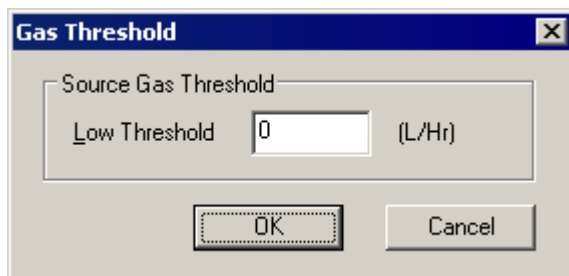


Figure 2-5 Gas Threshold Dialog

Clicking **Gas Thresholds** on the Configure Shutdown on Error dialog (Figure 2-3) invokes the Gas Threshold dialog; an MS error is flagged if the source gas flow rate falls below the value entered in the **Low Threshold... (L/Hr)** text box.

Note: The Gas Thresholds button is only available when the MS error option is selected.

2.3.4 The Auto Control Tasks Page

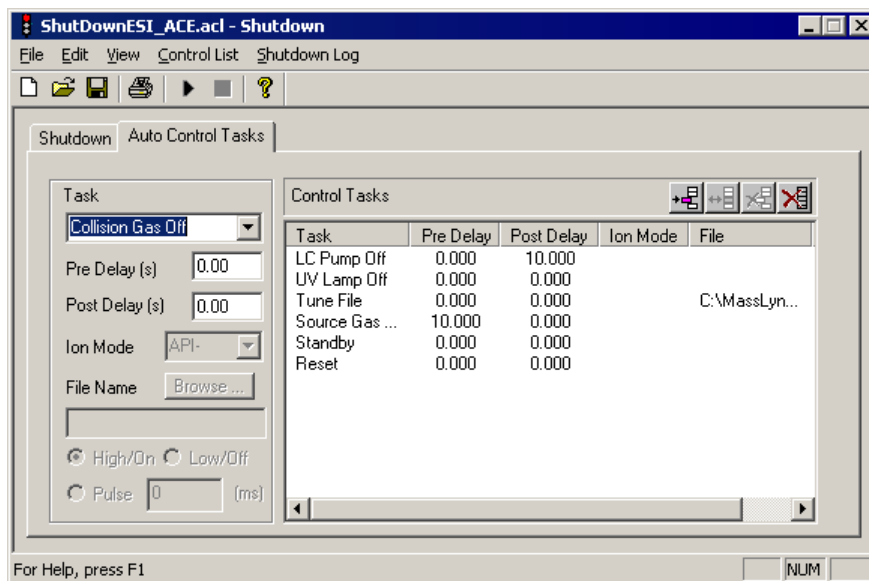


Figure 2-6 Auto Control Tasks Page

Task
Pre-Delay

This is a dropdown list box with all the available tasks.
This is the length of time that will elapse before the current task is performed

Post-Delay

This is the length of time that will elapse after the current task has been completed and before the next task is started.

E.g. a Post delay of 60s, in the Tune File task above, means that there will be a delay of 60 seconds before the next task is started, to allow the machine to stabilize with the new tune page settings.


Ion Mode

This is a dropdown list box with all the available ionization modes.

File Name

This is the name of the Tune file to be used. The file name can be typed in, including the full path name, or selected from the browser displayed when the Browse button is pressed.

To Add a Task


1. Select a task from the dropdown Task list box.
2. Enter the required parameters.
3. Press the add  button.

Note: If this is a new task timetable the task will be added to the end of the list. If a task has been inserted into the task timetable then all subsequent tasks will be added after the inserted task. To add a task to the end of the timetable after inserting a task, double click below the last entry in the timetable and then add the new task.

To Insert a Task


1. Click on the entry in the task timetable before which you want to insert the new task.
2. Select a task from the dropdown Task list box.
3. Enter the required parameters.
4. Press the add button. The task will be inserted before the selected entry.

To Modify a Task


1. Click on the entry in the task timetable. The details for the task will be displayed in the fields on the left of the screen.
2. Change the required parameters.
3. Press the modify  button. The details will change in the task timetable.

To Delete a Task


1. Click on the entry in the task timetable. The details for the task will be displayed in the fields on the left of the screen.

2. Press the add  button. The task selected will be deleted from the task timetable.

To Delete All Tasks


Press the add  button. All tasks will be deleted from the task timetable.

To Change the Width of a Column

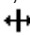
The width of the columns can be changed, by positioning the mouse pointer on the heading between two columns until the  symbol appears, and then click and drag until the column is the required width.

Saving/Loading Startup and Shutdown Files


To Open a Startup or Shutdown file

1. Press  or select Open from the File menu. This displays the Open file dialog.
2. Select a data file and press the Open button.

To Save a Startup or Shutdown file

1. Press  or select Save or Save As from the File menu. If this is a new file, or the Save As option has been selected, the Save As dialog is displayed
2. Type a name into the File Name box and press the Save button.

Printing Startup and Shutdown Files

1. Press  or select **Print** from the File menu. This displays the Print dialog.
2. Select the printer, print range and number of copies and press the OK button.


Creating Startup and Shutdown Files

Press or select **New** from the **File** menu.

Running Startup and Shutdown Files

If Startup or Shutdown is selected from the MassLynx Short cut bar or from the Shutdown editor Control List menu then the automatic startup and shutdown files are run.

To Run a Different Startup or Shutdown file: -

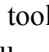1. Open the required file in the Shutdown editor and press the  toolbar button or select Run List from the Shutdown editor Control List menu.

2. Press the toolbar button or select Stop List from the Shutdown editor Control List menu if you wish to stop running this file.

Alternatively if the Enable startup before batch or Enable shutdown after batch options are checked on the Shutdown Page (Section 2.3.3) the files will be run before or after a Sample List run.

2.3.1 Shutdown Log

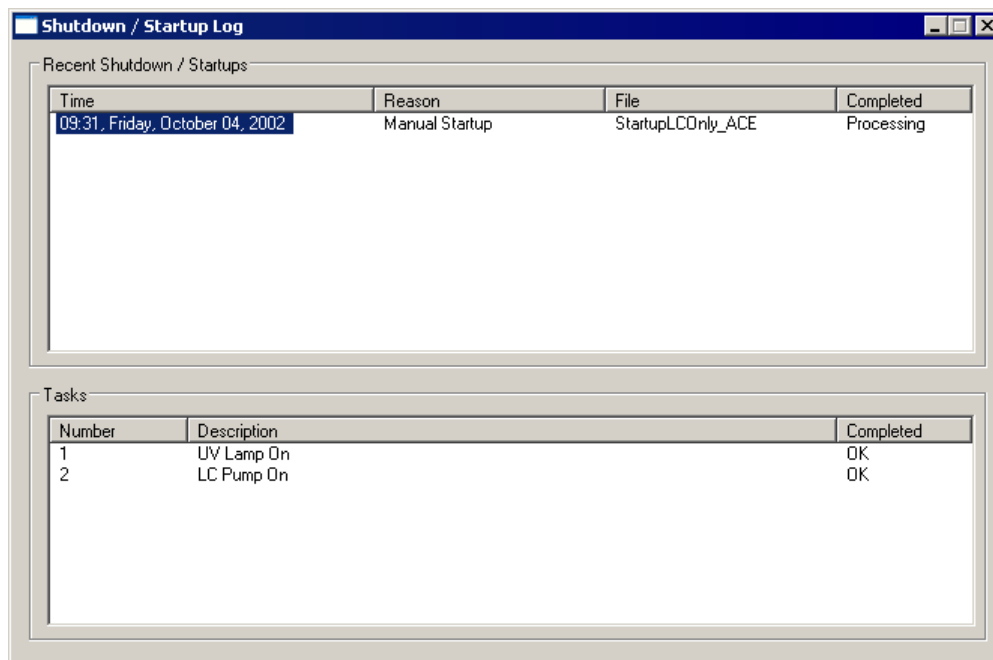


Figure 2-7 Shutdown Startup Log

Shutdown / Startup Log

The Shutdown / Startup Log keeps a record of the most recent startups and shutdowns. Select **Shutdown Log > Recent Shutdowns and Startups** invokes the **Shutdown Startup Log** (Figure 2-7). Selecting the most **Recent Shutdowns and Startups** in the top pane will show the tasks carried out in the bottom pane.

Log Parameters

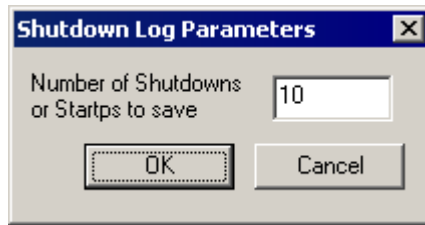


Figure 2-1 Log Parameters Dialog

Selecting **Shutdown Log > Log Parameters** from the main menu bar invokes the Shutdown Log Parameters dialog. The number of startups and shutdowns recorded in the Startup / Shutdown Log can be altered by changing the number in the text box.

Chapter 3

Electrospray

3.1 Introduction

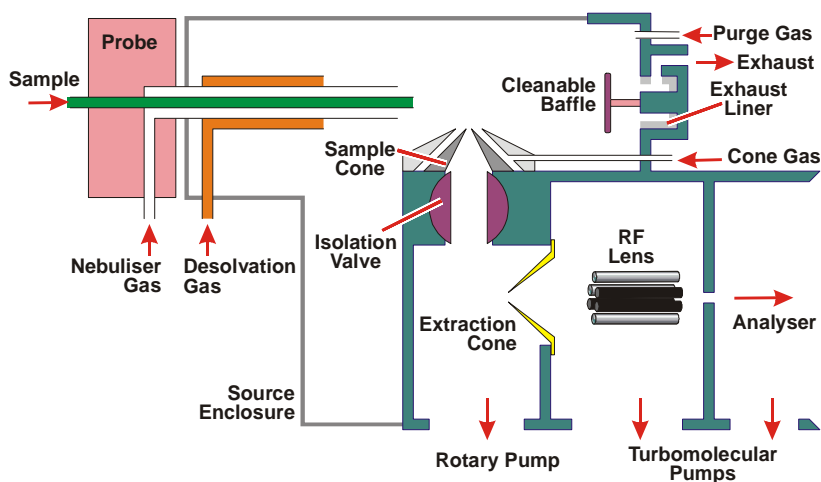


Figure 3-1 Electrospray Scheme

The ESI interface consists of the standard Z-spray source fitted with an electrospray probe. See the following chapter for additional information concerning the optional nanoflow interface.

Mobile phase from the LC column or infusion pump enters through the probe and is pneumatically converted to an electrostatically charged aerosol spray. The solvent is evaporated from the spray by means of the desolvation heater. The resulting analyte and solvent ions are then drawn through the sample cone aperture into the ion block, from where they are then extracted into the analyzer.

The electrospray ionization technique allows rapid, accurate and sensitive analysis of a wide range of analytes from low molecular weight (less than 200 Da) polar compounds to biopolymers larger than 100 kDa.

Generally, compounds of less than 1000 Da produce singly charged protonated molecules ($[M+H]^+$) in positive ion mode. Likewise, these low molecular weight analytes yield ($[M-H]^-$) ions in negative ion mode, although this is dependent upon compound structure.

High mass biopolymers, for example peptides, proteins and oligonucleotides, produce a series of multiply charged ions. The acquired data can be transformed by the data system to give a molecular weight profile of the biopolymer.

The source can be tuned to fragment ions within the ion block. This can provide valuable structural information for low molecular weight analytes.

The most common methods of delivering sample to the electrospray source are:

Syringe Pump and Injection Valve

A flow of mobile phase solvent passes through an injection valve to the electrospray source. This is continuous until the pump syringes empty and need to be refilled. Sample is introduced through the valve injection loop (usually 10 or 20 μ l capacity) switching the sample plug into the mobile phase flow. Tuning and acquisition are carried out as the sample plug enters the source. (At a flow rate of 10 μ l/min a 20 μ l injection lasts 2 minutes.)

Reciprocating Pump and Injection Valve

A flow of mobile phase solvent passes through an injection valve to the electrospray source. Sample injection and analysis procedure is the same as for the syringe pump. The pump reservoirs are simply topped up for continuous operation. The most suitable reciprocating pumps for this purpose are those which are specified to deliver a flow between 1 μ l/min and 1 ml/min. A constant flow at such rates is more important than the actual flow rate. The injection valve on reciprocating pumps may be replaced by an autosampler for unattended, overnight operation.

Infusion Pump

The pump syringe is filled with sample in solution. The infusion pump then delivers the contents of the syringe to the source at a constant flow rate. This arrangement allows optimization and analysis while the sample flows to the source at typically 5-30 μ l/min. Further samples require the syringe to be removed, washed, refilled with the next sample, and replumbed.

A 50:50 mixture of acetonitrile and water is a suitable mobile phase for the syringe pump system and the reciprocating pump systems. This is appropriate for positive and negative ion operation.

Positive ion operation may be enhanced by 0.1 to 1% formic acid in the sample solution.

Negative ion operation may be enhanced by 0.1 to 1% ammonia in the sample solution. Acid should not be added in this mode.

These additives should not be used in the mobile phase for flow injection analysis (FIA) studies, to allow easy change over between positive and negative ion analysis.

Degassed solvents are recommended for the syringe and reciprocating pumps. Degassing can be achieved by sonification or helium sparging. The solvents should be filtered, and stored under cover at all times.

It is wise periodically to check the flow rate from the solvent delivery system. This can be carried out by filling a syringe barrel or a graduated glass capillary with the liquid emerging from the probe tip and timing a known volume, say 10 μ l. Once the rate has been measured and set, a note should be made of the back pressure readout on the pump as fluctuation of this reading can indicate problems with the solvent flow.

3.1.1 Post-column Splitting

Although the electrospray source can accommodate flow rates up to 1 ml/min, it is recommended that the flow is split post-column to approximately 200 μ l/min. Also, even at lower flow rates, a split may be required to save valuable samples.

The post-column split consists of a zero dead-volume tee piece connected as shown in Figure 3-2.

The split ratio is adjusted by increasing or decreasing the back pressure created in the waste line, by changing either the length or the diameter of the waste tube. A UV cell may also be incorporated in the waste line, avoiding the requirement for in-line, low volume “Z cells”. As the back pressure is varied, the flow rate at the probe tip should be checked as described above.

These principles apply to splitting for both megafLOW and normal flow electrospray.

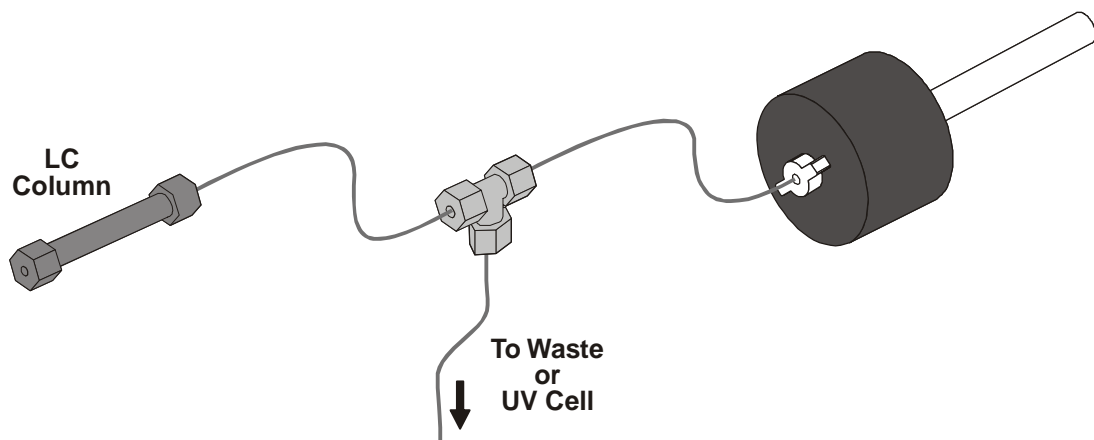


Figure 3-2 Post Column Split

3

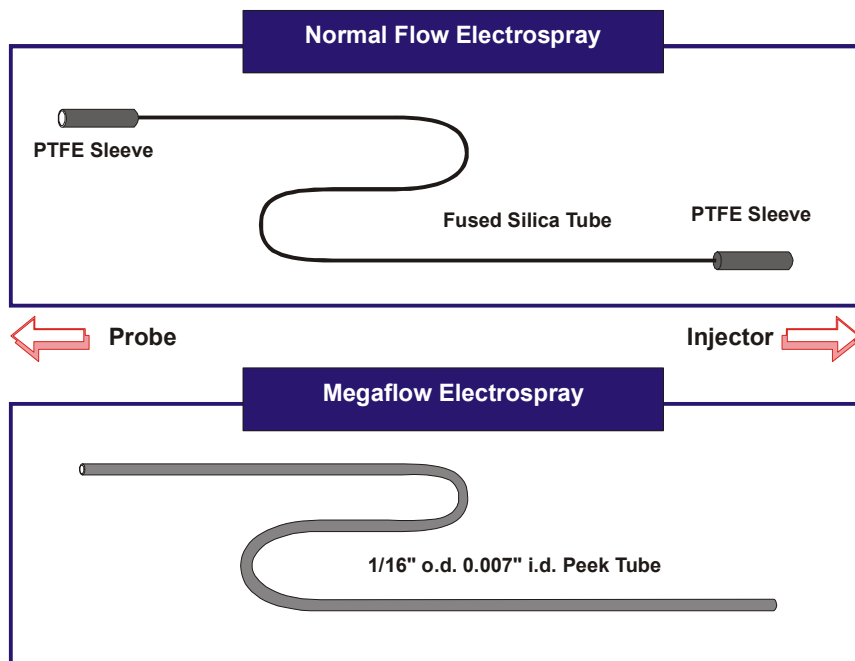


Figure 3-3 Comparison of Normal and Megaflow Electrospray

3.1.2 Megaflow

Megaflow electrospray (Figure 3-3) enables flow rates from 200 $\mu\text{L}/\text{min}$ to 1 mL/min to be accommodated. This allows microbore (2.1mm) or 4.6mm diameter columns to be interfaced without splitting.

Changing Between Flow Modes

When changing between megaflow and standard electrospray operation, it is essential that the correct tubing is used to connect the probe to the sample injector. For megaflow operation 1/16" o.d., 0.007" i.d. peek tubing, easily identified by its yellow stripe, is used. This replaces the standard fused silica tube, together with the PTFE sleeves.

3.2 Operation

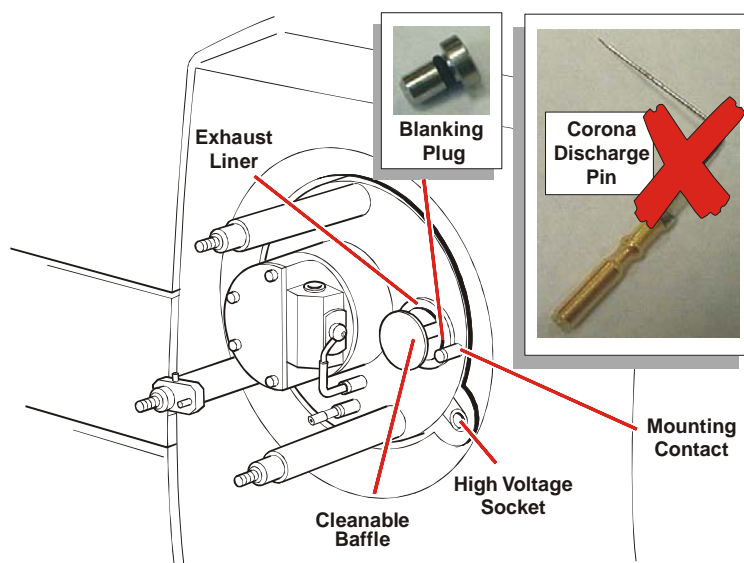


Figure 3-4 ESI Source



Caution: Operating the source without the source enclosure will result in solvent vapor escape and the exposure of hot surfaces and high voltages.

Ensure that the source is assembled as described in Maintenance and Fault Finding, and that the instrument is pumped down and prepared for electrospray operation as described in Routine Procedures.

Ensure that a supply of nitrogen has been connected to the gas inlet at the rear of the instrument and that the head pressure is between 6 and 7 bar (90-100 psi).

Ensure that the exhaust liner and the cleanable baffle are fitted to the source.

This is important for optimum electrospray intensity and stability when operating at low flow rates.

Checking the ESI Probe

1. Connect the electrospray probe to a pulse free pump.
2. Solvent should be degassed to prevent beam instabilities caused by bubbles.
3. Connect the PTFE tubing of the electrospray probe to Nebuliser Gas on the front panel. Secure with the nut provided.
4. With the probe removed from the source turn on the liquid flow at 10 $\mu\text{l}/\text{min}$ and check that liquid flow is observed at the tip of the capillary.
5. To avoid unwanted capillary action effects, do not allow liquid to flow to the probe for long periods without the nitrogen switched on.
6. Turn on **Nitrogen** by selecting **API Gas**, and check that a nebuliser flow of less than 100 litres/hour is registered.
7. To monitor the flow rate, select Window then Gas Flow on the tune page and observe the readback window.
8. Check that there is gas flow at the probe tip and ensure that there is no significant leakage of nitrogen elsewhere.
9. Adjust the probe tip to ensure complete nebulization of the liquid.

There should be approximately 0.5 mm of sample capillary protruding from the nebulizing capillary (Figure 3-5).

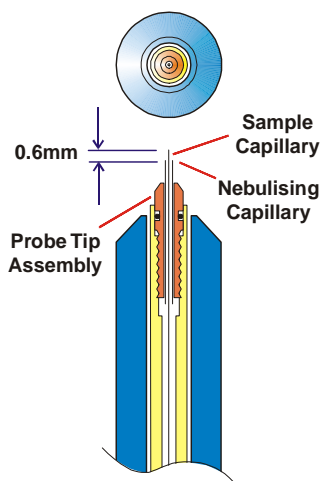


Figure 3-5 Electro spray Tip

The tip of the electro spray probe can influence the intensity and stability of the ion beam. A damaged or incorrectly adjusted probe tip will lead to poor electrospray performance.

10. Using a magnifying glass ensure that both inner and outer stainless steel capillaries are straight and circular in cross-section.
11. Ensure that the inner stainless steel capillary is coaxial to the outer capillary.

If the two capillaries are not coaxial, it is possible to bend the outer capillary slightly using thumbnail pressure.

12. Insert the probe into the source and tighten the two thumb screws.
13. Plug the probe high voltage cable into **Capillary / Corona** on the front panel.

Obtaining an Ion Beam

1. If necessary, change the ionization mode using the **Ion Mode** command.
2. Using the needle valve on the front panel, set the desolvation gas flow rate to 300 litres/hour.
3. Turn on the liquid flow at 10 µl/min and set **Desolvation Temp** to 100°C.

Tuning and Optimization

The following parameters, after initial tuning, should be optimized using a sample representative of the analyte to be studied. It will usually be found, with the exception of the sample cone voltage, that settings will vary little from one analyte to another.

Probe Position

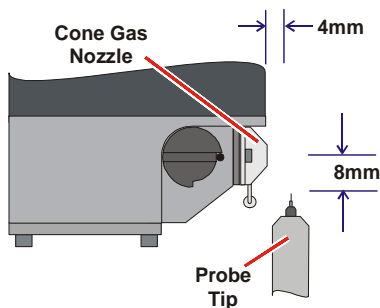


Figure 3-6 Optimum Probe Position

The position of the probe (Figure 3-6) is adjusted using the probe adjustment collar (in and out) and the adjustment knob (sideways) located to the left of the probe (Figure 3-7). The two screws can be adjusted singly or simultaneously to optimize the beam. The position for optimum sensitivity and stability for low flow rate work (10 $\mu\text{l}/\text{min}$) is shown.

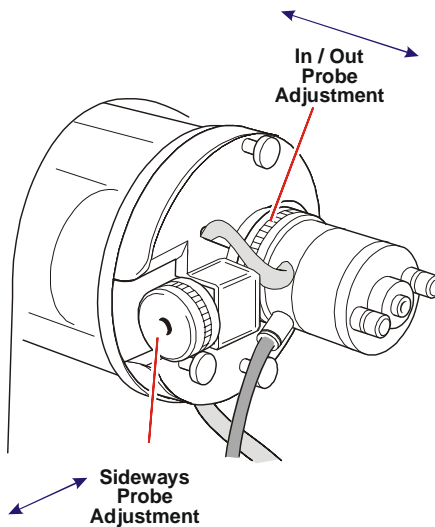


Figure 3-7 Probe Adjustment

Small improvements may be gained by varying the position using the sample and solvent system under investigation. The following information should be considered when setting the probe position:

- 10 mm of movement is provided in each direction, with 1.25mm of travel per revolution of the probe positioning controls.
- At higher liquid flow rates the probe tip should be positioned further away from the sample cone to achieve optimum stability and sensitivity. The position is less critical than at lower flow rates.

Nebuliser Gas

Optimum nebulization for electrospray performance is achieved with a nitrogen flow between 70 and 90 litres per hour. This can be achieved by fully opening the Nebuliser flow control valve, which is situated on the instrument's front panel.

Desolvation Gas

The desolvation gas, also nitrogen, is heated and delivered as a coaxial sheath to the nebulized liquid spray by the desolvation nozzle.

Note: The position of the desolvation nozzle heater is fixed relative to the probe tip and requires no adjustment.

The **Desolvation Gas** flow rate is adjusted by the control value situated on the instrument's front panel. The optimum **Desolvation Temp** and flow rate is dependent on mobile phase composition and flow rate. A guide to suitable settings is given below.

The **Desolvation Gas** flow rate indicated on the MassLynx tune page represents total drying flow, that is desolvation gas + cone gas (nanoflow only) + purge gas (if enabled).

Solvent Flow Rate $\mu\text{l/min}$	Desolvation Temp $^{\circ}\text{C}$	Desolvation Gas Flow Rate Litres/hour
<10	100 - 120	200 - 250
10 - 20	120 - 250	200 - 400
20 - 50	250 - 350	200 - 400
>50	350 - 400	500 - 750

Higher desolvation temperatures give increased sensitivity. However increasing the temperature above the range suggested reduces beam stability. Increasing the gas flow rate higher than the quoted values leads to unnecessarily high nitrogen consumption.



Attention: Do not operate the desolvation heater for long periods of time without a gas flow. To do so could damage the source.

Cone Gas

The cone gas reduces the intensity of solvent cluster ions and solvent adduct ions. The cone gas flow rate should be optimized by increasing until solvent cluster ions and / or adduct ions are reduced as much as possible without diminishing the intensity of the ion of interest, normally (M+H)⁺.

Typical cone gas flow rates are in the range 100 to 300 litres per hour.

Purge Gas

The purge gas is not necessary for most ESI applications. It may be useful for megafLOW operation where an analyte is susceptible to acetonitrile adducting.

Purge gas is enabled simply by removing the blanking plug from the outlet situated within the source enclosure.

Purge gas flow rate is a constant fraction (30%) of the total desolvation gas flow.

Source temperature

100°C is typical for 50:50 CH₃CN:H₂O at solvent flow rates up to 50 µl/min. Higher source temperatures, up to 150°C, are necessary for solvents at higher flow rates and higher water content.



Attention: The maximum operating temperature for the source heater is 150°C. Do not set **Source Temp** higher than 150°C

Capillary Voltage

Capillary usually optimizes at 3.0kV, although some samples may tune at values above or below this, within the range 2.5 - 4.0kV for positive electrospray. For negative ion operation a lower voltage is necessary, typically between 2.0 - 3.5kV.

At high flow rates this parameter may optimize at a value as low as 1kV.

Sample Cone Voltage

A **Cone** setting between 25V and 70V will produce ions for most samples, although solvent ions prefer the lower end and proteins the higher end of this range. Whenever sample quantity and time permit, Cone should be optimized for maximum sensitivity, within the range 15V to 150V. Increasing Cone will increase ion fragmentation within the source.

Extraction Cone Voltage

Extractor optimizes at 0 - 5V. Higher values may induce ion fragmentation of low molecular weight samples.

3.2.3 Megaflow Hints

With this high flow rate technique the setup procedure involves making the following adjustments:

- Increase **Desolvation Gas** flow to 500 litres/hour.
- Increase **Desolvation Temp** to 400°C.
- Increase **Source Block Temp** to 150°C.
- Move the probe further away from the sample cone.

Note: When changing from electrospray to megaflow operation it is not necessary to adjust any source voltages.



Attention: The maximum operating temperature for the source heater is 150°C. Do not set **Source Block Temp** higher than 150°C.

Cluster ions are rarely observed with Z-spray. However solvent droplets may form within the source enclosure if the source and desolvation temperatures are too low.

Refer to the previous section on operating parameters for typical desolvation gas flow rates.

Purge gas can be used during megaflow operation to stop the source enclosure from overheating. This is also beneficial when the analyte is susceptible to acetonitrile adducting. Purge gas is enabled by removing the blanking plug from the outlet situated within the source enclosure.

If the sample is contained within a 'dirty matrix' the probe may be moved away from the sample cone to extend time between source cleaning operations. This may incur a small loss in sensitivity.



Caution: It is normal for the source enclosure, the glass tube and parts of the probe mounting flange, to get hot during prolonged megaflo operation. Care should be taken when handling source components during and immediately after operation.

The source enclosure will run cooler if purge gas is used.



Attention: For health and safety reasons always ensure the exhaust line is vented outside the building or to a fume hood.



Attention: Ensure that a plastic bottle is connected in the exhaust line to collect any condensed solvents.

3

3.2.4 Removing the Probe

To remove the probe from the source proceed as follows:

1. On the tune page deselect Operate to put the instrument into standby mode.
2. Switch off the liquid flow and disconnect from the probe.
3. Deselect API Gas and turn off Nitrogen.
4. Disconnect the probe cable from the instrument.
5. Disconnect the nebulizing gas supply from the instrument.

3.3 Sample Analysis and Calibration

3.3.5 General Information

Care should be taken to ensure that samples are fully dissolved in a suitable solvent. Any particulates must be filtered to avoid blockage of the transfer line or the probe's capillary. A centrifuge can often be used to separate solid particles from the sample liquid.

There is usually no benefit in using concentrations greater than 20 pmol/μl for biopolymers or 10 ng/μl for low molecular weight compounds.

Higher concentrations will not usually improve analytical performance. Conversely, for biopolymers, lower concentrations often yield better electrospray results. Higher levels require more frequent source cleaning and risk blocking the transfer capillary.

Optimization for low molecular weight compounds may usually be achieved using a concentration of 1 ng/μl.

Samples with phosphate buffers and high levels of salts should be avoided. Alternatively, at the expense of a small drop in sensitivity, the probe can be pulled away from the sample cone to minimize the deposit of involatile material on the cone.

To gain experience in sample analysis, it is advisable to start with the qualitative analysis of known standards. A good example of a high molecular weight sample is horse heart myoglobin (molecular weight 16951.48) which produces a series of multiply charged ions that can be used to calibrate the m scale from 800-1600 in either positive ion or negative ion mode.

Polyethylene glycol mixtures, for example 300 / 600 / 1000, are low molecular weight samples suitable for calibrating the m scale from approximately 100 to 1200 in positive ion mode. A mixture of sugars covers the same range in negative ion mode.

Alternatively, sodium iodide or caesium iodide can be used for calibration.

Typical ES Positive Ion Samples

- Peptides and proteins.
- Small polar compounds.
- Drugs and their metabolites.
- Environmental contaminants (e.g. pesticides, pollutants).
- Dye compounds.
- Some organometallics.
- Small saccharides.

Typical ES Negative Ion Samples

- Some proteins.
- Some drug metabolites (e.g. glucuronide conjugates).
- Oligonucleotides.
- Some saccharides and polysaccharides.

3.4 Chromatographic Interfacing

Electrospray ionization can be routinely interfaced to reversed phase and normal phase chromatographic separations. Depending on the LC pumping system, chromatography column and setup, there are some basic options:

- Microbore and capillary chromatography separations employing 1mm diameter (and smaller) columns can be interfaced directly to the electrospray probe. Typical flow rates for such columns may be in the region of 3-50 $\mu\text{l}/\text{min}$. It is suggested that a syringe pump is used to deliver these constant low flow rates through a capillary column. Alternatively, accurate pre-column splitting of higher flow rates from reciprocating pumps can be investigated.

In all cases, efficient solvent mixing is necessary for gradient elution separations. This is of paramount importance with regard to low flow rates encountered with capillary columns. HPLC pump manufacturers' recommendations should be heeded.

- 2.1mm diameter reversed phase columns are gaining popularity for many separations previously addressed by 4.6mm columns. Typically flow rates of 200 $\mu\text{l}/\text{min}$ are

used, allowing direct coupling to the electrospray source. The increased sample flow rate requires increased source temperature and drying gas flow rate.

A UV detector may be placed in-line to the probe, provided that the volume of the detector does not significantly reduce the chromatographic resolution. Whenever a UV detector is used, the analog output may be input to MassLynx NT for chromatographic processing.

- The interfacing of 4.6mm columns to the electrospray source can be achieved either by flow splitting or by direct coupling. In both cases an elevated source temperature and drying gas flow rate are required. In general, the best results are obtained by splitting after the column using a zero dead volume tee piece so that 200-300 $\mu\text{l}/\text{min}$ is transferred to the source.



Attention: The maximum operating temperature for the source heater is 150°C. Do not set Source Block Temp higher than 150°C.

Conventional reverse phase and normal phase solvent systems are appropriate for LC-electrospray.

Involatile buffers may be used but prolonged periods of operation are not recommended. When using involatile buffers the probe should be moved as far away from the sample cone as possible. This may reduce sensitivity slightly, but will reduce the rate at which involatile material will be deposited on the sample cone.

Trifluoroacetic acid (TFA) and triethylamine (TEA) may be used up to a level of 0.05%. If solvents of high aqueous content are to be used then tuning conditions should be appropriate for the solvent composition entering the source.

Higher source temperatures (150°C) are also recommended for high aqueous content solvents. Tetrahydrofuran (THF) should not be used with peek tubing.

3.4.1 LC-MS Sensitivity Enhancement

The sensitivity of a LC-MS analysis can be increased or optimized in a number of ways, by alterations to both the LC operation and the MS operation.

In the LC area some examples include the use of high resolution columns and columns with fully end capped pickings. For target compound analysis, techniques such as trace enrichment, coupled column chromatography, or phase system switching can have enormous benefits.

Careful choice of the solvent, and solvent additives or modifiers, may also prove important.

Chapter 4

Nanoflow Electrospray

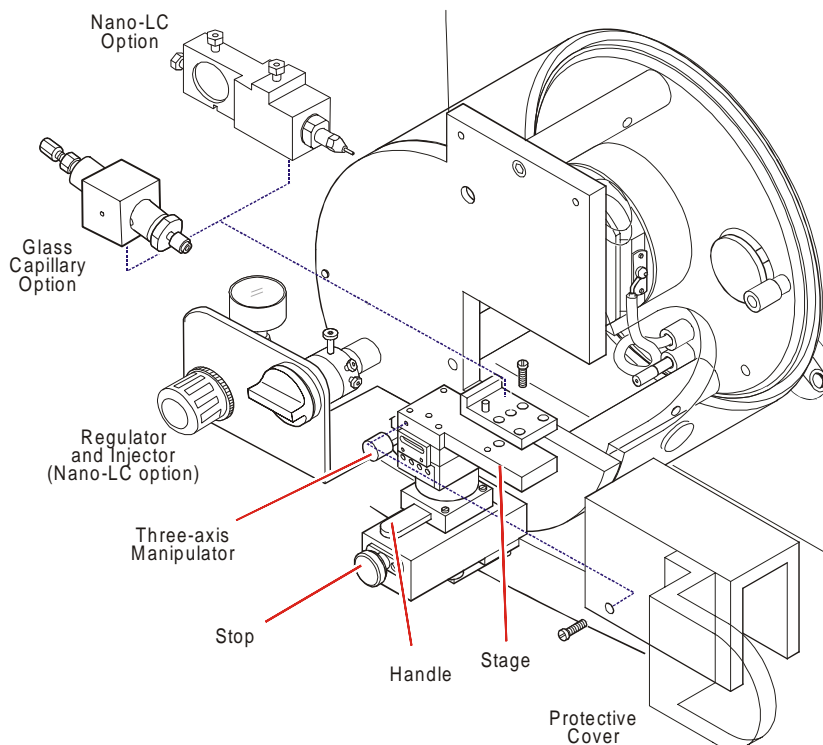


Figure 4-1 The Nanoflow Interface

The optional nanoflow interface allows electrospray ionization to be performed in the flow rate range 5 to 1000 nL/min. There are two options for the spraying capillary, which can be alternately fitted to the interface:

- Borosilicate metal coated glass capillary.

Metal coated glass capillaries allow the lowest flow rates to be obtained although they are used for one sample only and must then be discarded.

- Nano-LC.

This option is suitable for flow injection analyses or for coupling to nano-HPLC, and uses a pump to regulate the flow rate down to 100 nl/min. If a syringe pump is to be used, a gas-tight syringe is necessary to obtain correct flow rates without leakage. A volume of 25µl is recommended.

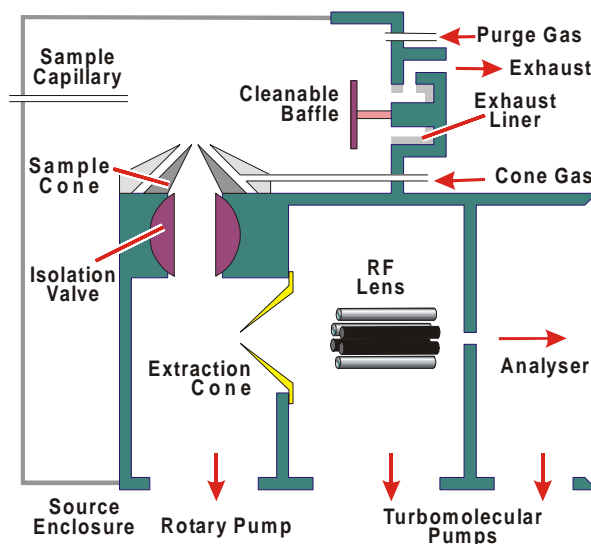


Figure 4-2 Nanoflow Electrospray Scheme

For a given sample concentration, the ion currents observed in nanoflow are comparable to those seen in normal flow rate electrospray. Great sensitivity gains are therefore observed when similar scan parameters are used, due to the great reductions in sample consumption.

The nanoflow end flange consists of a three-axis manipulator, a stage, a protective cover and a stop / handle arrangement for rotation of the manipulator and stage.

The manipulator and stage are rotated by 90 degrees to change option or, in the glass capillary option, to load a new nanovial.



Attention: Failure to use the stop and handle to rotate the stage can result in permanent damage to the three-axis manipulator.

4.1 Installing the Interface

To change from the normal electrospray interface and install the nanoflow interface:

1. If fitted, remove the probe.
2. Remove the moulded cover from around the source.
3. Undo the three thumb screws and withdraw the probe adjustment flange assembly and glass tube.
4. Place the glass tube, end on, on a flat surface and place the probe support flange assembly on top of the glass tube.
5. Remove the PTFE encapsulated source O ring.



Caution: When the source enclosure has been removed the ion block heater is exposed. Ensure that the source block heater has been switched off and has cooled before proceeding. Observe the **Source Block Temp** readback on the Tune Page

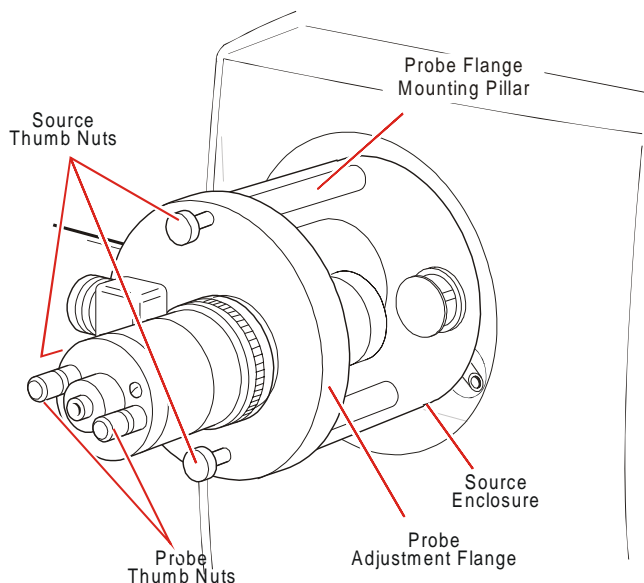


Figure 4-3 Electrospray Source Being Removed

6. Unscrew the three probe flange mounting pillars, using the holes to obtain the necessary leverage.

7. If the cone gas nozzle is not in place, close the sample cone isolation valve. Remove the two screws that secure the sample cone and fit the cone gas nozzle.
8. Replace the two screws.
9. Connect the cone gas outlet to the cone gas nozzle using the PTFE tubing provided (Figure 4-4). Open the sample cone isolation valve.

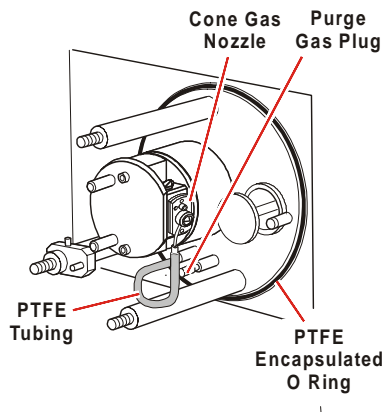


Figure 4-4 Nanoflow Source Without Cover

The cone gas flow rate is set at 30% of the total desolvation gas flow.

10. Ensure that the purge gas is plugged (disabled).
11. Ensure that the cleanable baffle, the exhaust liner and the corona discharge pin blanking plug are fitted.
12. Fit a viton O ring and the three shorter nanoflow pillars.
13. Install the perspex cover and the nanoflow end flange, securing this with socket head screws.

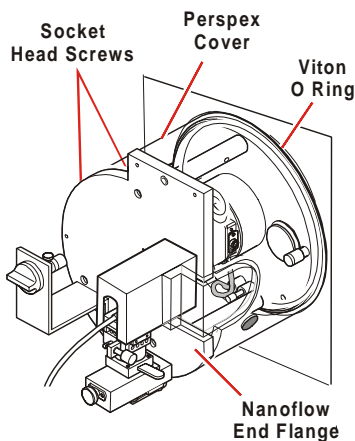


Figure 4-5 Nanoflow Source with Cover and End Flange

Do not attempt to refit the moulded cover.

14. If not already in place, attach the microscope or camera brackets using the screw hole and dowels at the top of the bracket.
15. Insert the flexible light guide into the grommet at the base of the perspex cover.
16. Set the light source to its brightest.
17. Block the **Desolvation Gas** outlet on the instrument's front panel.
18. Close the nebuliser needle valve.

The cone gas is split from the desolvation gas internally.

19. Attach the two cables to the sockets marked Capillary / Corona and Probes on the front panel of the instrument.
20. Set **Source Block Temp** to approximately 80°C.



Attention: The maximum operating temperature for the source heater is 150°C. Do not set **Source Block Temp** higher than 150°C.

4.1.1 Operation of the Camera System

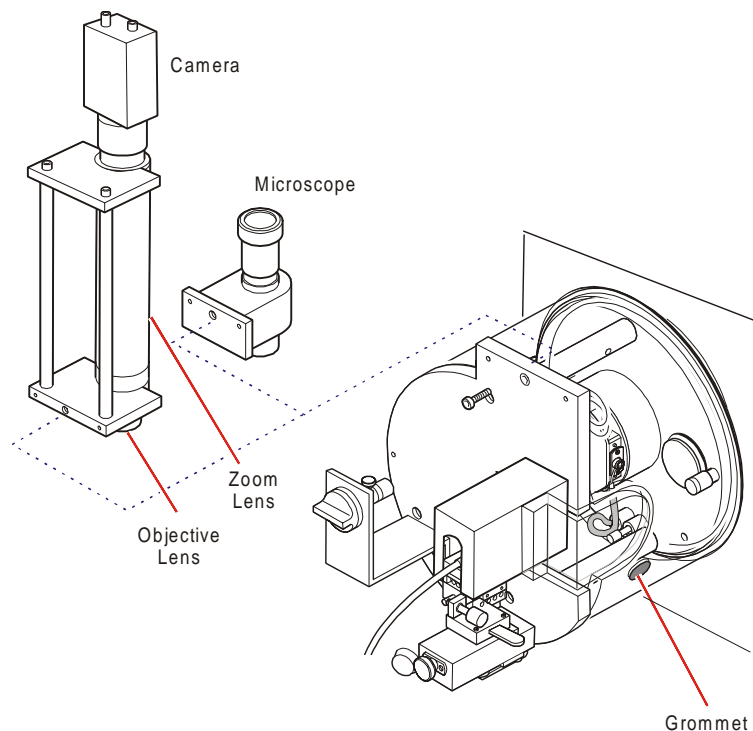


Figure 4-6 Nanoflow Camera and Microscope System

Magnification is controlled by the zoom lens. A fine focus can be achieved by rotating the objective lens.

Using the Microscope

Focusing is adjusted by rotating the top of the microscope.

4.1.2 Glass Capillary Option

Installation



Caution: Do not touch the sharp end of the capillary. As well as the risk of injury by a sliver of glass, the capillary may contain toxic samples.



Caution: The capillaries are extremely fragile and must be handled with great care. Always handle using the square end of the capillary. The needle may become inoperable if the sharp end is touched.

1. With the stage rotated outwards, unscrew the union from the end of the assembly.
2. Carefully remove the capillary from its case by lifting vertically while pressing down on the foam with two fingers.
3. Over the blunt end of the capillary, pass the knurled nut, approximately 5mm of conductive elastomer and finally the union.
4. Tighten the nut (finger tight is sufficient) so that 5mm of glass capillary is protruding from the end of it. This distance is measured from the end of the nut to the shoulder of the glass capillary.
5. Load sample into the capillary using either a fused silica syringe needle or a GELoader tip.
6. Screw the holder back into the assembly - finger tight is sufficient.
7. Ensure that Capillary is set to 0V on the tune page.
8. Rotate the stage back into the interface using the stop and handle.

When using a GELoader tip, break the nanovial in half, by scoring with a fused silica cutter. This enables the GELoader to reach the tip of the nanovial

Nanovial Tip Position

The position of the nanovial is adjusted as shown (Figure 4-7). The tip is in line with the centre of the sample cone, at a distance between two and three times the diameter of the cone orifice, as observed through the microscope.

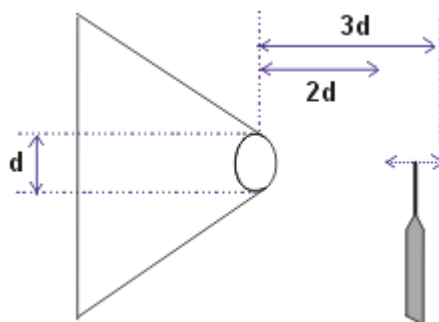


Figure 4-7 Nanovial Tip Position

Operation

Manoeuvre the stage so that the microscope or camera can view the capillary tip.

Using the nanoflow regulator, apply pressure to the back of the tip until a drop of liquid is seen.

On the tune page, select APIGas to turn on Nitrogen.

Select Operate.

Set Capillary between 600 and 1.0kV.

Adjust Desolvation Gas flow to 100 litres/hour using the knob on the front panel of the instrument.

An ion beam should now be visible on the tune page.

Tune the source voltages, adjust the gas flow and adjust the three-axis manipulator for maximum ion current.

The ion current may change dramatically with very slight changes of position but the high resolution of the threads in the manipulator allows very fine tuning.

Restarting the Spray

Should the spray stop, it is possible to restart it by adjusting the three-axis manipulator so that, viewed under magnification, the capillary tip touches the sample cone and a small piece of the glass hair shears off. Set the capillary to zero when doing this.

It may also be necessary to apply some back pressure to the holder to force a drop of liquid from the capillary. Up to 1.4 bar (20 psi) can be applied and, with this pressure, a drop should be visible unless the capillary is blocked.

4.2 Nano-LC Option

4.2.1 Installation

With the sprayer assembly removed from the stage and with reference to Figure 4-8:

1. Cut approximately 25mm of the red stripe peek tubing and, using the plug cap and a Valco nut, set a ferrule to the correct position on the tubing.

At this stage the ferrule is required only to grip the tubing lightly, and should not be too tight.

2. Cut the peek such that 10mm of the peek protrudes from the back of the ferrule.
3. Thread approximately 70mm of the 90 micron o.d. fused silica through the new fitting.
4. Ensure that the fused silica is flush with the peek sleeve.
5. Again using the plug cap, tighten the nut further to ensure that the fused silica is gripped. Some force may be required to do this.
6. Remove the sleeved fused silica from the plug cap and remove the Valco nut.
7. Place an O ring onto the peek tube, using tweezers if necessary.

The O ring is required to seal the region between the ferrule and the end of the thread on the nano-LC chamber.

8. Thread the sleeved fused silica through the nano-LC chamber.
9. Rotate the microvolume union in the body such that the ferrule seat is aligned correctly.
10. Insert the chamber into the nano-LC body and tighten using a pair of spanners.

The capillary can now be checked for flow by connecting the output from a Harvard syringe pump to the other side of the union and setting the flow to 1 $\mu\text{L}/\text{min}$, using a micropipette to measure the flow. It is recommended that a syringe with a volume of no more than 50 millilitres is used.

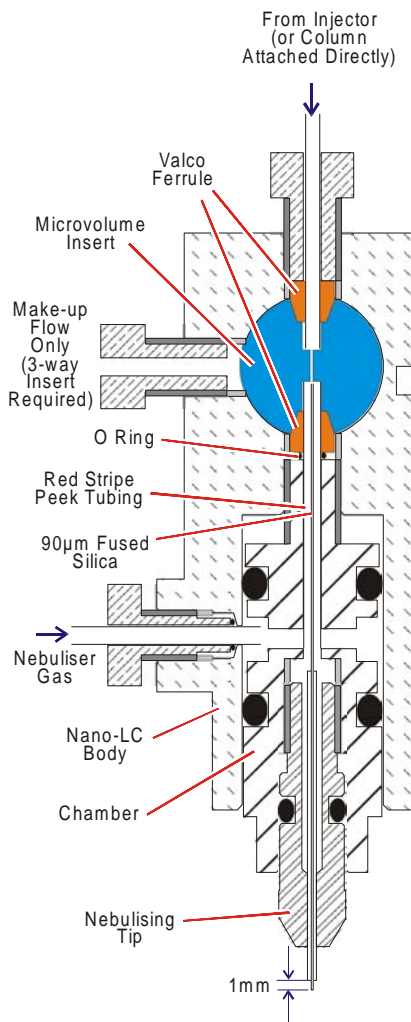


Figure 4-8 Nano-LC Assembly

11. Thread the fused silica through the nebulizing tip and screw in the nano-LC chamber such that it is screwed in approximately half way.
12. Cut the fused silica using a tile cutter and adjust the nebulizing tip further, such that 1mm of fused silica protrudes from the tip.
13. Attach the nebulizing gas tubing to the sprayer using an O ring and the special screw.
14. Attach the sprayer assembly to the stage.

It may be necessary to alter the position of the thumbscrew underneath the baseplate to attach the sprayer correctly.

15. Swing the stage into the interface using the stop and handle.

4.2.2 Operation

For tuning purposes it may be useful to infuse a known sample in 95% water using a Harvard syringe pump.

1. Set the liquid flow to about 200 nl/min.
2. Switch on Gas at the MassLynx tune page.
3. Set the pressure of the gas on the regulator to approximately 0.5 bar (7 psi).
4. Ensure there are no leaks of gas at the sprayer, particularly where the PTFE tubing is connected to it.

Note: *By viewing under magnification, the spray emanating from the capillary may be examined and tuned by altering the nebulizing tip such that a fine spray is observed. Altering the gas slightly may also help in this tuning process.*

5. Swing the stage back out of the source and place the cover over the sprayer ensuring that the tubing coming from the sprayer is threaded correctly through it.
6. Lock the cover in place with two screws.
7. Swing the stage back into the source and alter the translation stage (in / out direction) such that the capillary is approximately 5mm from the cone.
8. Select **Operate** and set **Capillary** to approximately 2.5kV.

An ion beam should now be present.

9. Optimize the ion beam by altering the position of the spray using the controls of the translation stage.
10. The sprayer can now be connected to the HPLC system. The injection valve is plumbed as follows:
 - P from the pump.
 - C to the column (or to the union).
 - S is the sample port, attach a VISF sleeve here.
 - W is a waste port.

Note: *A short tail of fused silica, attached to the entrance port of the union, and the use of low pressure PTFE connectors will remove the need to move the stage. This will prevent*

accidental alteration of the sprayer's position when changing between tuning and HPLC operation.

4.3 Changing Options

To change between the glass capillary and the nano-LC options:

Rotate the stage outwards.:



Attention: *Failure to use the stop and handle to rotate the stage can result in permanent damage to the three-axis manipulator.*

Remove the protective cover and release the captive screw located underneath the stage.

Lift off the holder and replace it with the alternative holder, securing it with the captive screw.

Replace the protective cover, ensuring that either the PTFE back pressure tubing (glass capillary option) or the fused silica transfer line is fed through the slot in the back of the protective cover along with the HV cabling.

Chapter 5

APCI

Atmospheric Pressure Chemical Ionization (APCI) is an easy to use LC-MS interface that produces singly-charged protonated or deprotonated molecules for a broad range of involatile analytes.

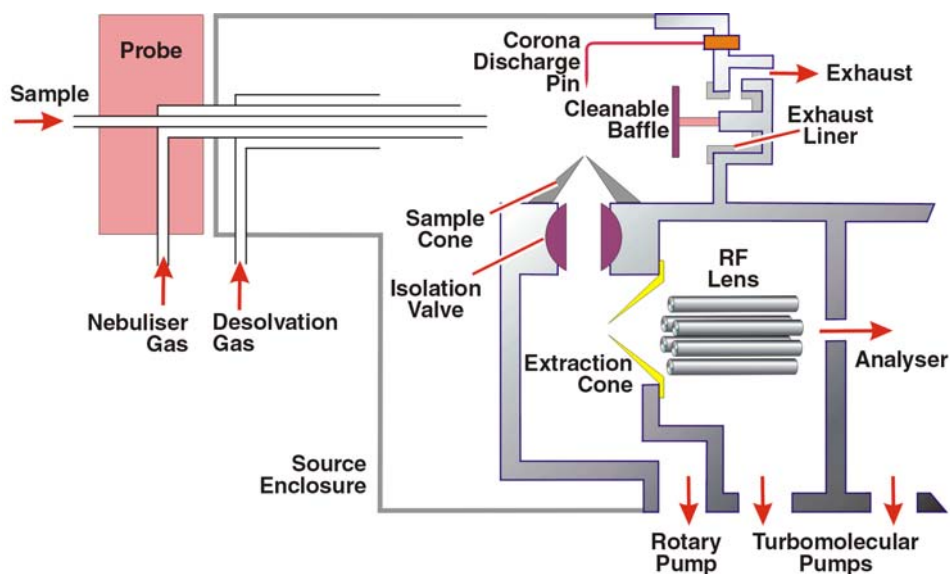


Figure 5-1 Schematic Representation of APCI

The ability to operate with 100% organic or 100% aqueous mobile phases at flow rates up to 2 ml/min makes APCI an ideal technique for standard analytical column (4.6 mm i.d.) normal phase and reverse phase LC-MS.

The APCI interface consists of the standard Z-spray source fitted with a corona discharge pin and a heated nebuliser probe. Mobile phase from the LC column enters the probe where it is pneumatically converted into an aerosol and is rapidly heated and converted to a vapor / gas at the probe tip. Hot gas from the probe passes between the sample cone and the corona discharge pin, which is typically maintained at 2.5kV. Mobile phase molecules rapidly react with ions generated by the corona discharge to produce stable reagents ions.

Analyte molecules introduced into the mobile phase react with the reagent ions at atmospheric pressure and typically become protonated (in positive ion mode) or deprotonated (in the negative ion mode). The sample and reagent ions pass through the sample cone into the ion block prior to being extracted via the extraction cone into the RF lens.

Changeover between electrospray and APcI operation is simply accomplished by changing the probe and installing the corona discharge pin within the source enclosure.

For APcI operation, the desolvation gas is not heated in the desolvation nozzle. However, it is important that desolvation gas is used throughout.

The background spectrum for 50:50 acetonitrile:water is dependent upon the settings of Cone and Extractor. The main reagent ions for typical sample cone and extraction cone voltages of 40V and 10V respectively are 42, 56, 83 and 101.

The transmission of these ions will be dependent on the setting of RF Lens. A lower RF Lens voltage is required for optimum transmission at lower m.

Acetonitrile adducting may be minimized by optimization of the probe position, as described in the chapter entitled Electrospray.

5.1 Preparation

Ensure that the source is assembled as described in Maintenance and Fault Finding, and that the instrument is pumped down and prepared for APcI operation as described in Routine Procedures.

APcI may be operated with or without the cleanable baffle fitted.

Ensure that a supply of nitrogen has been connected to the gas inlet at the rear of the instrument and that the head pressure is between 6 and 7 bar (90-100 psi).

5.1.1 Checking the Probe

1. Ensure that the probe heater is off.
2. Unplug the probe from the instrument's front panel and remove the probe from the source.
3. Connect the PTFE tube to the **Nebuliser Gas** outlet on the front panel.
4. Remove the probe tip assembly by carefully loosening the two grub screws.

5. Disconnect the heater from the probe body by pulling parallel to the axis of the probe.
6. Ensure that 0.5 to 1 mm of fused silica is protruding from the stainless steel nebuliser tube.
7. Connect the LC pump to the probe with a flow of 50:50 acetonitrile:water at 1 ml/min.
8. Check that the liquid jet flows freely from the end of the capillary and that the LC pump back pressure reads 250 to 400 psi.
9. Check that the nitrogen supply pressure is 6 to 7 bar (90 to 100 psi).
10. Select API Gas and turn on Nitrogen.
11. Check that the liquid jet converts to a fine uniform aerosol.
12. Switch off the liquid flow.
13. Deselect **API Gas** and turn off **Nitrogen**.
14. Reconnect the probe tip assembly.
15. Insert the APcI probe into the source and secure it by tightening the two thumb screws.
16. Connect the probe cable to **Probes** on the instrument's front panel.

5.2 Obtaining an Ion Beam

1. Ensure that the corona discharge pin is fitted as described in Section 2.2.2 on page 26 and that the pin is connected using the APcI HV cable.
2. Ensure that the APcI probe is fitted as described above, that the desolvation gas tube is connected to the front panel, and that the cone gas and purge gas outlets are plugged.
3. If necessary, change the ionization mode using the Ion Mode command.
4. Set **Source Block Temp** to 150°C.



Attention: The maximum operating temperature for the source heater is 150°C. Do not set Source Block Temp higher than 150°C

5. Set **APcI Probe Temp** to 20°C with no liquid flow and Nitrogen off.
6. Initially set **Corona** to 2.5kV and Cone to 30V.

When Source Block Temp reaches 150°C:

1. Select API Gas to switch on the nitrogen gas.
2. Using the valves on the front of the instrument, adjust Desolvation Gas to 150 litres/hour and set Nebuliser Gas to its maximum setting.
3. Set one of the peak display boxes to show masses down to at least 100 Da.
4. Select Operate.
5. Set APcI Probe Temp to 350°C.

When APcI Probe Temp reaches 350°C:

1. Start the LC pump at a flow of 1 ml/min.
2. Optimize Corona so that the peaks reach maximum intensity.
3. Optimize the probe position for intensity and stability.

The two screws can be adjusted singly or simultaneously to optimize the beam.

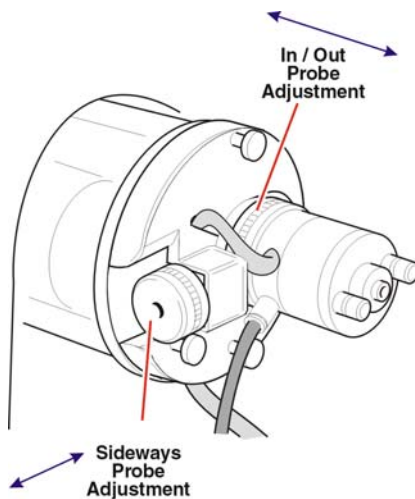


Figure 5-2 Illustration of Probe Adjustment

The position of the probe will affect sensitivity. However, if the sample is contained in a 'biological matrix' or is contained in an involatile solvent the probe should be moved away from the sample cone and towards the corona discharge pin



Caution: It is normal for the source enclosure, the glass tube and parts of the probe adjustment flange to reach temperatures of up to 60°C during prolonged APcI operation. Care should be exercised when handling source components immediately after operation.



Caution: Switch off the liquid flow and allow the probe to cool (<100°C) before removing it from the source.



Attention: Failure to employ a desolvation gas flow during APcI operation may lead to heat damage to the source.

5.3 Hints for Sample Analysis

Tuning

- Start by tuning on the solvent ions.
- It is generally found that the most significant analyte tuning parameter to adjust following tuning on the solvent ions is **Cone**.
- Fine tuning on the analyte of interest can be performed either by large loop injections (100µl) or by constant infusion in the mobile phase typically at analyte concentrations of a few ng/µl.
- 10µl loop injections can be monitored using real time chromatogram updates.

Mobile Phase

- The choice of mobile phase is an important compound specific factor in APcI. For example, steroids prefer methanol:water mixtures as opposed to acetonitrile:water.
- Analyte sensitivity is also dependent on mobile phase composition, which can be varied from 100% aqueous to 100% organic for any particular mixture.

Probe Temperature

This can be a critical factor for some analytes.

- Involatile samples (for example steroids) generally require high probe temperatures ($>400^{\circ}\text{C}$).
- Volatile samples (for example pesticides) can be analyzed with low probe temperatures ($<400^{\circ}\text{C}$).
- In some cases, too high a probe temperature can lead to thermal degradation of labile samples.

Desolvation Gas

Although a **Desolvation Gas** flow of approximately 150 litres/hour is typical for most samples, this flow rate should be tuned for maximum sensitivity while ensuring that the flow rate is not decreased below 100 litres/hour.

5.4 Removing the Probe

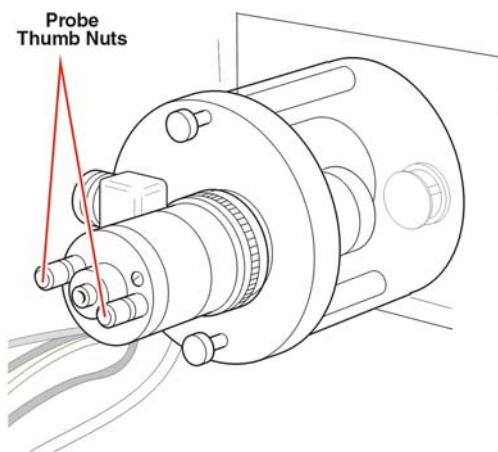


Figure 5-3 APcI Probe showing twin Thumb Nuts used for Removal

After a session of APcI operation:

Turn off the LC flow.

Set **APcI Probe Temp** to 20°C .

Deselect **Operate** to put the instrument in standby mode.

When the probe temperature falls below 100°C :

Deselect **API Gas** and turn off **Nitrogen**.

Undo the two thumb nuts and remove the probe from the source (Figure 5-3).



Caution: Take care when removing the APcl probe. There is a risk of burns to the operator.



Attention: Removal of the APcl probe when hot will shorten the life of the probe heater.

If the instrument is not to be used for a long period of time the source temperature should be reduced to 60°C.

Chapter 6

Tuning

Before sample data are acquired, the instrument should be tuned and, for highest mass accuracy, calibrated using a suitable reference compound (see Chapter 7).

Consult the relevant sections of this manual for information concerning source tuning procedures in the chosen mode of operation. See:

- Chapter 3, Electrospray
- Chapter 4, Nanoflow Electrospray
- Chapter 5, APcI

For information on optimizing peak shape, intensity and resolution and adjusting the tuning parameters see the following sections:

- Section 6.2.4, Source Tuning Page
- Section 6.2.5, Quadrupole Tuning Page (Analyzer)
- Section 6.2.6, Time Of Flight Page
- Section 6.3, Obtaining an Ion Beam

Note: Care should be taken to optimize the value of the collision energy.

6.1 Tuning

The preceding sections have outlined the software controls and connections to establish the physical conditions required as a prerequisite to mass spectrometer operation.

Within the source and analyzer enclosures, electric fields controlled via MassLynx are applied to the components to manage the ion beam generated, according to the sample to be analyzed, the ionization mode and the type of information required. The ion optical system elements are indicated in Figure 6-1 on page 76.

Considering the variable nature of the beam with different samples, the instrument should be tuned for signal strength and calibrated for accurate m using suitable reference compounds, prior to the acquisition of sample data.

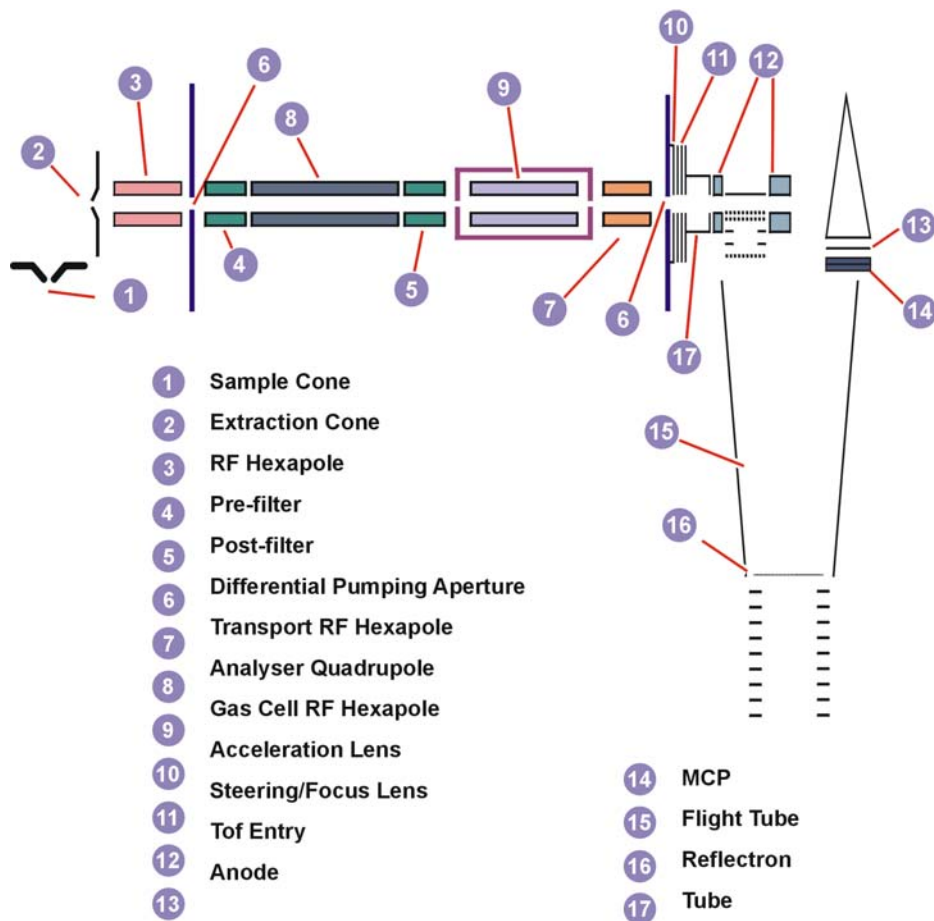


Figure 6-1 The Ion Optical System of the Q-ToF *micro*

Tuning parameters have been grouped in MassLynx into 3 pages as in Section 6.3 on page 89. Full details of source tuning procedures for electrospray, APcI and nanoflow electrospray are given in the relevant chapter of this manual.

6.2 Tune Page

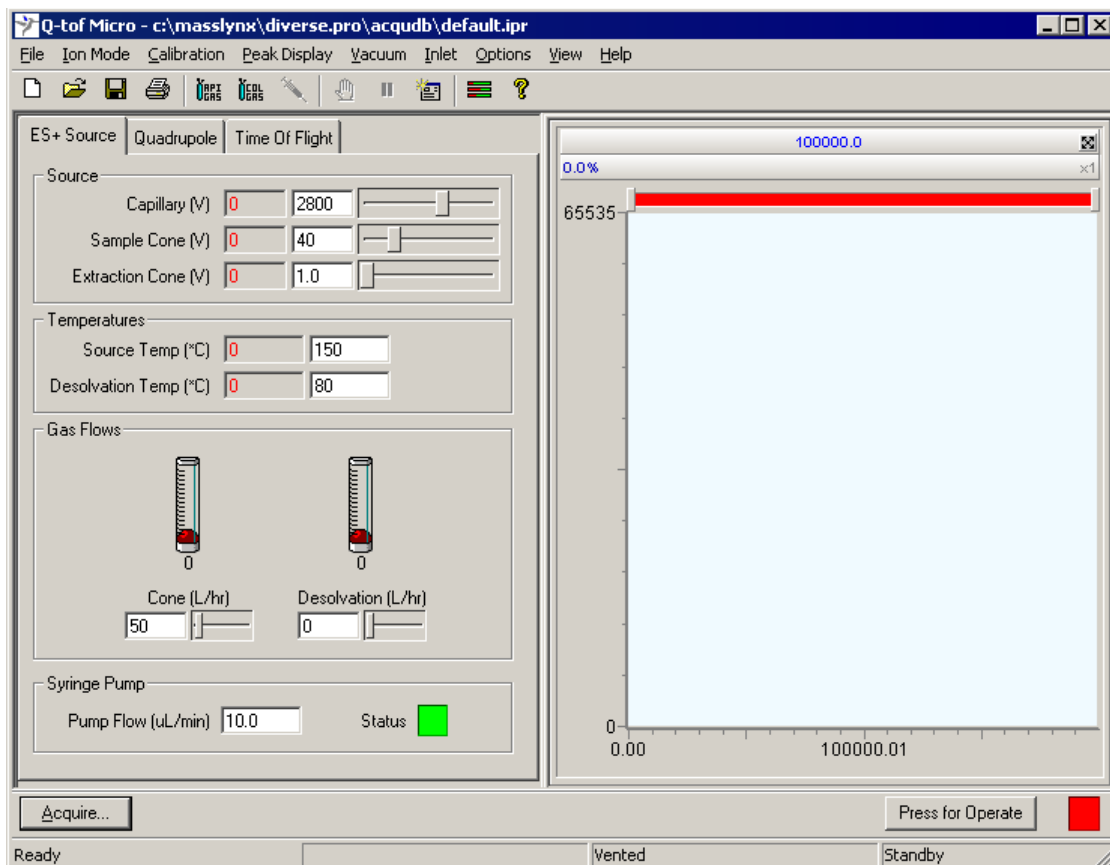


Figure 6-2 The QTof *micro* Tune Page

The Tune Page Consists of a peak display window and 3 pages of parameters viewable by selecting the relevant tab. These are:

- The Source Page (see Section 6.2.4)
- The Quadrupole Page (see Section 6.2.5)
- The time of Flight Page (see Section 6.2.6)


6.2.1 Tune Page Basics

To Display the Tune Page:

Select MS Tune from the Masslynx Shortcut bar (Figure 6-2).

Printing Tune Information

To print a report, containing a copy of the tune peak information displayed on the screen along with a record of each parameter setting:

Press , or select Print from the Tune Page File menu.

This report is not configurable by the user.

Experimental Record


Tuning parameters are stored with every data file as part of the experimental record. The tuning parameters for a particular data file can be viewed or printed from the data browser, see the MassLynx User Guide, Selecting and Viewing Data, for more information.

6.2.2 Saving and Restoring Parameter Settings

Whole sets of instrument tuning parameters can be saved to disk as a named file and then recalled at a future date.

A tune parameter file contains the latest settings for the source controls for all supported ionization modes not just the ionization mode currently selected. Tune parameter files also contain settings for the analyzer, inlet set points and peak display.

To save the current tune parameters with the existing file name:


1. Press , or choose **Save** from the tune page **File** menu.
2. Press **Save**.

To save the current tune parameters with a new file name:

1. Select **Save As** from the tune page **File** menu.
2. Enter a new file name or select an existing file from the list displayed.
3. Press **Save**.

If the selected file already exists on disk a warning is displayed. Press Yes to overwrite the existing information or No to enter a different file name.

To restore a saved set of parameters:

1. Press , or choose **Open** from the tune page **File** menu.
2. Select the required tuning parameter file, either by typing its name or by selecting from the list displayed.
3. Press **Open**.

6.2.3 Changing the Peak Display

To change the display using the mouse:



Figure 6-3 The Customize Display Menu

Click in the peak display area with the right mouse button to display the pop up menu (Figure 6-3).

The display area for each peak can be individually changed, e.g. the peak color for peak 1 can be red and for peak 2 green etc.

Customize Plot Appearance

To change the color of the background and traces and to change the number of traces displayed:

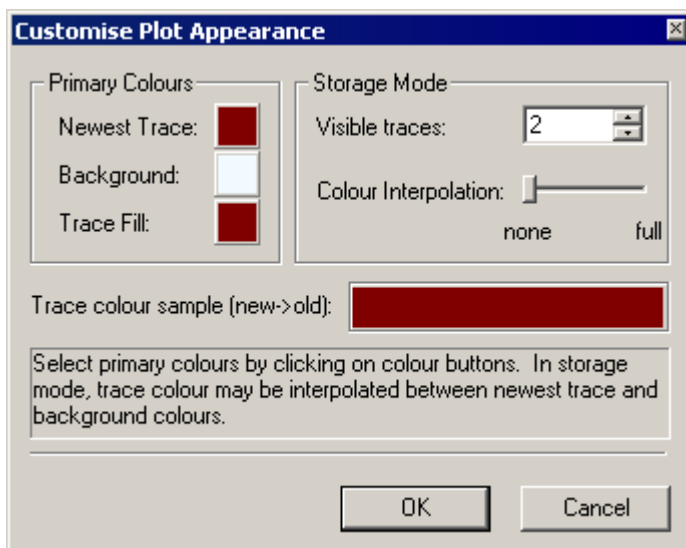


Figure 6-4 The Customise Plot Appearance Dialog

Select Customize, Plot Appearance.

The Customize Plot Appearance dialog is displayed (Figure 6-4).

To change the colours on the display:

Press **Newest Trace**, **Background** or **Trace Fill** and select a new color from the dialog displayed.

To change the number of traces:

Use to change the number, or enter a new value in the **Visible Traces** box, within the range 2 to 20.

If more than one trace is displayed then the older traces can be displayed in a different shade to the new ones:

Drag the **Color Interpolation** slider toward the full position. The color of the old traces is shown in the **Trace color sample (new->old)** field.

Trace

From the pop-up menu do one of the following:

- Select the **Trace, Outline** option to display the peak outline only.
- Select the **Trace, Fill** option to fill the trace with the trace fill color.
- Select the **Trace, Min/Max** option to show the minimum and maximum data points only.

The option selected has a tick next to it.

Intensity

Select either **Intensity**, **Relative Intensity** or **Intensity, Absolute Intensity** as required.

Select Intensity, **Normalize Data** to display normalized data.

The options selected each have a tick next to them.

Grid

The options allow vertical and horizontal grid lines to be independently displayed or hidden.

Selected Options have ticks next to them. Selecting an option a second time deselects the option.

6.2.4 Source Tuning Page

The main parameters that are changed during tuning can be found on this page. The two that have the most effect are the **Capillary Voltage** and **Sample Cone Voltage**.

ES+ Source | Quadrupole | Time Of Flight

Source

Capillary (V) 2958 3000

Sample Cone (V) 153 30

Extraction Cone (V) 121.0 2.0

Temperatures

Desolvation Temp (°C) 121 120

Source Temp (°C) 83 80

Gas Flows

Cone (L/hr) 20

Desolvation (L/hr) 347 360

Syringe Pump

Pump Flow (uL/min) 0.0 5.0 us

Figure 6-5 The Source Tuning Page

The positive ion electrospray (ESP+) Source Tuning page is shown. Suggested tuning parameters are as follows:

Capillary	This sets the absolute voltage on the electrospray probe or APCI corona needle, and is typically adjusted to 3000V . The value shown in the display is for illustration only.
Sample Cone	This sets the voltage on the sampling cone relative to the extraction lens. It is dependent on compound and charge state. For multiply charged species this is set to 30 - 50eV , and higher for singly charged species. In general, higher cone voltages are needed for larger mass ions.
Extraction Cone	This sets the voltage on the extraction lens, and is normally set from 0 - 2V .
Desolvation Temp	This sets the temperature of the desolvation gas heater. It is usually set to 150°C , and increased for higher solvent flow rates.
Source Temp	The source block temperature is usually set to 80°C , but is increased for higher solvent flow rates.
Gas Flows	Gas Flow will normally be set at 20 l/h for cone and 380 - 400 l/h for Desolvation
Syringe Pump	Enter the desired flow rate for the Analyte.
Lock Spray	Choose from Reference or Analyte . Note: This box only appears if LockSpray is selected on the MassLynx Options dialog

6.2.5 Quadrupole Tuning Page (Analyzer)

The screenshot displays the Quadrupole Tuning Page (Analyzer) with three tabs: ES+ Source, Quadrupole (selected), and Time Of Flight. The page is organized into three main sections: Quad, Gas Cell, and Automatic MS Profiling.

Quad Section:

- Ion Energy (V): 1.6
- Collision Energy (V): 10
- Low Mass Res: 5.0
- High Mass Res: 5.0
- RF 1 DC Offset (V): 119.3 (red), 1.5
- Pre/Post Filter (V): 113.7 (red), 4.9
- RF 2 DC Offset (V): 96.9 (red), 11.2
- Pump Aperture 2 (V): 102.9 (red), 5.0

Gas Cell Section:

- Plate 1 (V): 110.7 (red), 3.0
- Entrance (V): 108.7 (red), 109.0
- Can (V): 0, 0
- Plate 2 (V): 105.7 (red), -3.0

Automatic MS Profiling Section:

- Auto MS Profiling: ☒
- Quad Mass (AMU): 0
- Gas Cell RF (V): 97

Figure 6-6 The Quadrupole Tuning Page

Note: The Parameters on this page do not routinely need changing.

Ion Energy

This sets the quadrupole offset DC with respect to the collision energy of the ions travelling between the quadrupoles. This should be set between 1 - 1.8V.

Collision Energy	This sets the collision energy of the ions when they reach the collision cell. The cell itself is grounded but the collision energy voltage is simultaneously applied to all the optical elements preceding it. i.e cone, extraction lens, source hexapole, differential pumping aperture, and quadrupole. In Q-ToF micro gas is always introduced into the cell, affording collisional cooling and consequently higher resolution. The collision energy is set to 10eV for MS mode to maximize ion transmission but produce little or no fragmentation. When in the MS/MS mode the collision energy is adjusted to give the best fragmentation pattern for the selected parent
Low Mass Res and High Mass Res	These set the resolving DC on the quadrupoles. The two sliders are set to give constant resolution across the mass range. When the quadrupoles have been set up, settings of 15 on both sliders should give unit resolution at 20% peak height.
RF1 DC Offset	Normal value is between 0 and 1V. The value shown in the display is for illustration only (item 3 in Figure 6-2).
Pre/Post Filter	Normal value is between 5 and 15V (items 4 and 5 in Figure 6-2).
RF2 DC Offset	Normal value is between 0 and 4V (item 7 in Figure 6-2).
Pump Aperture	Normal value is between 2 and 10V (item 6 in Figure 6-2).
Plate 1	This is the voltage applied to the entrance plate of the gas cell. Normal value is 2V. The value shown in the display is for illustration only.
Entrance	This sets the voltage on the pusher entrance and exit, and defines the axial speed of the ions through the TOF. The theoretical value is 105 eV when the ion beam should be central to the detector. It should be optimized when looking at a TOF beam, not on the first detector. It should be possible to correctly tune the instrument using this value. The value shown in the display is for illustration only.
Can	This is the voltage applied to the hexapole gas cell surrounding shield. Normal value is 0.
Plate 2	This is the voltage applied to the exit plate of the gas cell. Normal value is -3V.
Auto MS Profiling	If the box is checked then automatic MS/Hex profile will be used. If the box is not checked then Quad Mass and Gas Cell must be set manually.

Quad Mass This is the mass in Da at which the quad mass is set in RF only mode

Gas Cell RF (V) This is the peak to peak amplitude of the applied RF voltage.

6.2.6 Time Of Flight Page

The screenshot shows the 'Time Of Flight' tab selected. The interface is organized into three main sections:

- Transfer Lens:**
 - Acceleration Lens (V): -91 (red), 200
 - Focus (V): 218 (red), 0
 - Steering (V): 5.5 (red), 2.6
 - Tube Lens (V): -91 (red), 90
- Tof Flight Tube:**
 - Grid 2 Offset (V): 4.25
 - Flight Tube (V): -5644 (red), 5630
 - Reflectron (V): 1721 (red), 1730
 - Pusher Offset (V): 0.80
 - Pusher (V): 812.7 (red), 810.0
 - Puller (V): 718.9 (red), 710.0
 - Manual Pusher: ☐
 - Cycle Time (us): 40
 - Width (us): 9
- MCP:**
 - MCP Detector (V): 2675 (red), 2700

Figure 6-7 The Time Of Flight Page

Note: The Parameters on this page do not routinely need changing.

Acceleration Lens	Usually set to maximum voltage (200V) in all modes of operation. The factory setting should be noted, before attempting to improve resolution by changing this voltage.
Focus	This adds an equal voltage to both the top and bottom steering/focus lens halfplates. For maximum TOF resolution it is set to zero.
Steering	This adjusts the voltage difference between the top and bottom half plates of the steering/focus lens. It acts as a y-deflector, directing the beam into the pusher. A setting close to zero should produce an optimum beam. Better sensitivity can be achieved, however, with a voltage (positive or negative) close to zero, e.g. +1.0, -0.5 etc., and not more than $\pm 2V$.
Tube Lens	This is set to its optimum value in the factory and by the engineer at installation. It helps shape the ion beam on entry into the pusher and so has a large effect on resolution and peak shape.
Grid 2 Offset	<p>This determines the voltage difference between the two plates of the first acceleration region of the ToF (Pusher and Grid3). It is used as a fine tuning control to optimize resolution.</p> <p>The optimum value for this element may change with pusher frequency. The pusher frequency is dependent on mass range when in automatic pusher mode and takes one of six discrete values. See pusher rates on Page 41.</p>
Flight Tube	<p>This sets the flight tube voltage. It is always set to 5630V.</p> <p>Adjusting this will change the peak position and resolution.</p>
Reflectron	This is set expressed as a percentage of the flight tube voltage and its value is 1780V. Adjusting this will change the peak position and resolution.
Pusher Offset	This is the DC value applied to the puller plate to collimate the beam. Normal value is -1V. The value shown in the display is for illustration only.
Pusher	This sets the amplitude of the pusher pulse and is normally set to 830V.
Puller	This sets the amplitude of the puller pulse and is normally set to 645V.

MCP	This sets the voltage on the TOF detector. MCP must be conditioned before applying high voltage (see Section 2.1.3 on page 23). The TOF analyzer is usually operated with this detector set at 2800V. The value shown in the display is for illustration only. <i>It is recommended that a record of these values is kept for future reference.</i>
Manual Pusher	If the Manual Pusher box is selected from the Other sub-menu then the repetition frequency of the pusher pulse is determined by the Cycle Time entered. See below.
Width	This is the width of the pusher pulse in μs , the default value is 9 μs but for low mass ($< 80\text{Da}$) acquisitions a lower value may be needed.

Manual Pusher

Cycle Time may be set between 30 μsec and 255 μsec . Setting Cycle Time at less than 30 μsec defaults to 30 μsec . Entering 0 μsec switches the pusher off.

Increasing the flight time reduces the duty cycle (sampling efficiency) of the TOF analyzer resulting in decreased sensitivity.

For the best mass accuracy the instrument should be re-calibrated if the flight time is changed.

If Manual Pusher is not selected then the repetition frequency of the pusher pulse is determined automatically according to the highest m requested in the acquisition range, as shown in the table below.

Maximum Flight Time	Highest
33 μsec	1000
47 μsec	1001 - 2000
66 μsec	2001 - 4000
94 μsec	4001 - 8000
132 μsec	8001 - 16000
188 μsec	16001 - 32000

6.2.7 Other Tune Page Settings

TDC Settings

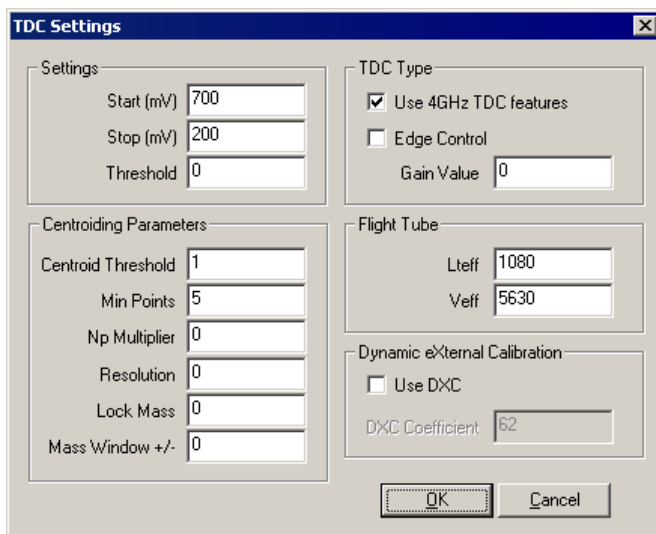


Figure 6-8 The TDC Settings Dialog

To access the TDC (time to digital converter) settings select **Options >TDC settings**.

Start	This is the size of the trigger signal that is necessary to trigger the TDC (start the clock). The start signal is derived from the pusher voltage itself, and a typical value is 500mV. This voltage may be different in negative ion mode.
Stop	This is the size of pulse needed to register as being an ion, so stopping the clock. It is usually set at 150mV, a value high enough to prevent electronic noise being detected as ions.
Threshold	This parameter should normally be set to zero. Setting to 1 will cause all peaks in the spectrum with one count to be thresholded out.
Centroid Threshold	Any peak above this value in height (counts) will be classed as ions and a centroid will be produced
Min Points	Any peak above this value in width (bins) will be classed as ions and a centroid will be produced.
NP Multiplier	Used in deadtime correction calculations. A normal value is 0.7.

Resolution	Used in deadtime correction calculations. A normal value is 6000, but will vary according to instrument resolution.
Lock Mass	The mass to be entered in accurate /exact mass corrections to ensure data drift is minimized.
Mass Window	The window searched around the set lock mass in order to establish a suitable peak to be the lock mass
Veff	Effective flight tube voltage, used for setting up the non-calibrated mass scale. A normal value is 5630.
Temp Correction	Enables the DXC corrections to be applied, if installed.
Lteff	This is used to make the TOF mass measurement nominally correct without a calibration. The value is usually around 1080. See Chapter 7 for further details.

6.3 Obtaining an Ion Beam

This section refers exclusively to the positive ion mode of operation using the ESI source. This is often the simplest mode of operation for getting and checking the ion beam.

Note: *This method should also be followed through when it is suspected there may be a problem with the instrument. It is likely that, should the user ever need to contact Micromass' Technical Support specialists, they will ask that these steps have been performed.*

It is standard practice to use a sample with which the user is familiar. This sample should be of reasonable concentration and have at least one readily recognized component or pattern of masses.

Commonly used compounds are listed below:

Compound	Mwt
Leucine Enkephalin	556Da
PEG 200/400/600 (mix)	100-1000Da
Sodium Iodide	100-2500Da

A normal concentration is 2ng/μl. The general method outlined here can be used in either mode - positive or negative.

The sample should be infused at a rate of **5 μ l/min** via a syringe pump (on the instrument or otherwise) or using a dedicated LC system. This ensures a uniform, continuous spray that will last for a relatively long period.

6.3.1 Standard Tune Parameters

1. Enter the Standard Tuning Parameters (or load from a file) as described in Section 6.2 and Section 6.3.

Ensure that the Peak Display window is open.

2. Switch the instrument into **OPERATE**.

Press the button on the tune page marked ‘Press For Operate’ to get the following (the indicator goes to green and the button changes to ‘Press For Standby’)

3. Ensure that the **API Gas** , **Col Gas** , and **Syringe** , buttons are selected and operating.

Note: *It may take a while for the sample to appear in the peak display window as any dead volume will need to be cleared first.*

4. Assuming that the probe is spraying the **Peak Display** window should be showing some response.

Note: *Adjustment of the Capillary and Sample Cone Voltage will vary the peak height and shape*

6.3.2 Checking the Resolution

While there may be a beam showing in the peak display, a single peak needs to be “Tuned” to give a resolution of at least **5000** (FWHM), for the QToF *micro* to perform to its specification.

The resolution is calculated by the following equation:

$$Resolution = (Mass) \div (PeakWidthHalfHeight)$$

In the following example Leucine Enkephalin (2ng/ml) has been infused at 5 μ l/min.

1. Press the Tune Page **Acquire** button. This invokes the Tune page Acquisition dialog (Figure 6-9)

Enter a **Data File Name** together with any necessary comments.

Ensure **TOF MS** is selected for the **Function** and that the **Data Format** is set to **Continuum**.

Set the mass range from **100 Da** to **1000 Da** and a Duration of 1 minutes.

The scan time can be any reasonable value, but a setting of 1s with an inter scan time of 0.1 s gives one second scans for each individual spectrum saved to disk.

2. Click **Start**.

Tune Page Acquisition

Data File Name: RESOLUTION

Text: leu Enk

Function: TOF MS

Data Format: Continuum

Masses

Set Mass (m/z): 50

Start Mass (m/z): 100

End Mass (m/z): 1000

Times

Run Duration (mins): 1

Scan Time (s): 1.00

Inter Scan Time (s): 0.1

Start Close Origin...

Figure 6-9 The Tune Page Acquisition Dialog

3. Open **Chromatogram** from the MassLynx Sample List menu and click on Real Time Update.
4. **Combine** data from at least 30 scans. On completion the **Spectrum** window will open

Note: See the MassLynx User's Guide for details on how to use Chromatogram, Spectrum and Combine.

5. In the Spectrum window Zoom in on the main peak at 556 Da.

Note: Ensure that a Grid is shown on the Spectrum. Select **Display > View** to display the **Spectrum Display** dialog

- With the aid of the grid right click and drag the mouse across the width of the peak at half height. The width is show in the bottom left corner of the Spectrum window (Figure 6-10)

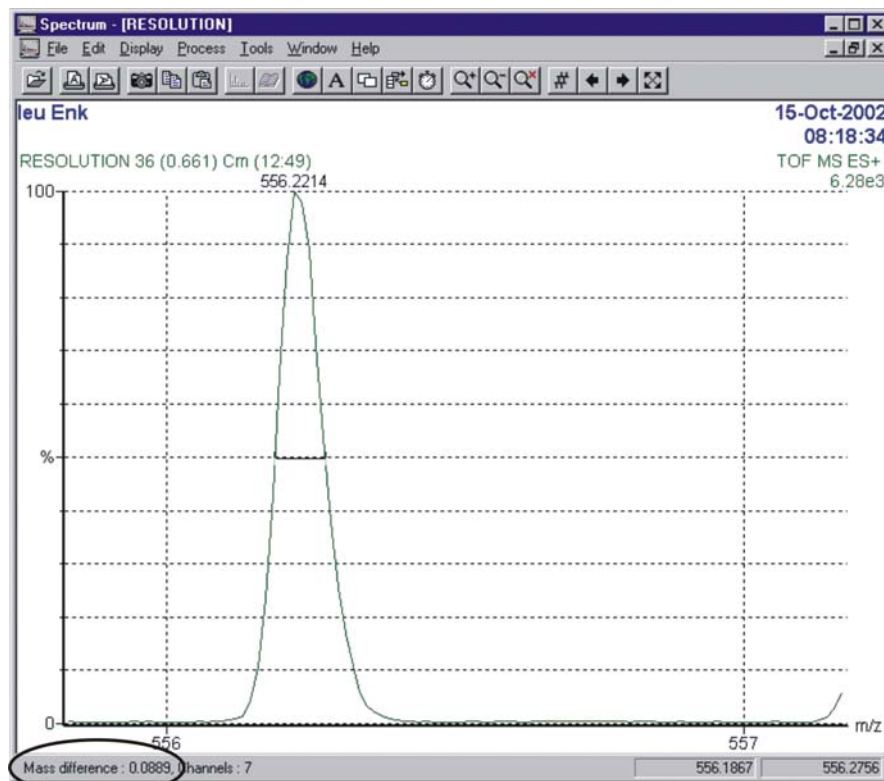


Figure 6-10 Combined Spectrum of Leucine Enkephalin Showing Peak Width at Half Height

- The resolution can now be calculated.

$$\text{Resolution} = 556.2214 \div 0.0899 = 6187$$

In this case it is quite clearly adequate and above 5000.

If the resolution is not above 5000 the process will have to be repeated with different tuning parameters. With experience it will possible estimate the resolution from the Peak Display on the Tune Page.

Note: If a resolution of 5000 can not be achieved it is likely that there is a problem with the instrument and Micromass Technical Support should be contacted.

Chapter 7

Calibration

There are several steps to carrying out a Calibration for Accurate Mass on the Micromass Q-ToF Micro MS:

1. Calculation of Nominal Mass Accuracy (Section 7.2.1)
2. Acquire data for Calibration. (Section 7.2.20)
3. Check for Deadtime distortion (Section 7.2.3)
4. Calibration (Section 7.2.4)
5. Deadtime Correction (Section 7.2.6)

You will need to refer to the following chapters of the MassLynx User's Guide.

- Chapter 6 - Chromatogram
- Chapter 7 - Spectrum
- Chapter 8 - Strip and Combine

Note: *To use the Calibration file for Accurate Mass refer to Chapter 15 "Accurate Mass Measure" of the MassLynx User's Guide.*

7

7.1 Introduction

Extremely accurate mass measurements can be performed with the Q-ToF micro, due to the elevated resolution and inherent stability of the calibration law of orthogonal TOF instruments.

The basic time of flight calibration from mass (m) to time (t) is of the form:

$$\sqrt{m/z} = Q + Pt$$

where:

The term P represents the resultant gain from the instrument geometry (pathlengths and voltages).

Q is an offset, arising from propagation delays through the electronics (detector rise time and delays of trigger signals through cables).

If a data file is acquired from the instrument with no calibration applied, then it is assumed that the offset is zero and the gain P is calculated from the instrument geometry.

7.1.1 Nominal Mass Accuracy - Lteff

Lteff is the effective length of the flight tube, it is important that this should be set up to give at least nominal mass accuracy.

Nominal mass measurement is achieved on the Q-ToF micro by adjustment of the **Lteff** factor, a term which quantifies the difference between the indicated and actual mass.

A TOF spectrum of a standard compound is acquired with **Lteff** set to its default value of 1078.

A new value of **Lteff** can be calculated from the following relation:

$$Lteff = 1078 \sqrt{m_{ind} \div m_{act}}$$

where:

m_{ind} = indicated m/z

m_{act} = actual m/z

The new value is entered on the Tune Page under **Options > TDC Parameters** (Figure 7-1).

Figure 7-1 TDC Settings Dialog

All Subsequent mass measurements will be nominally correct. This value should only ever need to be set once.

7.1.2 Calibration

With no calibration applied, the spectral data in MassLynx is merely a set of mass intensity pairs $\{M_n, I_n\}$ based upon instrument geometry.

The inherent relationship between mass and time shown above makes it prudent to generate higher order calibration coefficients that are applied to the square root of the nominal masses $\{M_n\}$:

$$\sqrt{M_c} = A + B\sqrt{M_n + CM_n + DM_n^{3/2}}$$

where:

The terms A, B, C, D.... are calculated by fitting a polynomial to the acquired mass spectral data.

M_c is the calibrated displayed mass.

- If a polynomial of order 1 is requested, the values for A & B are calculated, and the higher terms are set to zero.

- For a polynomial of order 5 (the highest supported in MassLynx) there will be six terms generated.
- When calibrating over a large mass range (>500 Da) it is advisable to use a higher order polynomial, as the deviations from the straight line fit become more appreciable.
- Once a calibration has been generated from a reference compound such as PEG it should be used as an ‘instrument calibration’ to be applied to all subsequently acquired data.

The procedure for this is described below.

7.1.3 Lock Mass

Temperature variations in the environment and in the instrument power supplies can cause drifts in measurements of a few hundred parts per million (ppm) over the course of a day. For accurate mass work, the instrument should be kept in OPERATE at all times to enable stabilization of the power supplies.

Users can compensate for instrument drift by applying a single point lock mass correction that recalculates the term B in the above equation. This can be done in two ways:

1. Using a Micromass LockSpray / NanoLockSpray Interface
2. Teeing in, post column, a Lock Mass reference compound.

Note: *There are inherent problems with the second approach which can lead to poor exact mass capability. Consequently we recommend that the first option is used when ever possible.*

7.1.4 Deadtime Correction

The data acquisition system for the instrument is a time to digital converter (TDC). This is an ion counting system which generates a mass spectrum by storing the arrival times of ions in a histogram memory.

After the arrival and registration of an ion by the TDC there is a minimum time interval before a subsequent ion arrival can be registered. This is called the ‘dead time’ of the TDC and is of the order of 5 nanoseconds.

At high ion currents some of the ions generated are not registered, leading to a shift to lower mass centroids, with consequently lower measured areas on reported peaks.

However, the MassLynx software incorporates a correction facility which allows for accurate mass measurements to be achieved over a large range of ion currents, and the use of dead time correction is described below.

7.2 Generation of an Instrument Calibration

Make up an analyte solution consisting of:

PEG 300 (10 nanolitres/millilitre),

PEG 600 (10 nanolitres/millilitre),

in a stock solution of 2mM ammonium acetate in 50/50 acetonitrile/water.

Dilute this by a factor of ten (with the stock solution) and introduce this to the instrument using a syringe pump operating at 5 microlitres per minute.

Note: The Quad ramp is predetermined for every mass range and can not be set by the user.

7.2.1 Calculation of Lteff

The screenshot shows the 'Tune Page Acquisition' dialog box. It has a title bar with a close button. The dialog contains several input fields and dropdown menus. The 'Data File Name' field is set to 'LEUCINE_ENKEPHALIN'. The 'Text' field is set to 'Lteff Nominal Mass Accuracy'. The 'Function' dropdown is set to 'TOF MS'. The 'Data Format' dropdown is set to 'Continuum'. There are two sections: 'Masses' and 'Times'. The 'Masses' section has three input fields: 'Set Mass (m/z)' set to 50, 'Start Mass (m/z)' set to 100, and 'End Mass (m/z)' set to 1000. The 'Times' section has three input fields: 'Run Duration (mins)' set to 2, 'Scan Time (s)' set to 0.52, and 'Inter Scan Time (s)' set to 0.1. At the bottom, there are three buttons: 'Start', 'Origin...', and 'Close'.

Figure 7-2 The Tune Page Acquisition Dialog

To calculate L_{teff} it is recommended that a solution such as Leucine Enkephalin is used that will produce a single peak.

1. Make up a solution of Leucine Enkephalin at 2 ng/ml in 50:50 methanol:water.
2. Introduce this into the instrument with the syringe pump at 5 $\mu\text{l}/\text{min}$.
3. Click **Aquire**, this will invoke the **Tune Page Acquisition** dialog (Figure 7-2).
4. Acquire for about 2 minutes.
5. **Combine** the data and then **Center** the Spectrum.

Note: For further information on Combining data see Chapter 8 of the MassLynx User's Guide and for Centering Spectral data, Chapter 6.

6. Note the indicated mass. In this case 561.2170.

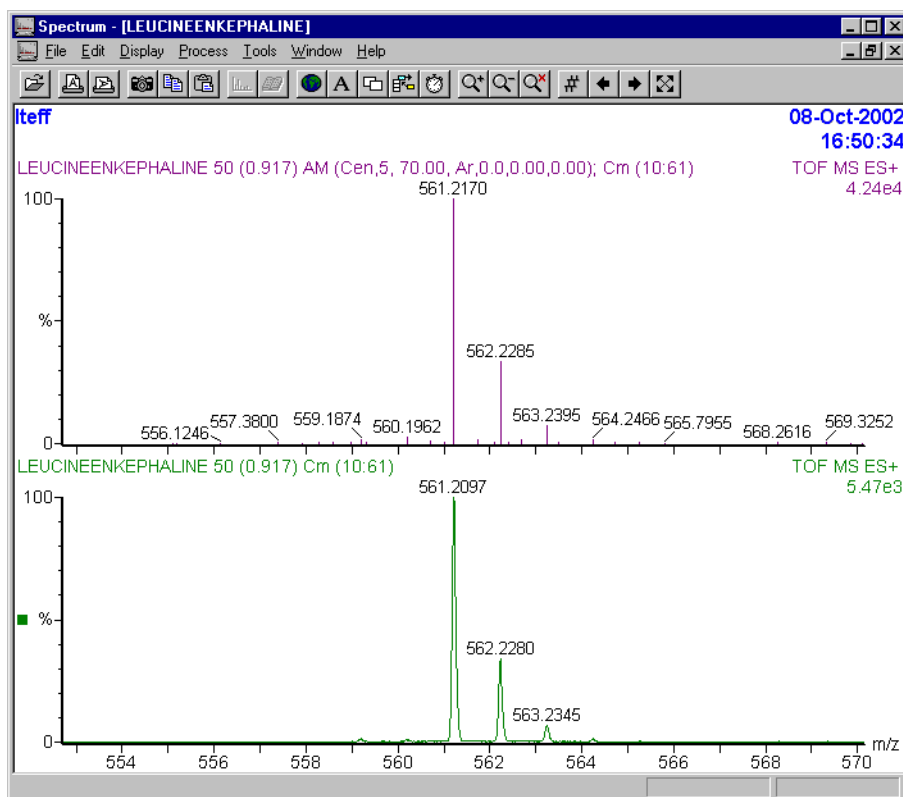


Figure 7-3 Indicated Mass of Leucine Enkephalin prior to L_{teff} Correction

7. Use this value and the actual mass of Leucine Enkephalin (556.2771) to calculate Lteff. In this case:

$$L_{\text{teff}} = 1082.78.$$

8. Enter the new Lteff value on the **TDC Settings** dialog.
9. Repeat steps 3 - 5. The new Lteff setting will give a nominal mass accuracy for Leucine Enkephalin.

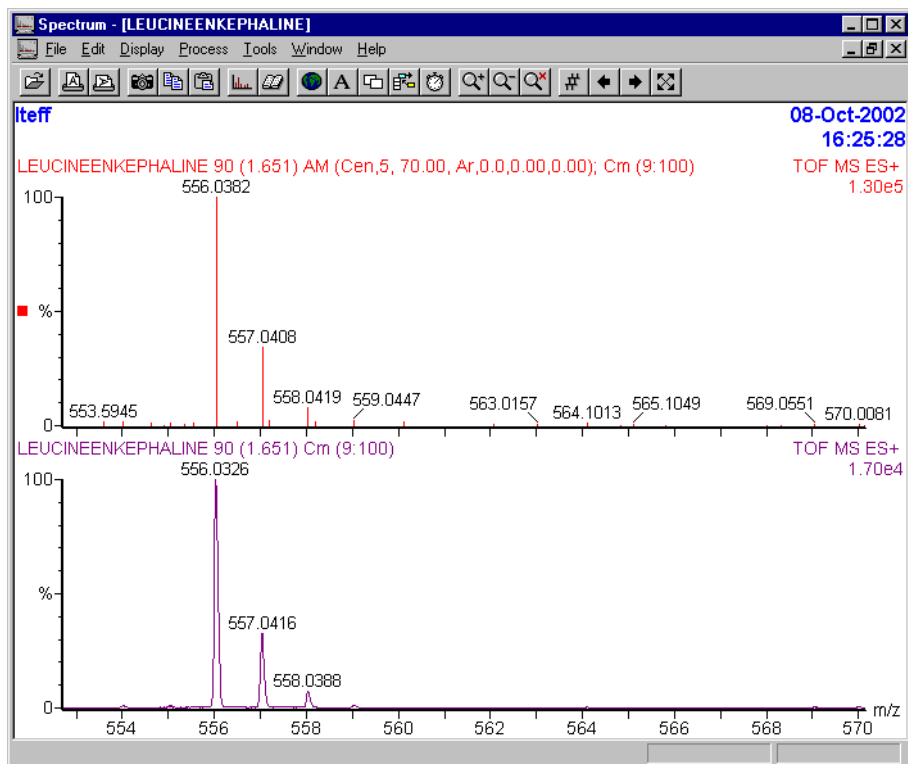


Figure 7-4 Indicated Mass of Leucine Enkephalin after Lteff Correction.

7.2.2 Acquire Data for Calibration

1. Set the Cone Voltage (**30 - 40** volts) so as to give a good distribution of PEG peaks.
2. Check the resolution as detailed in the previous chapter.
3. Acquire data for one minute over the range **100 - 1000** Da, with scan time of 1 second (for details on Acquiring Data see Section 7.2.1 on page 97).
4. Combine at least **30** scans of data.

Check that the signal is not too intense and that the data acquired does not have an ion count above **200** counts/sec.

7.2.3 Check for Deadtime Distortion

With the data acquired in the previous section:

1. Select **Process > Center** from Spectrum, this invokes the **TOF Spectrum Center** dialog (Figure 7-6). Select the **TOF** button from this dialog to invoke the **QTOF Accurate Mass** dialog. Check that both the resolution and Np multiplier are set to 0 (Figure 7-5).

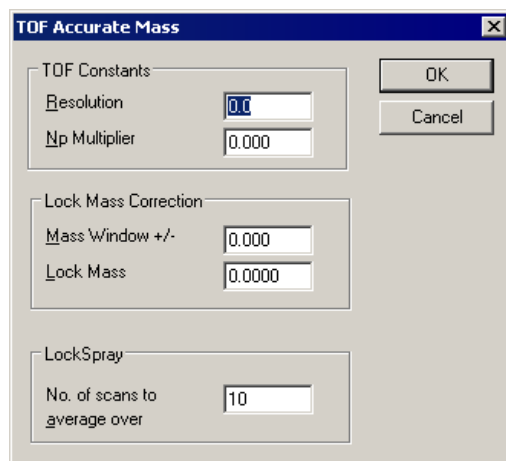


Figure 7-5 The TOF Accurate Mass dialog

2. Center the Spectrum with the values shown in (Figure 7-6).
For details on how to center spectra see the MassLynx User's Guide.

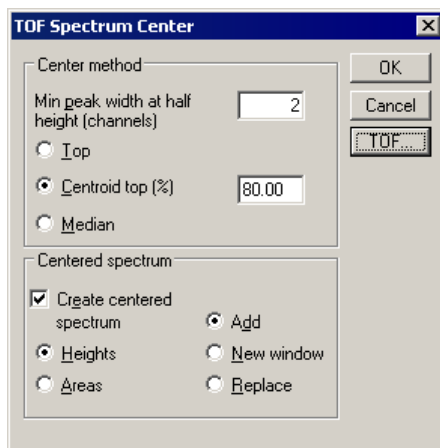


Figure 7-6 The TOF Spectrum Center Dialog

3. Re-center the data using a **Resolution** of **5000** and **Np Multiplier** of **0.7**. Check that any differences in the centroided masses of the two centered spectra are less than 1 mDa

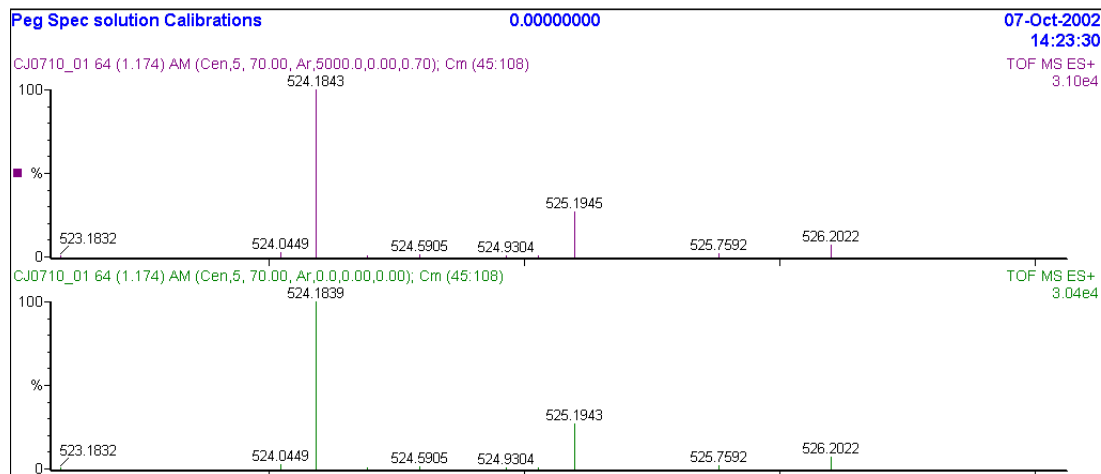


Figure 7-7 Comparison of Centered Spectra With (top) and Without (bottom) Deadtime Correction

If the difference between the two peaks is more than 1 mDa then repeat the acquisition with less intense peaks, either by dilution or by reducing the capillary voltage.

Once a satisfactory centered spectrum has been acquired, save it in the spectrum history.

7.2.4 Calibration

1. Select **Calibration > Calibrate Instrument** from the Tune Page. This invokes the Calibration Window (Figure 7-8).
2. Select the Reference file to use from the drop down list of the Calibration Window. (**pehnh4**, ref will give the correct masses with the sample used here).
3. Select **Calibrate > From File**. This will invoke the **Display Calibration Graphs** dialog (Figure 7-9).
4. Use the **Browse** option to select the centered spectrum which has been previously saved. In this example it is from the data file CJ0710_02.
5. Once the spectrum has been selected Click **OK**.

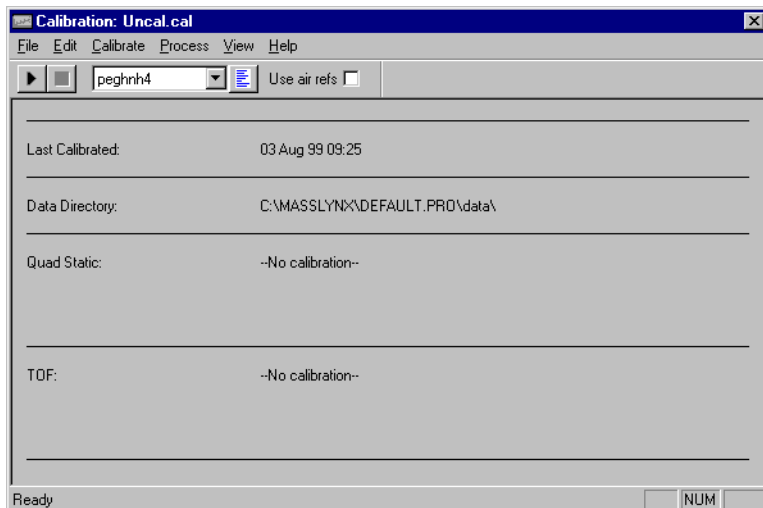


Figure 7-8 The Calibration Window.

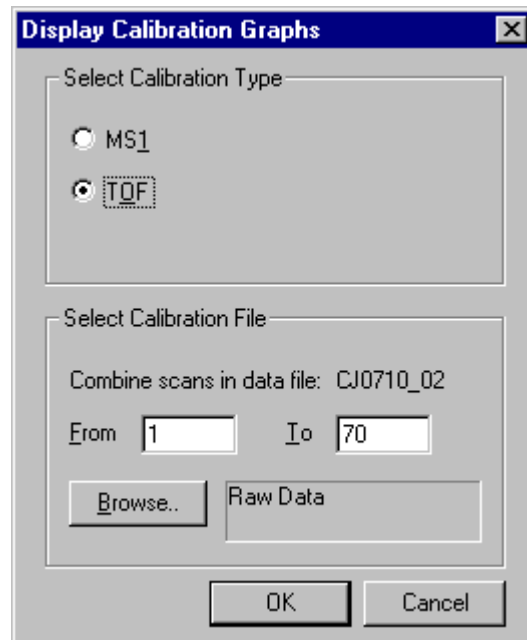


Figure 7-9 The Display Calibration Graphs Dialog

6. The Centered Spectrum will now be processed and compared against the Calibration File. When the processing is complete the results are displayed in a new calibration window.

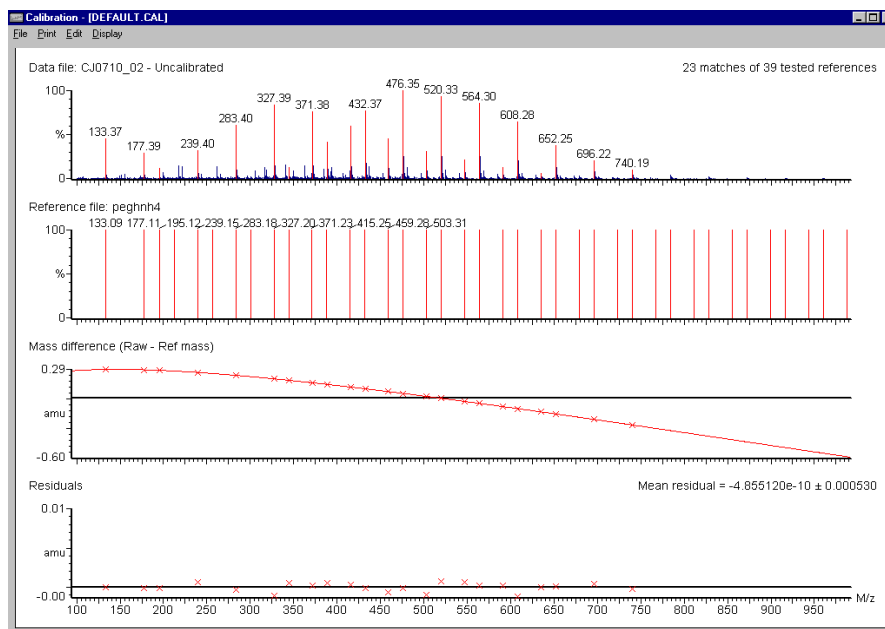


Figure 7-10 Calibration Window Showing Residual Errors

Note: The residual errors should all be less than 3 mDa.

7. Outliers can be removed from the calibration curve by right clicking the outlying peak in the reference file spectrum and clicking the associated peak on the data file spectrum.
8. The calibration parameters can be altered by selecting **Edit** to invoke the **Calibration Parameters** dialog (Figure 7-11). Set the **Intensity Threshold** to **5%**. This will mean that lower intensity peaks will not be included in the calibration and will generally give a better fit.

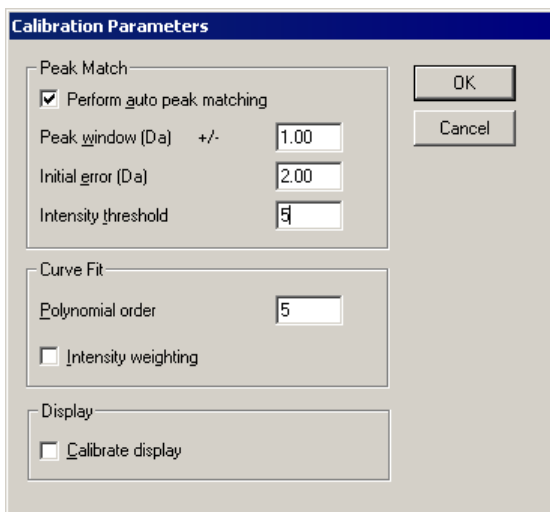


Figure 7-11 The Calibration Parameters Dialog

A fifth order polynomial has been chosen here, as a large mass range is being covered by the calibration. It can be useful to do a first order calibration initially so that outliers are spotted more easily. When these have been removed reset the curve to a fifth order.

9. When the parameters have been set **Exit** the window and Accept the Calibration (below).

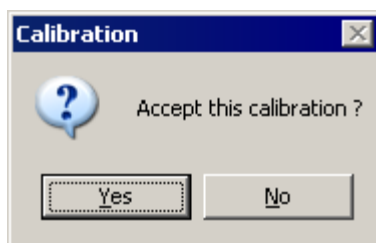


Figure 7-12 The Accept Calibration Dialog

10. The calibration will be updated with the time of the update displayed in the **Calibration Window** (Figure 7-8).

Note: Always Select **File > Save As** before exiting the calibration window and save the calibration with a unique filename. Exiting directly may result in the **Uncal** default file being overwritten.

7.2.5 LockMass Correction

The application of a single point lock mass correction will now correct for subsequent instrument drift and bring masses back to within 5 ppm RMS, on the condition that there is no isobaric chemical interference with either lock mass or analyte peaks.

The lock mass is found under the TOF Spectrum Center menu when QTOF is clicked:

Centering the spectrum with these parameters will force the peak at 520 to be 520.333 exactly and recalibrate the entire mass spectrum.

7.2.6 Deadtime Correction

If you wish to Acquire Accurate Mass data at above 200 counts/sec then a deadtime correction needs to be applied to the calibration.

1. Acquire a spectrum of PEG at the higher concentration of **10 nl/ml** and create a centered spectrum with no lock mass and no correction. (All parameters set to zero.)(Figure 7-5)

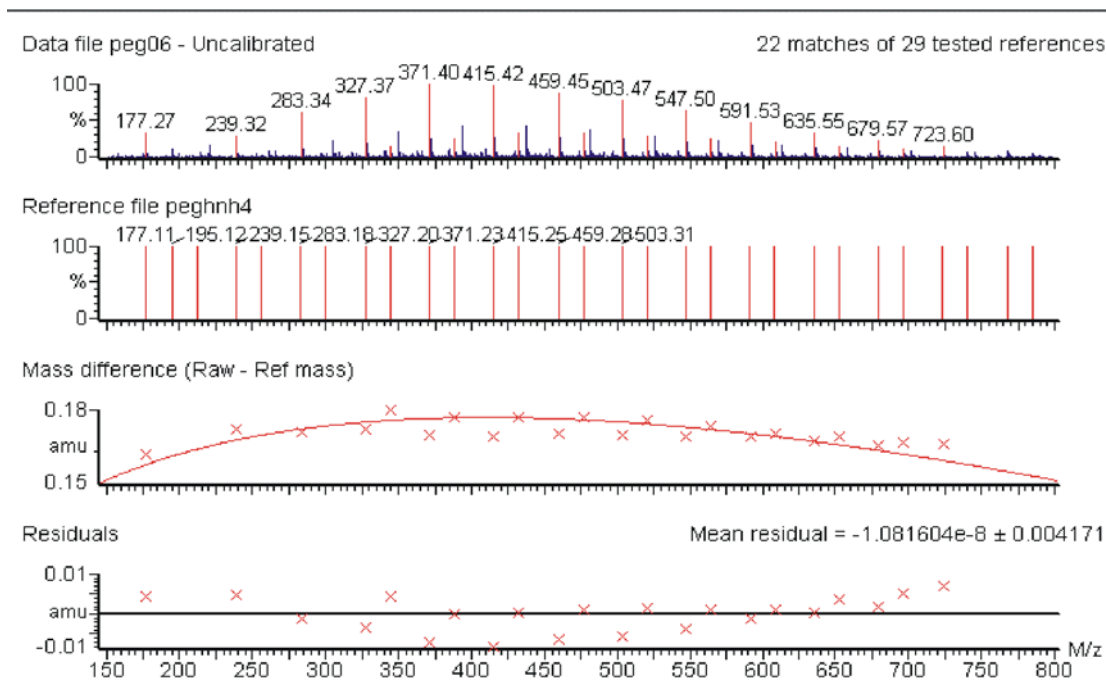


Figure 7-13 Calibration Window Deadtime Distortion

2. Follow steps for calibration in Section 7.2.4 on page 102, to when the Calibration results window is open and the residual errors are revealed.

There should now be a reduction in the mass deviation observed between adjacent peaks of high and low intensity.

3. Repeat this process at subsequently lower values of Np until the deviation is minimized.

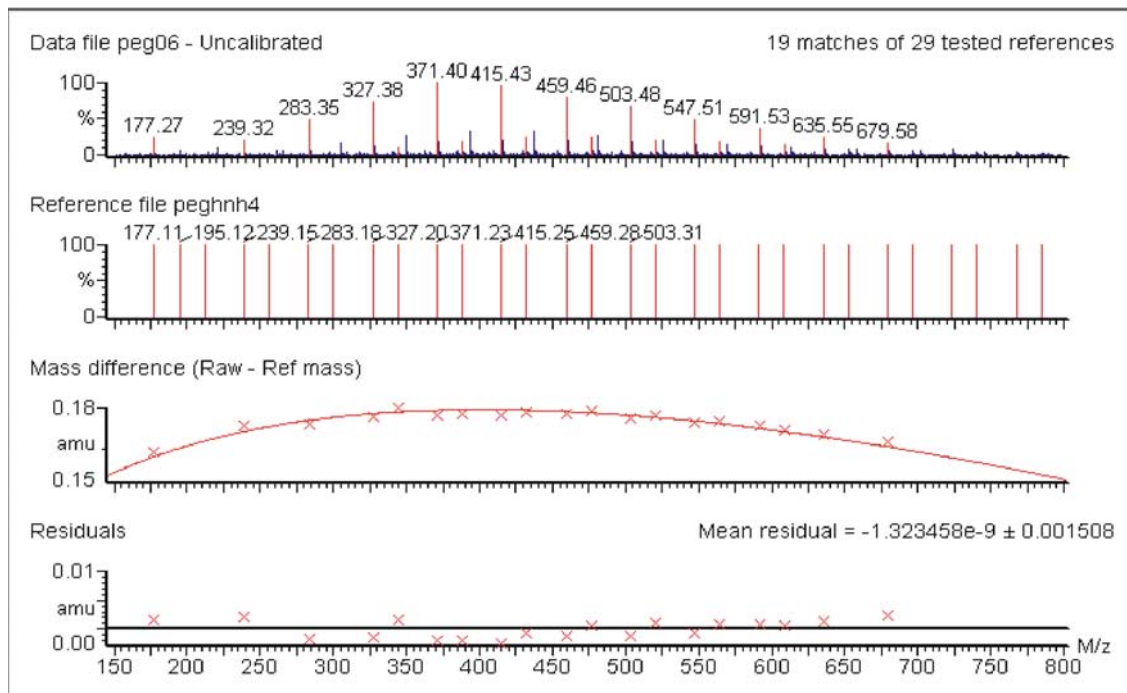


Figure 7-14 Same data as Figure 7-13 but with Np multiplier of 0.7

At this point the peaks should lie close to the axis - the line of best fit:

4. Typical values for **Np Multiplier** will be about **0.7**, and the **Resolution** used should be as measured on the instrument at mass 500.
5. Now that the deadtime correction has been set up, the difference in reported areas and masses between corrected and uncorrected centered data can be seen. In the example below, the 371 peak has been shifted by 8.7 mDa.

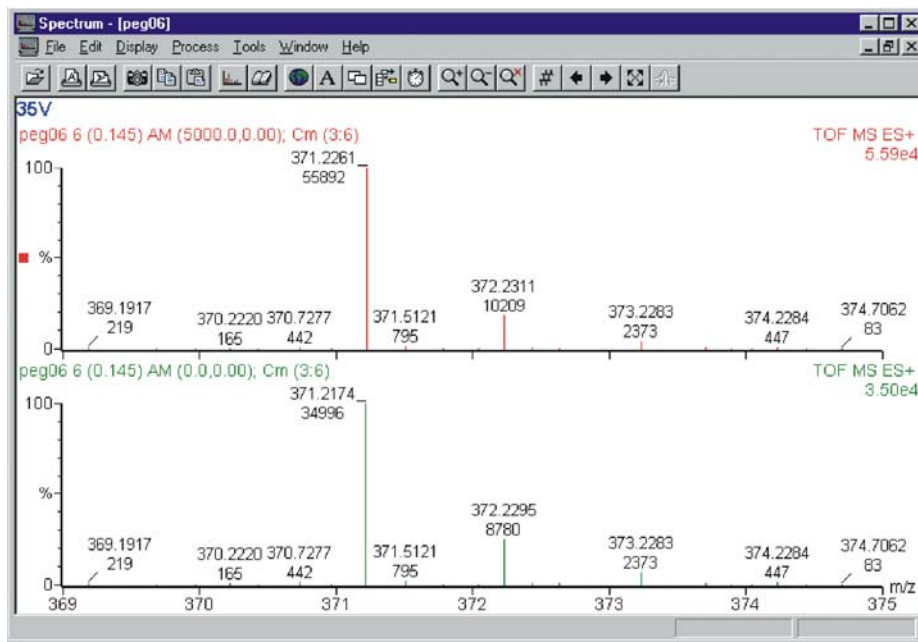


Figure 7-15 Two Spectra One (bottom) with no Deadtime Correction and the Other (top) with.

Further Notes

The model successfully corrects for deviations of up to about 15 mDa at mass 500, at which point the limits of the model are reached, with no correction applied at higher ion current.

Users should familiarize themselves with the ion current range over which successful mass measurements can be made. Once the Resolution and Np Multiplier figures have been evaluated they can be left active in the menu without affecting centered data in any adverse way.

*The procedure involved using the Spectrum, Make Calibration commands is a post-processing calibration, producing a data file of a particular format and with the extension **.scl**.*

*The file resulting from the commands Instrument Calibrate, Save Calibration has the same file extension **.cal**, but it should be noted that because it is an instrument calibration file, it is of an incompatible format to that described above.*

Folders and naming conventions should be appropriately assigned.

7.3 Exact Mass Measurement: Additional Hints

- Best results are obtained if the lock mass gives an intensity of approximately 100 to 200 counts per second as shown on the real time tune display.
- The lock mass should be chosen to be at the upper end of the mass range used.
- Do not change transfer lens voltages without re-calibrating the instrument.
- Changing Cone will not change the calibration.
- Always be aware of possible chemical interference problems, either on the sample or the lock mass peak.
- When performing a base calibration (using PEG, for example) better results may be obtained by using a stronger solution and moving the probe off axis to limit the ion current, rather than using a lower concentration.
- This has the effect of minimizing any possible chemical interference from background ions.
- Always check the stability of the spray (see the instructions for the setting of the probe tip in Maintenance and Fault Finding).
- Short term variations in the spray produce fluctuations in the number of ions per peak per pusher pulse, giving rise to errors in the deadtime correction calculations. The number of ions per peak per pusher pulse is calculated from the ion current integrated over a period of time thus only giving an average value.
- To obtain the best deadtime correction, only combine scans of a similar intensity, either at the top of a chromatographic peak or in the tail of a chromatographic peak.
- This ensures the number of ions per peak per push is calculated correctly.
- The deadtime correction algorithm can only correct for ion intensities and m shifts up to a limit. If the ion current approaches approximately 10,000 counts per peak (without correction) in one second then the limits of the model are being reached.
- If the limits of the model are exceeded no correction is applied - the same result will be obtained by centring the data with Resolution and Np Multiplier set to zero.
- If the limits of the deadtime correction algorithm are exceeded it may be possible to use the C₁₃ isotope instead.
- For the best mass accuracy when mass measuring doubly charged ions, it is advisable to use a doubly charged lock mass peak.
- To obtain the true number of ions per peak, areas must be selected on the peak centre menu.

- The standard deviation in the determination of the mass (strictly m/z) centroid of a triangular-shaped peak (σ_{ppm}) due to ion statistics alone is given by the equation below.

A triangle is assumed to be a close enough approximation to the shape of the mass spectrometer peak for the equation to be valid.

$$\sigma_{\text{ppm}} = 10^6 \Delta M / M(24N_p)^{0.5} \text{ ppm}$$

where:

ΔM is the width (m) of a triangular peak across the base.

M is the m value of the peak.

N_p is the number of ions per peak.

Using the above equation we can calculate the number of ions per peak required to give a standard deviation of 5 ppm when measuring a peak at 500 m.

Assuming 5000 (FWHM) resolution, so $\Delta M = 0.2 m/z$ (width at base = twice width at half height), then:

$$N_p = ((0.2^2 \times 10^{12}) / (24 \times 500^2 \times 5^2)) = 267$$

Thus standard deviations of less than 5 ppm cannot be expected unless the number of ions per peak is greater than 267.

Chapter 8

Data Acquisition

The following sections of this chapter will guide you through the basic acquisition methods used in MassLynx.

- Section 8.1, Starting an Acquisition
- Section 8.2, Monitoring an Acquisition
- Section 8.3, The Experiment Setup Editor
- Section 8.4, Automated Data Dependent Acquisition (DDA)

Before starting an acquisition the instrument will need to be tuned and calibrated (refer to Chapter 6, Tuning and Chapter 7, Calibration). You will also need to refer to the following chapters of the MassLynx User's Guide.

- Chapter 4 - Sample Lists
- Chapter 6 - Chromatogram
- Chapter 7 - Spectrum

8.1 Starting an Acquisition

There are two ways of starting an acquisition:

- single sample acquisition from the tune page.
- multiple sample acquisition from the MassLynx sample list.

8.1.1 Starting an Acquisition from the Tune Page

The easiest way to acquire data is directly from the tune page.

- Acquisitions can be started and stopped.
- Inlet programs cannot be used.
- Analog data cannot be acquired.
- Multiple sample sequences cannot be acquired.

To start a single sample acquisition:

1. Press Acquire on the tune page, or choose Acquire. This will invoke the **Tune Page Acquisition Dialog**.
2. Make any required changes to the settings.
3. Press Start

Figure 8-1 Tune Page Acquisition Dialog

Parameters

Data File Name	The can be up to 128 characters. If the file already exists on disk, a prompt is given to rename the file or to overwrite the existing one. The file is written to the data directory of the current project.
Set Mass	Specifies the precursor mass that is used for the TOF MSMS acquisition. This control is disabled if the function selected is TOF MS.
Start Mass and End Mass	These specify the masses at which the scan starts and stops. Start Mass must be lower than End Mass.

Run Duration	The length of the acquisition, measured in minutes.
Scan Time	Specifies the duration of each scan in seconds.
Inter Scan Time	Specifies the time in seconds between a scan finishing and the next one starting. During this period no data are stored.
Calibration	Allows an appropriate, previously generated calibration file for either positive or negative ion to be selected - see chapter entitled Calibration and Exact Mass.

Pressing **Origin** allows additional information about the sample to be analyzed to be entered into the following fields:

- Submitter
- Job
- Task
- Conditions

To change the directory into which data are acquired:

1. Cancel the acquisition.
2. Create a new project by choosing MassLynx top level file menu, Project Wizard or open an existing one by choosing Open Project, from the MassLynx top level file menu.
3. The Text area is used to enter the sample description. The description can be displayed on any output of the acquired data and has a maximum length of 74 characters. To display text on more than one line press CTRL+Return at the end of a line.
4. The type of acquisition Function used to collect the data can be either of the following:
 - ToF MS
 - ToF MS/MS

More information is given in Function List Editor later in this chapter.

The Data Format that are collected and stored on disk can be any of the following:

- Centroid
- Continuum
- MCA

More information is given on data formats later on in this chapter.

8.1.2 Multiple Samples

The MassLynx top level screen contains a sample list editor for defining multiple experiments. The list of samples is set up using a spreadsheet style editor, which can be tailored to suit different requirements.

For more information see the MassLynx User's Guide, Chapter 4.

8.2 Monitoring an Acquisition

Acquisition status is shown on the MassLynx screen. The run time is shown on the MS panel and the scan status, sample number and scan number are shown on the Status bar at the bottom of the page.

8.2.1 The Acquisition Status Window

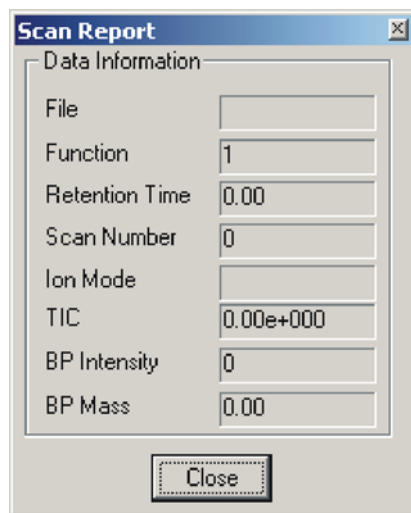


Figure 8-2 The Scan Report Dialog

The acquisition status window, or scan report (Figure 8-2), provides a scan by scan statistical report of the progress of an acquisition.


To display the scan report dialog, from the tune page:

Select **Options > Scan Status**.

This shows details of the scan currently being acquired.

8.2.2 Chromatogram Real-Time Update


To view in real time the chromatogram that is currently being acquired:

1. Open the data file using the MassLynx data browser.
2. Press  or select **Real-Time Update** from the **Display** menu. The chromatogram display is updated as the acquisition proceeds.

For Further Details see the MassLynx User's Guide

8.2.3 Spectrum Real-Time Update

To view in real time the spectrum that is currently being acquired:

1. Open the data file using the MassLynx data browser.
2. Press  or select **Real-Time Update** from the **Display** menu.
3. Select **Enable Real-Time update**.

For Further Details see the MassLynx User's Guide

8.2.4 System Manager

To check the communications between the MassLynx software and the embedded PC:

Select **Options > Communications Status** from the Tune page menu bar.

8.2.5 Stopping an Acquisition

To halt the acquisition:

1. From the Tune page, press .
2. From the MassLynx screen choose Stop from the Run menu, or press .

Data acquired up to this point is saved.

8.3 The Experiment Setup Editor

The experiment setup editor is used to set up the function(s) that the mass spectrometer uses during an acquisition. A function list can be a mixture of MS or MSMS experiments

that can be arranged to run either sequentially or concurrently during an acquisition. Survey or parent functions must be run in isolation.

Typical uses for mixed function acquisitions are the acquisition of several MSMS product ion spectra over different retention windows.

A function list is produced, saved on disk and then referenced by name when an acquisition is started.

A simple function list is shown above, containing only one function: a centroided mode TOF MS scan, between 100 and 1500 amu using ES+ ionization. Immediately above the function bar display is a time scale that shows from when the function is active, and for how long it runs. In this case the acquisition runs from 0 to 60 minutes.

To access this dialog:

Select MS Method from the MassLynx **Instrument Shortcut** bar.

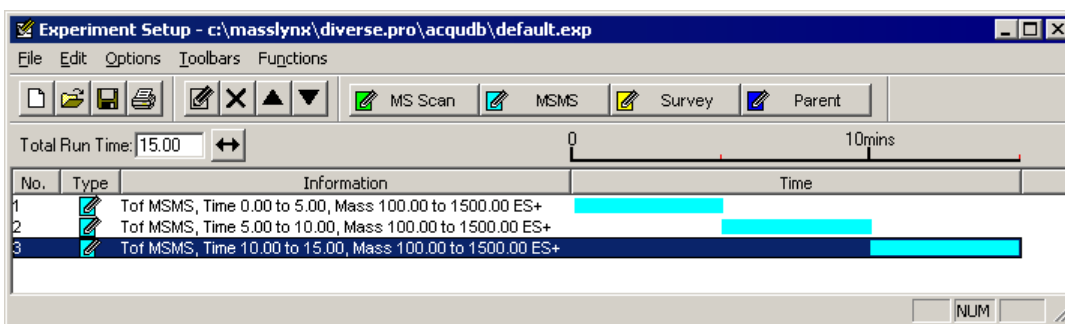


Figure 8-3 The Experiment Setup Editor

One example of a function list shown in Figure 8-3 shows 3 different MSMS functions in sequential time windows is shown above. The currently selected function is highlighted and enclosed in a rectangular frame. If the display shows more than one function a new function can be selected either by clicking with the mouse, or by using the arrow keys on the keyboard.

8.3.1 The Experiment Setup Editor Toolbar

The toolbar is displayed at the top of the tune window and allows some common operations to be performed with a single click.



Create a new function list.



Edit the selected function.



Open an existing function list.



Delete the selected function.



Save the current function list to disk.



Move the selected function up the list of functions.



Print the current window in portrait format.



Move the selected function down the list of functions.



Create a new function of the Indicated type.

8.3.2 Adding a New Function

To add a new function to the list:

1. Click one of the toolbar buttons, or select the required function from the **Functions** menu.


The editor for the function type selected is displayed showing default values.

2. Make any changes required to the parameters and press **OK** to add the new function.

The function editors for each scan type is discussed in detail later on in this chapter.


8.3.3 Modifying an Existing Function

To modify an existing function:

1. Select the function in the function list.
2. Press , or double click on the function.

This displays the appropriate editor for the function type and allows changes to be made.

The function list display is updated to show any changes.

*Entering a new value in **Total Run Time** and pressing  sets the maximum retention time for the experiment. The ratio of the functions defined is maintained. For example, if two functions are defined one from 0 to 5 minutes and the other 5 to 10 minutes then a **Total Run Time** of 10 minutes is displayed. If this value is changed to 20 then the first function now runs from 0 to 10 minutes and the second from 10 to 20 minutes.*


8.3.4 Copying an Existing Function

To copy an existing function:

1. Select the function in the function list.
2. Select **Copy** and then **Paste** from the **Edit** menu.
3. Modify the parameters as described above.

8.3.5 Removing a Function



To remove a function:

1. Select the function in the function list.
2. Press , choose **Delete** from the **Edit** menu, or press **Del** on the keyboard.
3. When asked to confirm the deletion, select **Yes**.

8.3.6 Changing the Order of Functions

Functions are displayed in ascending **Start Time** and **End Time** order and this order cannot be changed. For functions that have the same start and end time the order in which they are performed can be changed as follows:

Highlight the required function.

Press  or  repeatedly until the function is in the required position.

8.3.7 Setting a Solvent Delay

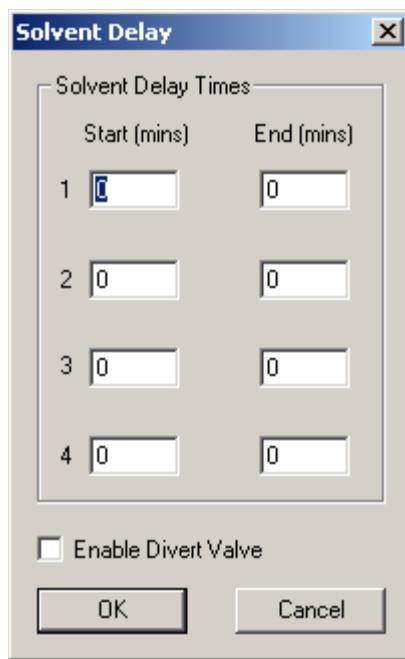


Figure 8-4 The Solvent Delay Dialog

No data is stored during the solvent delay period, which means that solvent peaks that would normally be seen eluting on the TIC chromatogram are no longer seen.

For APcI functions the APcI probe temperature is set to the value specified in the **APcI Probe Temp** control for the period of the solvent delay.

To set a solvent delay for a function list:

Select **Options > Solvent Delay** from the menu bar.

8.3.8 Analog Channels

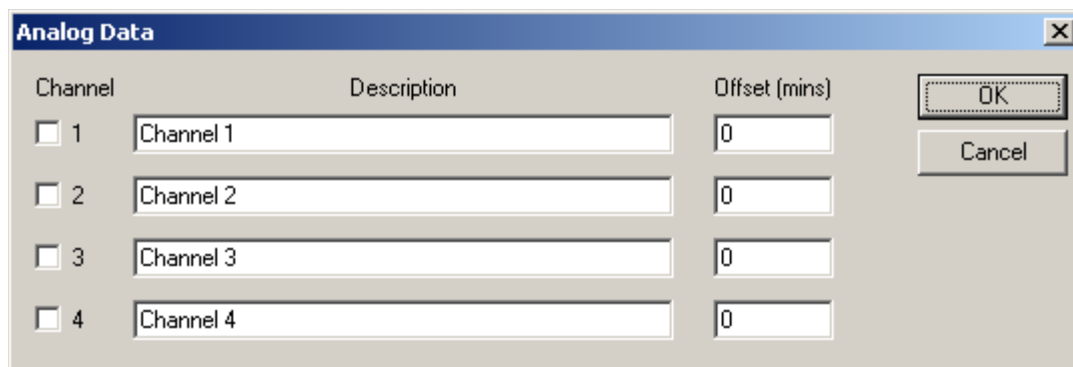


Figure 8-5 The Analog Data Dialog

If an analog channel's hardware option is fitted, up to 4 channels of analog data can be acquired, which are stored with the data acquired from the mass spectrometer. Analog channels are typically used to collect data from external units such as UV detectors or FID detectors. A reading is made from the external channel at the end of each scan and stored with the data for that scan. The resolution of the chromatography for an analog channel is therefore dependent on the scan speed used to acquire the mass spectrometry data.

To access this dialog:

Select **Options > Analog Data** from the Method Editor dialog.

To store data for an analog channel:

1. Check the box(es) for the channel required.
2. Enter a textual description for each of the selected analog channels.

This description is used on the analog chromatogram dialog as the channel description. See "Chromatogram" in the MassLynx User's Guide.

3. Enter an **Offset** to align the external unit with the mass spectrometer.
4. Press **OK**.

8.3.9 Saving and Restoring a Function List

To save a function list:

1. Choose Save As from the function list File menu.
2. Enter a new file name, or select an existing file from the list displayed.

3. Press Save.

If the file already exists on disk, confirmation is requested to overwrite the existing information.

Press Yes to overwrite the file, or No to select a different name.


When the editor is closed a prompt is issued to save any changed function lists.

To restore a saved function list:

1. Choose Open from the function list File menu.
2. Select the name of the function list to open, either by typing its name or by selecting it from the displayed list.
3. Press Open.

8.3.10 Setting up an MS Scan Function

The full scan function editor, activated by pressing or by selecting TOFMS Scan from the Functions menu, is used to set up centroid, continuum and MCA functions.

The full scan function editor, activated by pressing  MS Scan or by selecting **TOFMS Scan** from the **Functions** menu, is used to set up centroid, continuum and MCA functions.

Function:4 TOFMS Scan

Mass (m/z)
 Low Mass (m/z) 100
 High Mass (m/z) 1500

Method
 Ionization Mode ES+
 Data Continuum
 Scans To Sum 100000

Time (Mins)
 Start Time (Mins) 0
 End Time (Mins) 60

Scan Duration (secs)
 Scan Time 0.50
 Inter-Scan Delay 0.1

Cone Voltage
 Use Tune Page Setting ☐
 Cone Voltage (V) 30

APCI Probe
 Use Tune Page Setting ☐
 Probe Temperature 20

Reference Scan
 Scans to average 1
 Frequency (secs) 0
 Cone Voltage (V) 0

OK Cancel

Figure 8-6 The TOFMS Scan Dialog

Mass (m/z)

Low Mass and **High Mass** specify the masses at which the scan starts and stops.

Time (mins)

Start Time and **End Time** specify the retention time in minutes during which this function becomes active, and data are acquired.

Cone Voltage

When **Use Tune Page** is checked, the cone voltage set on the tune page at the start of the acquisition is used.

The cone voltage value cannot be altered during acquisition by typing new values into the tune page, since the new values are not downloaded during acquisition. This can only be done by acquiring from the tune page.

Method

Ionization Mode, specifies the ionization mode and polarity to be used during acquisition.

Data, specifies the type of data to be collected and stored on disk.

There are three options:

- **Centroid** stores data as centroided, intensity and mass assigned peaks. Data are stored for every scan.
- **Continuum**. The signal received by the interface electronics is stored regularly to give an analog intensity picture of the data being acquired. Data are not centroided into peaks, but are stored for every scan.

Due to the fact that data are acquired to disk at all times, even when no peaks are being acquired, continuum data acquisition places some extra burden on the acquisition system as compared to centroided acquisition. Data files tend to be significantly larger than centroided ones and the absolute scanning speed (amu/sec) is slower.

It is possible, however, to set a threshold below which the data are not stored. Depending on the nature of the data acquired, this can greatly reduce these effects. The threshold can be set so that data considered to be ‘noise’ can be discarded, thus improving data acquisition speed and reducing data file sizes.

- **Multi Channel Analysis (MCA)**. MCA data can be thought of as ‘summed continuum’, with only one intensity accumulated scan being stored for a given experiment. As each scan is acquired, its intensity data is added to the accumulated summed data of previous scans.

An advantage of MCA is that random noise does not accumulate as rapidly as real data and therefore effectively averages out over a number of scans. This emphasizes the real data and improves the signal to noise ratio.

The disadvantage of MCA is that, as there is only one scan, it cannot be used for time resolved data.

For MCA, **Scans to Sum** defines the number of scans to sum to create a spectrum.

Scan Duration (secs)

Scan Time specifies the duration of each scan in seconds while **Inter-Scan Delay** specifies the time in seconds between a scan finishing and the next one starting. During this period no data are stored.

APCI Probe

Probe Temp, in degrees centigrade, is enabled when **Ionization Mode** is set to APCI.

When **Use Tune Page Settings** is selected the APCI probe temperature set on the tune page at the start of the acquisition is used. This control is enabled when the ionization mode is set to APCI.

The APCI probe temperature value cannot be altered by typing new values into tune page during the acquisition since the new values are not downloaded during the acquisition. This can be done by acquiring from the tune page.

Setting up MS-MS Scanning Functions

The screenshot shows the 'Function:4 MSMS Scan' dialog box. It contains several sections with input fields and checkboxes:

- Masses:** Set Mass (m/z) 50, Low Mass (m/z) 100, High Mass (m/z) 1500.
- Method:** Ionization Mode ES+, Data Continuum, Scans To Sum 100000.
- Time (Mins):** Start Time (Mins) 0, End Time (Mins) 60.
- Scan Duration (secs):** Scan Time 0.50, Inter-Scan Delay 0.1.
- Cone Voltage:** Use Tune Page Setting ☐, Cone Voltage (V) 30.
- Collision Energy:** Use Tune Page Setting ☒, Collision Energy (V) 30.
- Reference Scan:** Scans to average 1, Frequency (secs) 0, Set Mass (m/z) 0, Cone Voltage (V) 0, Collision Energy (V) 0.
- APCI Probe:** Use Tune Page Setting ☐, Probe Temperature 20.

At the bottom right are 'OK' and 'Cancel' buttons.

Figure 8-7 The MSMS Scan Dialog

The **MSMS** scan function editor is activated by pressing the MSMS button or by selecting **TOF MSMS** from the Functions menu.

Many of the fields in the MSMS scan editor are similar to those described in Section 8.3.10 on page 121. Only fields that differ significantly are described below.

Masses

Set Mass specifies the precursor mass to be selected by the quadrupole for fragmentation in the gas cell.

When multiply charged species are being analyzed care should be taken when selecting the **High Mass** to ensure that the mass range is sufficient to include product ions with a lower charge but a higher m/z value than the precursor.

Time

The product ion spectra are acquired in the time window defined by the **Start Time** and **End Time**. Several MSMS scanning functions can be set up to either run sequentially or concurrently in overlapping time windows.

Collision Energy

When **Use Tune Page Settings** is checked, the collision energy set on the tune page at the start of the acquisition is used.

8.4 Automated Data Dependent Acquisition (DDA)

For the automated MS/MS analysis of unknown compounds, the Micromass Q-Tof *micro* has a powerful software control system to enable the instrument to perform data dependent acquisition (DDA), switching from the MS to MS/MS mode and then returning to the MS mode using data dependent criteria. The advantage of this method is that it removes the requirement to analyze the sample in MS mode to identify the target precursor ions and then re-run the sample in MS/MS mode to acquire the MS/MS data from each of these precursors. This is particularly valuable in the analysis of unknown samples using on-line chromatography where the target precursor ions and their retention times may well be quite different for each sample.

During acquisition the instrument is controlled by the MS FILE in the MassLynx Sample List. DDA can only be acquired using the Sample List and may not be started by acquiring from the Tune Page.

The MS FILE file is created using the method editor and then selected in the sample list in the usual way.

8.4.1 Creating a Data Dependent MS to MS/MS Switching Experiment

Accessing the Method Editor

In the MassLynx sample list, using the cursor to select the MS file, perform a right mouse button click and in the drop down menu select the Open command to access the Method editor for that MS file.

Creating a New Experiment

In the top banner of the Method editor window select File and in the drop down menu select New to create a new template with no predefined functions. Select the Survey button to create a data dependant acquisition.

Clicking the Survey button will launch Survey Scan Dialog (Figure 8-8). There are 8 tabbed pages available. These are used to access all the relevant parameters for the data dependent function switching experiment.

It is advisable to work through each page sequentially and fill in the appropriate values.

8.4.2 Survey

Function:1 Survey Scan

Include Exclude Collision Energy Adduct
Q-Tof Micro Survey MS to MSMS MSMS to MS MSMS Template

Masses
Low Mass (m/z) 100
High Mass (m/z) 1500

Method
Ionization Mode ES+
Data Continuum

Times
Start Time (Mins) 0
End Time (Mins) 60

Scan Duration
Scan Time (secs) 0.5
Inter-Scan (secs) 0.1

Cone Voltage
Use Tune Page Setting ☐
Cone Voltage (V) 30

Reference Scan
Scans to average 1
Frequency (s) 0
Cone Voltage (V) 0

OK Cancel Apply

Figure 8-8 Survey Scan Dialog: Qtof MicroSurvey Page

The mass range is used to define the start and end mass range for the MS data.

Time (Mins), defines the start and end times for the complete experiment. The spectrum integration period (**Scan Duration**) and **Inter-Scan Delay** period for the MS data are selected. In the MS mode the minimum time required to transfer data from the embedded PC to the data file on the host PC is about 100 mSec. Scan times of less than 1 sec may be employed if required, but because the Inter-Scan Delay is fixed (by the data transfer rate) the duty cycle will decline at very fast acquisition rates.

The **Cone Voltage** must either be set to **Use Tune Page** or set in this window to a suitable value.

8.4.3 MS to MSMS

Function:1 Survey Scan

Include Exclude Collision Energy Adduct

Q-ToF Micro Survey MS to MSMS MSMS to MS MSMS Template

MS to MSMS Switch Criteria

☐ TIC

☒ Intensity

Threshold 10

Detection Window (Da) 0.5

Number of Components 1

Retention Time Window (s) 1

Charge State

Tolerance Window +/- (Da) 3

Extraction Window (Da) 2

Precursor Selection

☒ Everything

☐ Included Masses only

☐ Included Masses Take Priority

Detected Precursor Inclusion

☒ Auto Exclude

☐ Always Include

☐ Include After Time (s) 10

Data

☐ Discard uninteresting survey scans

OK Cancel Apply

Figure 8-9 Survey Scan Dialog: MS to MSMS Page

This tab has all the appropriate values that dictate when to switch into the MS/MS mode from the MS mode. They are as follows: -

MS to MSMS Switch Criteria

The decision can be made using **TIC** or individual ion **Intensity**. Since multiple precursors may be selected for concurrent MS/MS experiments it is usual to set this to **Intensity**. It is useful to acquire some MS data to be able to measure the typical background level and the peak intensities from the sample peaks so that the detection Threshold can be set appropriately.

Within the mass range of the MS survey scan any peaks that rise above the Threshold (in ion counts/sec) will be identified as candidate precursor ions for selection.

To prevent the selection of the P+1 (^{13}C) or P+2 ($^{13}\text{C}_2$) isotope peaks as additional candidate precursors a Detection Window is set. The window is an absolute mass value, i.e:- a setting of 1Da will extract a 1Da range of the m scale and identify only the most intense peak contained within that specific m range. A suitable value for 1+ and 2+ ions would be 2Da.

Number of Components defines the number of precursor masses that can be selected for concurrent MS/MS data acquisition from one MS survey scan. Once the candidate precursor are identified the N most intense will be used for MS/MS. The maximum that can be specified is eight, but often a more suitable number for LC peaks about 10-20 sec wide at half height is 3 - 4.

The detection criteria set here may be subsequently modified by charge state considerations (see later)

Retention Time Window (seconds) This parameter is only active when an include or exclude file has been generated and is selected (see later). If an include or exclude file is utilized where a retention time is included (e.g. 18 mins), then this window is a time value in seconds which is added to the retention time during which the peak must elute to be selected for fragmentation. For example, a value of 30secs would mean the peak expected to elute at 18 mins must elute between 18.00 and 18.50 mins).

Charge State

These settings define the parameters used for the charge state recognition algorithm.

Tolerance Window. The window used when measuring between adjacent peaks (isotopes). This should be a tight tolerance, and as a default set to 0.1Da.

Mass Extraction Window. This value in Da is the mass range around the peak identified that is extracted to perform the CS recognition on. For peptide work a value of 2Da is recommended as a default.

(%) Intensity of Main Peak. This threshold value is utilized to remove low intensity noise peaks from being identified by the CS algorithm. A default value of 2% is recommended.

Precursor Selection

Everything. No user-predefined masses are used and all peaks will be selected automatically according to the criteria defined in the Automatic Function Switching parameters.

Included Masses Only will perform MS/MS solely on those masses appearing in the Include Masses list, or file, (see later section page 7) provided that their intensities rise above the predefined threshold set in the MS to MS/MS switch criteria.

Included Masses and Everything Else allows masses in the Include Masses list to take precedence. These will be selected even if there are more intense peaks that would have been selected automatically.

Detected Precursor Inclusion

Auto Exclude will exclude a m value and its associated detection window (which is now defined under the exclude tab) from subsequent selection.

Always Include allows the immediate selection of a previously selected m value.

Include After Time excludes a previously selected m value and its associated detection window, (which is now defined under the exclude tab) from subsequent selection for the period specified. This prevents multiple selection of the same precursor in the event of tailing from strong LC peaks, but allows the fragmentation of near isobaric components eluting at other retention times. The maximum value is 200 secs.

Data

Discard uninteresting survey scans stores only those survey scans which triggered an MS/MS experiment to disk. Depending on the exact nature of the experiment, discarding these scans may significantly reduce the data file size.

8.4.4 MSMS to MS

Q-ToF Micro Survey | MS to MSMS | MSMS to MS | MSMS Template

MSMS to MS Switch Criteria

☒ TIC falling below threshold

☐ TIC rising above threshold

☐ Intensity falling below threshold

Threshold (Counts/s) 0

MSMS to MS Switch Method

☒ Default

☐ After time (s) 10

Figure 8-10 Survey Scan Dialog: MSMS to MS Page

If the MSMS to MS switch method is set to **Default** then the mode change will be triggered by one of the mutually exclusive options selected in the MSMS to MS Switch Criteria. The **Intensity falling below threshold** option is measured on the base peak in the MS/MS spectrum. A value of 2-3 counts/sec is recommended for peptide work on Q-ToF micro.

If the MSMS to MS switch method is set to **After time** then the mode change will be as for the default setting, except that if these criteria are not met then a maximum time can be set. (For example in the above case the switch would occur after 7 seconds assuming the intensity of the base peak in the MS/MS spectrum had not fallen below 2 counts/sec).

These rules apply on a per precursor basis. For example if MS/MS were concurrently running on 3 precursors and one MS/MS spectrum fell below 2 counts/sec then MS/MS would be terminated on that precursor. The remaining 2 precursors would continue in MS/MS mode until either they fell below the 2 counts/sec threshold or each MS/MS precursor ion had been acquired for 7 seconds. If the intensity **Threshold** were set to 0 counts/sec then MS/MS acquisition would continue on all three functions until a total of 21 seconds had elapsed.

It is useful to set this to be equal to an integral number of cycles of MS/MS spectrum acquisition. For example if the MS/MS **Scan Time** is 1 sec with an **Inter Scan Delay** of

0.1 sec and CE profiling is in operation then to get 3 scans at each of three different collision energies would require the **After time** parameter to be set to 9.9 sec.

8.4.5 MSMS Template

Function:1 Survey Scan

Include Exclude Collision Energy Adduct
Q-ToF Micro Survey MS to MSMS MSMS to MS MSMS Template

Instrument Parameter File
File c:\masslynx\default.pro\acqddb\default.ipr Browse

Mass (m/z)
Start 50
End 1500

Method
Ionization Mode ES+
Data Continuum

Cone Voltage
☐ Use Tune Page
Cone Voltage (V) 30

Scan Duration (secs)
Scan Time 1
Inter-Scan Delay 0.1

Collision Energy
☐ Use Tune Page
Collision Energy (V) 30

OK Cancel Apply

Figure 8-11 Survey Scan Dialog: MSMS Template Page

This page is used to set-up the parameters for the acquisition of the MS/MS data from each of the precursors.

It is important that the correct Q-ToF tune file (.ipr file) is selected in the instrument parameter window,

The mass range for the acquisition of the MS/MS data is independent of the mass range selected for the MS survey data as in the example shown here. The spectrum integration (**Scan Duration**) and **Inter-Scan Delay** period can also be different to those values chosen in the MS survey mode.

The **Collision Energy** parameter will define a default collision energy that will be used if, for any reason, a value cannot be selected from the parameters defined later under the collision energy settings.

8.4.6 Include Masses

Function:1 Survey Scan

Q-ToF Micro Survey MS to MSMS MSMS to MS MSMS Template

Include Exclude Collision Energy Adduct

Include Mass

☐ Range

☐ File

Mass	Retention Time	Collision Energy	Cone Voltage	Charge State	Scan

Include Window +/- (mDa)

Charge State

☒ Use Include By Charge State

Charge State(s): Number of Include Components:

Figure 8-12 Survey Scan Dialog: Include Masses Page

The **Include Mass** option allows a number of predefined m values to be included in the decision making process for switching from MS to MS/MS mode. There are two distinct modes of operation, which are defined in the MS to MSMS criteria (Section 8.4.4).

Included Masses Only will restrict MS/MS data acquisition only to peaks having the same m as those defined in the Include Masses List, provided that these peaks rise above the threshold.

Included Masses & Everything else will preferentially perform MS/MS on masses in the **Included Masses** List, even though more intense masses may have been automatically detected, provided that the included masses are above the detection threshold.

The **Listed** option allows a list of m values separated by commas to be entered directly into the window. MS/MS experiments on these masses will be carried out using the parameters defined elsewhere.

In contrast, the **File** option allows included masses to be read from a file and specific parameters (**Retention Time, Collision Energy, Cone Voltage, Charge State**) may be associated with these m values to override those parameter set elsewhere. The files may be created or modified using the **New, Add, Delete, Save As and Save** buttons as necessary.

Charge State

The **Use Include By Charge State** option allows the selection of only those charge states specified to be selected for MS/MS analysis. This may be useful for example where singly charged ions are usually from the background and it is desirable, therefore, to preclude them from analysis.

The required charge states are specified and separated by commas as shown. It is possible to enter values greater than 4, but at present CE control is available only up to and including charge state 4.

The **Number of Include Components** option allows for the maximum number of components above the detection threshold from which the candidate precursors will be selected. For low level studies where singly charged ions are not to be selected for MS/MS it is clear that there may be many candidates above the threshold, among which the multiply charged ions may be the minor components. The range for this parameter is 1 - 100 and a value of 60 or 80 would be a reasonable default value. However it should be noted that higher values slow down data processing and hence acquisition times.

8.4.7 Collision Energy

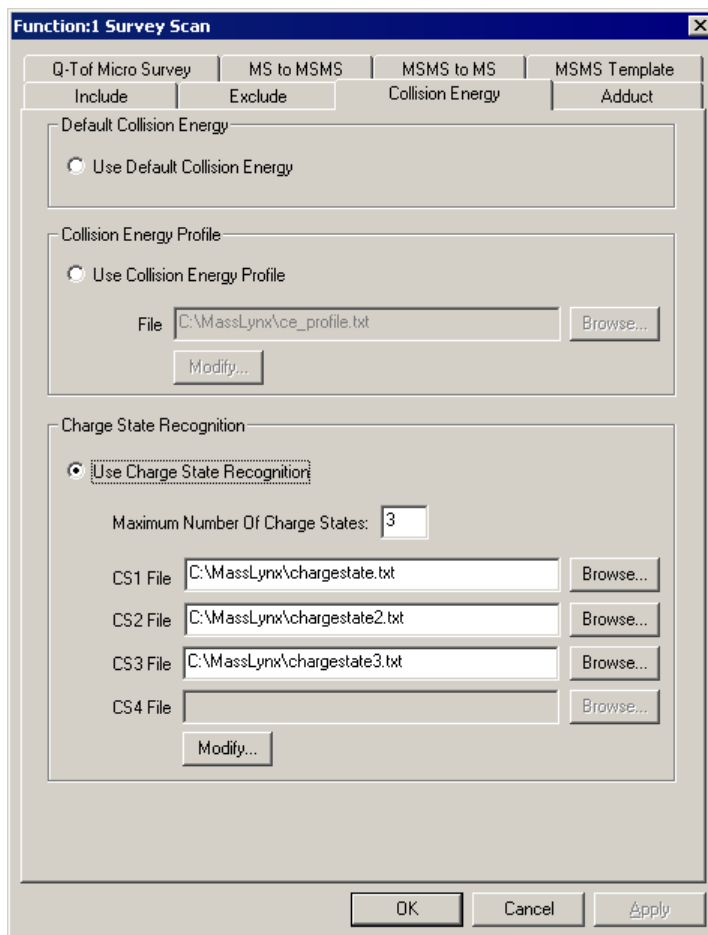


Figure 8-13 Survey Scan Dialog: Collision Energy Page

Use Default Collision Energy

When selected, all MS/MS data are acquired using the default collision energy specified in the MSMS template window (see Section 8.4.5).

Use Collision Energy Profile

This option allows a range of different collision energies to be used depending on the *m* value of the selected precursor. In addition, from 1 up to 5 collision energies may be used

for each precursor m to attempt to produce optimal fragmentation. Using the **Browse** button, files may be selected which have previously created collision energy profiles.

Using the **Modify** button new profiles can be created using the **Modify CE Profile** window and its associated **Add**, **Delete**, **New**, **Save** and **Save As** options.

The example shown is a typical set of collision energy profiles for peptides or protein digest applications. During the MS/MS experiment data are collected with the collision energy cycling through CE1, CE2....CEN on a per spectrum basis.

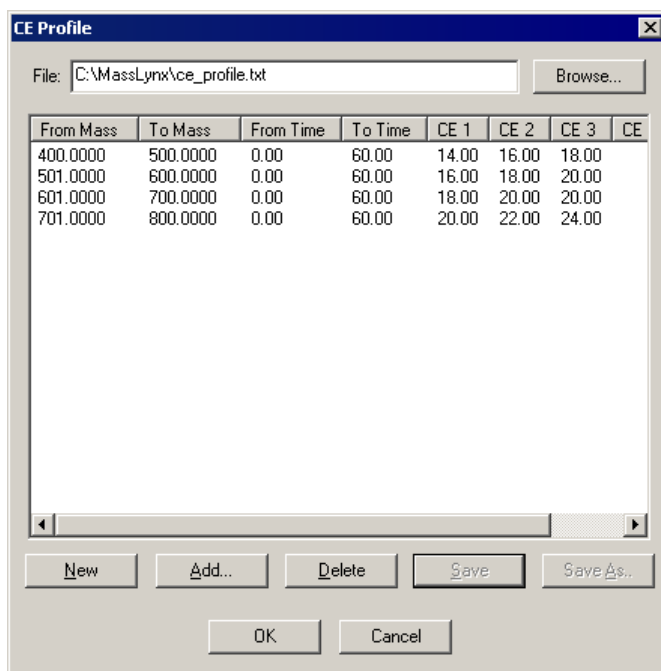


Figure 8-14 The CE Profile Dialog

Use Charge State Recognition

This mode will allow the collision energy to be set according the charge state (z) and m of the precursor as determined from the charge state recognition algorithm. Collision energy files for charge states from 1, 2, 3 and 4 may be created and selected using the Browse button. In this mode collision energy profiling by cycling around a number of different energies is less likely to be required for peptide analysis.

Using the **Modify** button new profiles can be created using the **Modify Charge State** window and its associated **Add**, **Delete**, **New**, **Save** and **Save As** options. The examples shown are a set of suggested collision energy profiles for peptides or protein digest applications when selecting 2+ ions for MS/MS analysis.

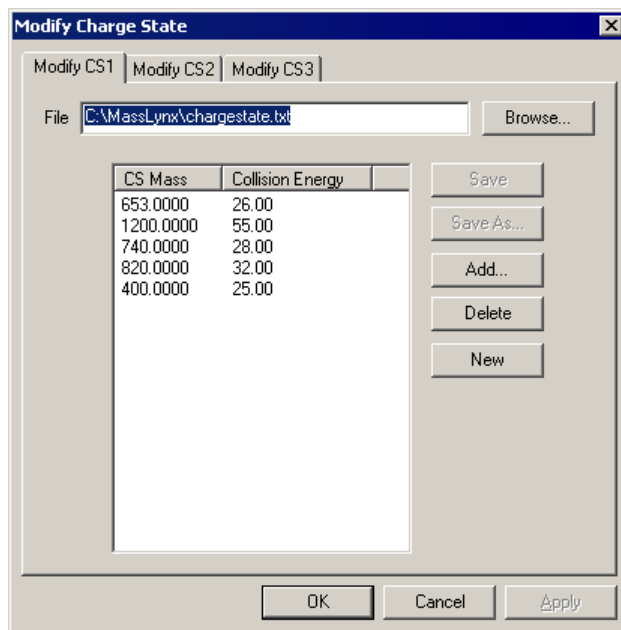


Figure 8-15 The Modify Charge State Dialog

8.4.8 Exclude Masses

The Exclude Mass option (Figure 8-16) allows predefined masses to be excluded from selection for MS/MS analysis. The Exclude Mass window is now defined on this page and a default value of at least 1000mDa should be selected for peptide work, otherwise C13 isotopes will still be selected for fragmentation.

The Range option allows a list of m values separated by commas, or a mass range to be entered into the location in the Exclude Masses window. This is suitable for a small number of components.

For a larger number of components and particularly for components or impurities not permanently present the File option can be used to specify m values and an associated retention time. The window around the retention time is specified in the MS to MSMS section above).

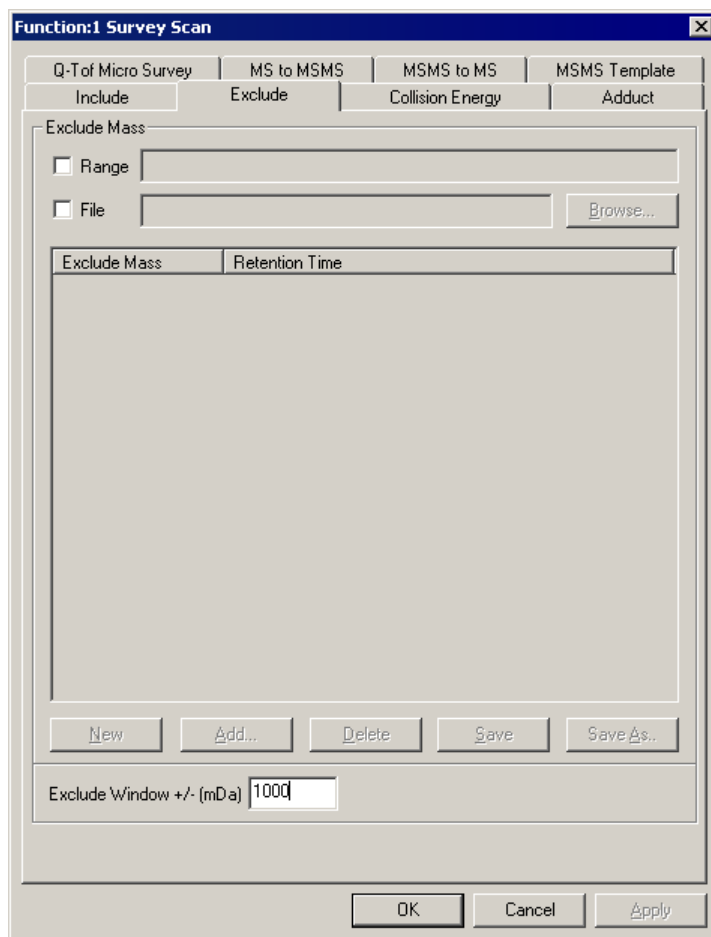


Figure 8-16 Survey Scan Dialog: Exclude Masses Page

Chapter 9


Parent Ion Discovery

Before starting a Parent Ion Discovery Experiment refer to Chapter 8, Data Acquisition. You will also need to refer to the following chapters of the **MassLynx User's Guide**.

- Chapter 4 - Sample Lists
- Chapter 6 - Chromatogram
- Chapter 7 - Spectrum

9.1 Creating a Parent Ion Discovery Experiment

In the MassLynx sample list use the cursor to select the MS file.

1. Select **MS Method** from the Instrument Shortcut bar
2. Select **File > New.** to create a new template with no predefined functions.
3. In the drop down menu select **New**
4. Select  **Parent** to create a parent ion discovery experiment. This invokes the **Parent Scan Dialog** (Figure 9-1). This consists of 10 tabbed pages.

9.1.1 Parent Survey Window

Function:1 Parent Scan

MS to MSMS | MSMS to MS | Include | Exclude | Adduct | Collision Energy

Parent Survey | Product Ions | Neutral Loss | MSMS Template

Masses (m/z)

Start Mass: 100

End Mass: 1500

Method

Ionization Mode: ES+

Data: Continuum

Times (Mins)

Start Time: 0

End Time: 60

Scan Duration (secs)

Scan Time: 0.5

Inter-Scan Delay: 0.1

Cone Voltage

☐ Use Tune Page

Cone Voltage (V): 30

Collision Energy

Low CE (V): 30

High CE (V): 30

Reference Scan

Scans to average: 1

Frequency (s): 0

Cone Voltage (V): 0

High/Low CE Switching

☒ Enable Automatic High/Low Collision Energy Switching

OK Cancel Apply

Figure 9-1 Parent Scan Dialog: Parent Survey Page

- Masses (m/z)** Specifies the start and end m/z range of the low and high-energy survey spectra
- Time (Mins)** Specifies the start and end times in minutes for the experiment.
- Cone Voltage** Allows the user to select either the cone voltage currently in use on the active tune page or to specify a user defined value by entering a suitable value into the text box

Method	Specifies the ionization mode and the mode of data collection. It is essential that centroid acquisition is selected as this technique relies on mass difference in the neutral loss experiment and absolute mass for the product ion experiment.
Scan Duration	Specifies the integration time for the low and high-energy survey scans. Typical values are 1 second for the scan and 0.1 second for the inter-scan delay.
Collision Energy	<p>These are the collision energies used during the low and high energy survey spectra. Typically the low energy is performed at 10eV or lower, whilst the high-energy data acquisition is performed at the appropriate value to produce the desired fragmentation.</p> <p>If a LockSpray interface is configured the Reference Scan options become available.</p>
Scans to Average	The parameter determines the number of reference scans to be averaged to provide a correction factor for exact mass measurements. Values between 5 and 10 are generally suitable.
Frequency	determines how often the reference will be selected and a value of 10 seconds is typical.
Cone Voltage	allows the use to select a suitable cone voltage for the reference compound being used.

9.1.2 Product Ions

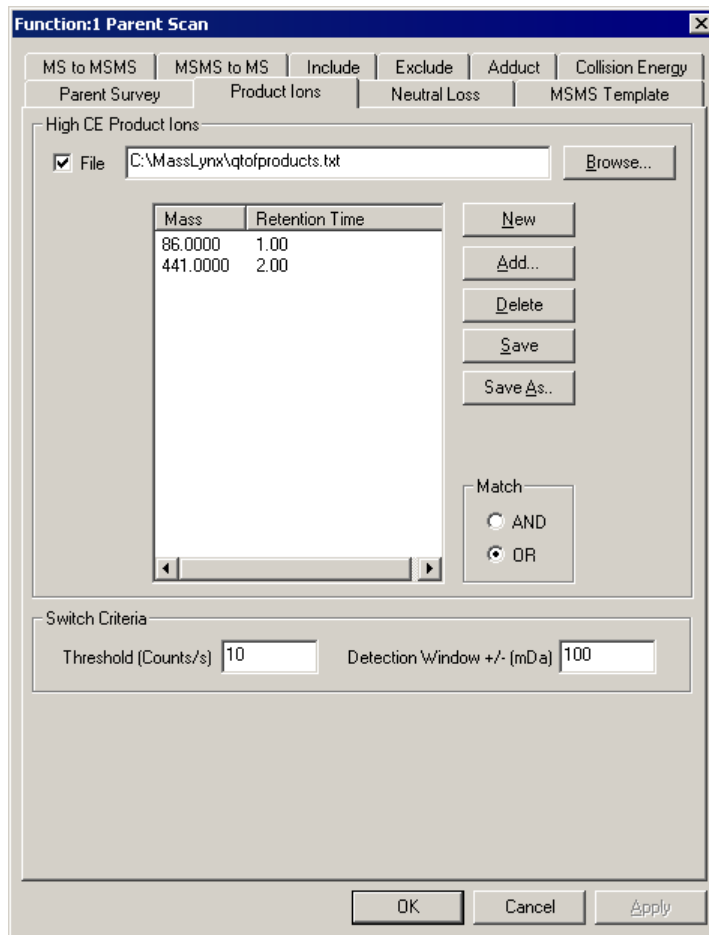


Figure 9-2 Example of a Screen Capture

If the parent ion discovery is to be triggered by the presence of product ions, the masses of these ions should be entered in the above list. Once created the list should be saved under a suitable name.

Switch Criteria

Any ion observed in the high-energy survey spectrum at the m specified above, +/- the detection window must exceed the Threshold (Counts/s) to be considered as the chosen product ion.

Detection Window +/- (mDa) is the mass window for detection placed around the product ion. Keeping this detection window low can increase the specificity of the experiment. Recommended value is 25 mDa. This value can be lowered to 10mDa in the case of an instrument operating with DXC or where a LockSpray interface is being used to apply a lock mass correction.

If more than one product ion is entered in the list then the user must specify the Match as either requiring the presence of one/or more product ions.

9.1.3 Neutral Loss

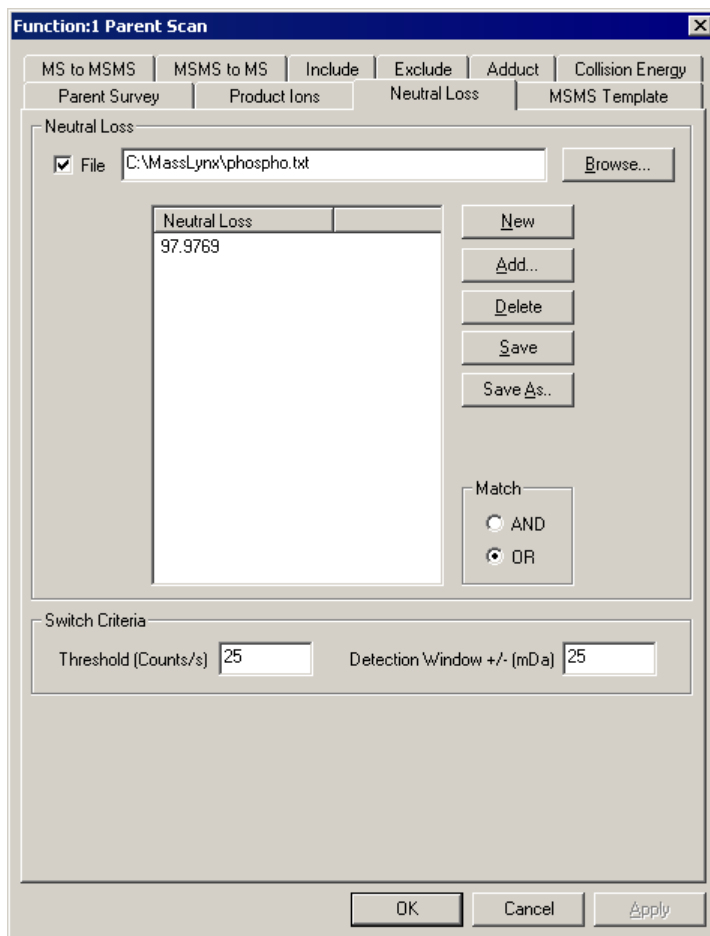


Figure 9-3 Parent Scan Dialog: Neutral Loss Page

If the parent ion discovery is to be triggered by the detection of a neutral loss then the mass of the loss should be entered into the above list. The value of 97.9769, as in the above example, is the loss of H_3PO_4 , which is commonly observed from phosphopeptides. If multiply charged parent ions are to be considered the software will automatically adjust the mass of the neutral loss accordingly.

9.1.4 MSMS Template

Function:1 Parent Scan

MS to MSMS | MSMS to MS | Include | Exclude | Adduct | Collision Energy
 Parent Survey | Product Ions | Neutral Loss | MSMS Template

Instrument Parameter File
 File: c:\masslynx\default.pro\acqddb\default.ipr [Browse]

Mass (m/z)
 Start: 50
 End: 1500

Method
 Ionization Mode: ES+
 Data: Continuum

Cone Voltage
☒ Use Tune Page
 Cone Voltage (V): 30

Scan Duration (secs)
 Scan Time: 1
 Inter-Scan Delay: 0.1

Collision Energy
☐ Use Tune Page
 Collision Energy (V): 30

[OK] [Cancel] [Apply]

Figure 9-4 Parent Scan Dialog: MSMS Template

This window is used to set-up the parameters for the MSMS data acquisition.

It is important that the correct Q-ToF tune file (.ipr file) is selected in the instrument parameter window

Mass (m)	Specifies the m range over which MS/MS data will be acquired. Currently it is vital that both the parent survey and MS/MS template interrogate data over the same m range during the experiment.
Cone Voltage	Allows the user to select the cone voltage currently in use on the active tune page or to specify a value by entering it into the text box.
Collision Energy	Will define a default collision energy that will be used in the MS/MS experiment if, for any reason, a value cannot be selected from the parameters defined later under the collision energy settings

9.1.5 MS to MSMS

MS to MSMS Switch Criteria

The decision can be made using **TIC** or individual ion **Intensity**. Since multiple precursors may be selected for concurrent MS/MS experiments it is usual to set this to **Intensity**. It is useful to acquire some MS data to be able to measure the typical background level and the peak intensities from the sample peaks so that the detection Threshold can be set appropriately.

Within the mass range of the MS survey scan any peaks that rise above the Threshold (in ion counts/sec) will be identified as candidate precursor ions for selection.

To prevent the selection of the P+1 (¹³C) or P+2 (¹³C₂) isotope peaks as additional candidate precursors a **Detection Window** is set. The window is an absolute mass value, i.e:- a setting of 1Da will extract a 1Da range of the m scale and identify only the most intense peak contained within that specific m range. A suitable value for 1+ and 2+ ions would be 2Da.

Number of Components defines the number of precursor masses that can be selected for concurrent MS/MS data acquisition from one MS survey scan. Once the candidate precursor are identified the N most intense will be used for MS/MS. The maximum that can be specified is eight, but often a more suitable number for LC peaks about 10-20 sec wide at half height is 3 - 4.

The detection criteria set here may be subsequently modified by charge state considerations (see later)

Retention Time Window (seconds) This parameter is only active when an include or exclude file has been generated and is selected (see later). If an include or exclude file is utilized where a retention time is included (e.g. 18 mins), then this window is a time value in seconds which is added to the retention time during which the peak must elute to be selected for fragmentation. For example, a value of 30secs would mean the peak expected to elute at 18 mins must elute between 18.00 and 18.50 mins).

Charge State

These settings define the parameters used for the charge state recognition algorithm.

Tolerance Window. The window used when measuring between adjacent peaks (isotopes). This should be a tight tolerance, and as a default set to 0.1Da.

Mass Extraction Window. This value in Da is the mass range around the peak identified that is extracted to perform the CS recognition on. For peptide work a value of 2Da is recommended as a default.

(%) Intensity of Main Peak. This threshold value is utilized to remove low intensity noise peaks from being identified by the CS algorithm. A default value of 2% is recommended.

Precursor Selection

Everything. No user-predefined masses are used and all peaks will be selected automatically according to the criteria defined in the Automatic Function Switching parameters.

Included Masses Only will perform MS/MS solely on those masses appearing in the Include Masses list, or file, (see later section page 7) provided that their intensities rise above the predefined threshold set in the MS to MS/MS switch criteria.

Included Masses and Everything Else allows masses in the Include Masses list to take precedence. These will be selected even if there are more intense peaks that would have been selected automatically.

Detected Precursor Inclusion

Auto Exclude will exclude a m value and its associated detection window (which is now defined under the exclude tab) from subsequent selection.

Always Include allows the immediate selection of a previously selected m value.

Include After Time excludes a previously selected m value and its associated detection window, (which is now defined under the exclude tab) from subsequent selection for the period specified. This prevents multiple selection of the same precursor in the event of tailing from strong LC peaks, but allows the fragmentation of near isobaric components eluting at other retention times. The maximum value is 200 secs.

Data

Discard uninteresting survey scans stores only those survey scans which triggered an MS/MS experiment to disk. Depending on the exact nature of the experiment, discarding these scans may significantly reduce the data file size.

Note: *It is not recommended to use this function with the neutral loss or product ion experiment.*

9.1.6 MSMS to MS

Function:1 Parent Scan

Parent Survey | Product Ions | Neutral Loss | MSMS Template

MS to MSMS | **MSMS to MS** | Include | Exclude | Adduct | Collision Energy

MSMS to MS Switch Criteria

☐ TIC falling below threshold

☐ TIC rising above threshold

☒ Intensity falling below threshold

Threshold (Counts/s) 3

MSMS to MS Switch Method

☐ Default

☒ After time (s) 15

Parent Switching

☒ Absence of Neutral Loss

OK Cancel Apply

Figure 9-5 Parent Scan Dialog: MSMS to MS Page

If the **MSMS to MS Switch Method** is set to **Default** then the mode change will be triggered by one of the mutually exclusive options selected in the **MSMS to MS Switch Criteria**. The Intensity falling below threshold option is measured on the base peak in the MS/MS spectrum. A value of 2-3 counts/sec is recommended.

If the **MS to MS Switch Method** is set to **After time** then the mode change will be as for the default setting, except that if these criteria are not met then a maximum time can be set. (For example in the above case the switch would occur after 15 seconds assuming the intensity of the base peak in the MS/MS spectrum had not fallen below 3 counts/sec).

If **Absence of Neutral Loss** is selected the instrument will switch back to the MS mode if the presence of the tentatively assigned neutral loss cannot be identified in the MS/MS mode. This is highly recommended for phosphopeptide neutral loss experiments.

9.1.7 Include

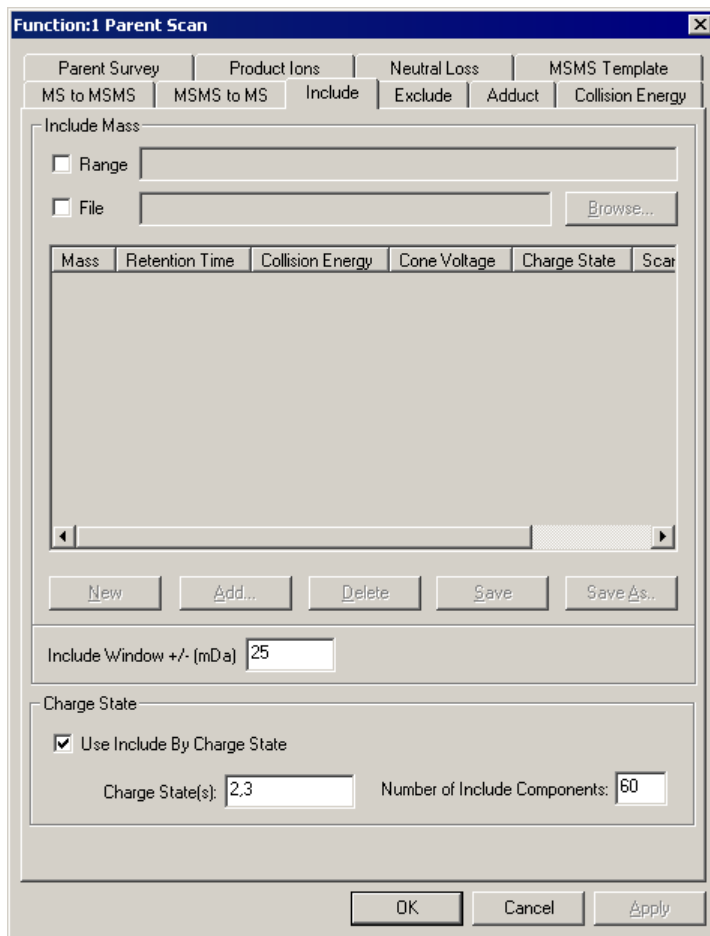


Figure 9-6 Parent Scan Dialog: Include Page

The **Include** option allows a number of predefined m/z values to be included as targeted precursor ions, for switching from MS to MS/MS mode. There are two distinct modes of operation, when operating with an include list and these are defined in the MS to MSMS criteria (see above).

Included masses only will restrict MS/MS data acquisition only to peaks having the same m/z as those defined in the Include Masses List, +/- include window in mDa, provided that these peaks rise above the threshold.

Included masses and Everything else mode will preferentially perform MS/MS on masses in the Include Masses List, even though more intense masses may have been automatically detected, provided that the included masses are above the detection threshold.

The **Listed** option allows a list of m/z values separated by commas to be entered directly into the window or an m/z range e.g. 500_1000. MS/MS experiments on these masses will be carried out using the parameters defined elsewhere.

In contrast, the **File** option allows included masses to be read from a file and specific parameters (**Retention Time, Collision Energy, Cone Voltage, Charge State**) may be associated with these m/z values to override those parameter set elsewhere. The files may be created or modified using the **New, Add, Delete, Save As** and **Save** buttons as necessary.

Charge State

The **Include By Charge State** option allows the selection of only those charge states specified to be selected for MS/MS analysis. This may be useful, for example, where singly charged ions are usually from the background and it is desirable to preclude them from analysis.

The required charge states are specified and separated by commas as shown. It is possible to enter values greater than 4, but at present CE control is available only up to and including charge state 4.

The **Number of Include Components** option window defines the maximum number of components above the detection threshold from which the candidate precursors will be selected. For low level studies where singly charged ions are not to be selected for MS/MS it is clear that there may be many candidates above the threshold, among which the multiply charged ions may be the minor components. The range for this parameter is 1 - 100 and a value of 60 to 80 would be a reasonable default value. It should be noted that raising this value will slow the acquisition process.

9.1.8 Exclude

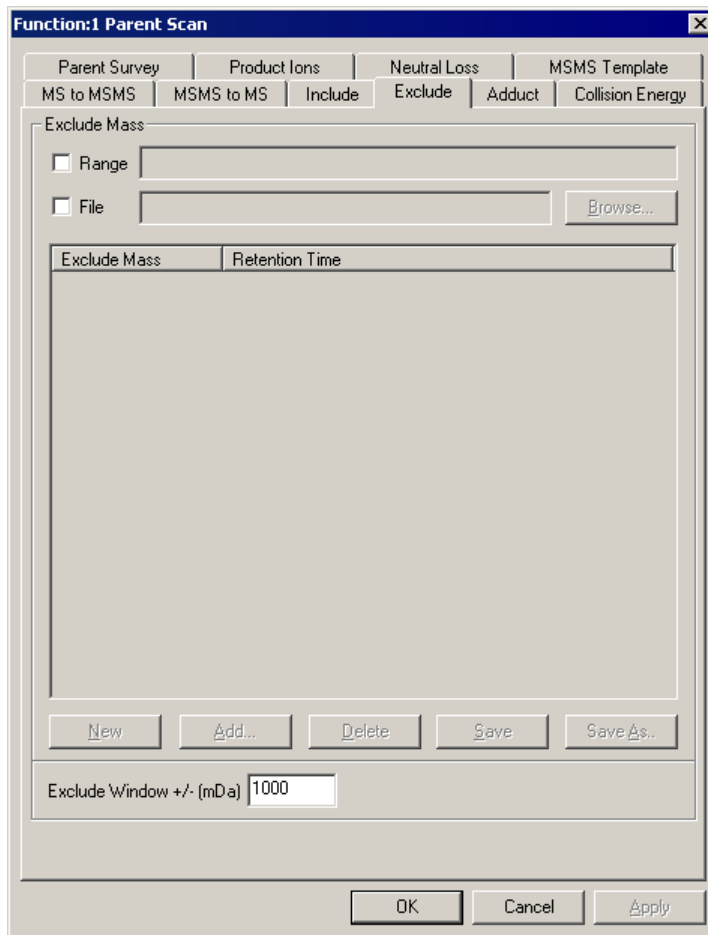


Figure 9-7 Parent Scan Dialog: Exclude Dialog

The **Exclude Mass** window allows predefined masses to be excluded from selection for MS/MS analyses. The **Exclude Mass** window is now defined on this page and a default value of at least 1000mDa should be selected for peptide work, otherwise C12 and C13 isotopes will still be selected for fragmentation as separate peaks.

An exclude **Range** e.g. 50_250 can be entered to prevent switching on solvent or background ions which are known not to be of interest.

Alternatively, the **Range** option allows a list of m/z values separated by commas to be entered into the location in the **Exclude Mass** window. This is suitable for a small number of components.

For a larger number of components and particularly for components or impurities not permanently present the **File** option can be used to specify m/z values and an associated retention time. The window around the retention time is specified in the MS to MSMS section above.

9.1.9 Collision Energy

Function:1 Parent Scan

Parent Survey | Product Ions | Neutral Loss | MSMS Template
MS to MSMS | MSMS to MS | Include | Exclude | Adduct | Collision Energy

Default Collision Energy

☐ Use Default Collision Energy

Collision Energy Profile

☒ Use Collision Energy Profile

File: C:\MassLynx\ce_profile.txt Browse... Modify...

Charge State Recognition

☐ Use Charge State Recognition

Maximum Number Of Charge States: 3

CS1 File: C:\MassLynx\chargestate.txt Browse...
CS2 File: C:\MassLynx\chargestate2.txt Browse...
CS3 File: C:\MassLynx\chargestate3.txt Browse...
CS4 File: Browse...
Modify...

OK Cancel Apply

Figure 9-8 Parent Scan Dialog: Collision Energy Page

Use Default Collision Energy

When selected all MS/MS data are acquired using the default collision energy specified in the MSMS template window.

Use Collision Energy Profile

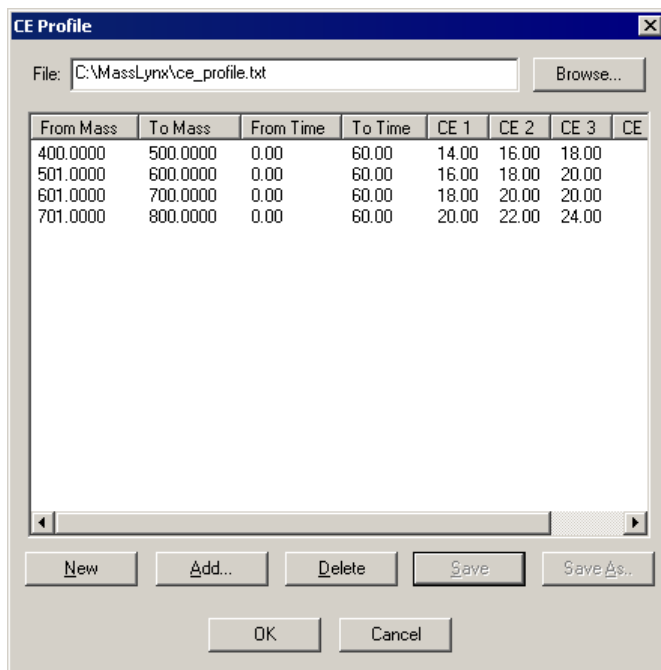


Figure 9-9 The CE Profile Dialog

The collision energy profile option allows a range of different collision energies to be used depending on the m/z value of the selected precursor. In addition, from 1 up to 5 collision energies may be used for each precursor m/z to attempt to produce optimal fragmentation. Using the **Browse** button, previously created files with collision energy profiles may be selected.

Using the **Modify** button new profiles can be created using the **Modify CE Profile** window and its associated **Add**, **Delete**, **New**, **Save** and **Save As** options.

During the MS/MS experiment data are collected with the collision energy cycling through CE1, CE2....CEN on a per spectrum basis.

The system will not interpolate values. MS/MS on any value of m/z not included in this table will be performed using the default value set in the MS/MS template. Below is a typical collision energy profile covering the range 400 to 1000 m/z .

Use Charge State Recognition

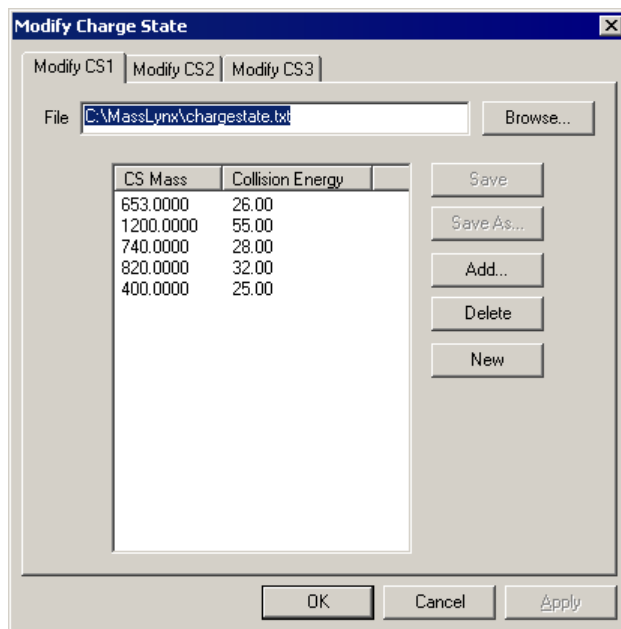


Figure 9-10 The Modify Charge State Dialog

This mode will allow the collision energy to be set according to the charge state and m/z of the precursor as determined from the charge state recognition algorithm. Collision energy profiles for charge states from 1, 2, 3 and 4 may be created and the text files selected using the **Browse** button. In this mode collision energy profiling by cycling around a number of different energies is not required. The system will linearly interpolate between values in the table.

Using the **Modify** button new profiles can be created using the **Modify Charge State** window and its associated **Add**, **Delete**, **New**, **Save** and **Save As** options. The examples shown are a set of suggested collision energy profiles for peptides or protein digest applications when selecting 2+ ions for MS/MS analysis.

Chapter 10

Maintenance and Fault Finding

Cleanliness and care are of the utmost importance whenever internal assemblies are removed from the instrument.

- Always prepare a clear clean area in which to work.
- Make sure that any tools or spare parts that may be required are close at hand.
- Obtain some small containers in which screws, washers, spacers etc. can be stored.
- Use tweezers and pliers whenever possible.
- If nylon or cotton gloves are used take care not to leave fibres in sensitive areas.
- Avoid touching sensitive parts with fingers.
- Do not use rubber gloves.
- Before reassembling and replacing dismantled components, inspect O rings and other vacuum seals for damage. Replace with new if in doubt.

Should a fault occur soon after a particular part of the system has been repaired or otherwise disturbed, it is advisable first of all to ensure that this part has been correctly refitted and / or adjusted and that adjacent components have not been inadvertently disturbed.



Caution: Many of the procedures described in this chapter involve the removal of possibly toxic contaminating deposits using flammable or caustic agents. Personnel performing these operations should be aware of the inherent risks, and should take the necessary precautions.

This Chapter describes maintenance for the following integral parts of the Q-ToF *micro*:

- Section 10.2, The Source
- Section 10.3, The Electrospray Probe
- Section 10.4, The APcI Probe
- Section 10.5, The Analyzer

In addition further advice and tips are given in the following sections:

- Section 10.6, Fault Finding
- Section 10.7, Cleaning Materials
- Section 10.8, Preventive Maintenance Check List

10.1 Electronics Maintenance

10.1.1 Cooling Fans and Air Filters

Always ensure that none of the cooling fans is obstructed. It is essential that the fan filter is checked and cleaned at regular intervals, and replaced if there is any doubt about its effectiveness.

10.1.2 The Vacuum System

The performance of the mass spectrometer will be severely impaired by the lack of a good vacuum in the ion transfer (hexapole) region or the analyzer.

- An excessive analyzer pressure results in a general loss in performance indicated by a loss of resolution and an increase in the background noise.
- As the vacuum deteriorates, the vacuum becomes insufficient to maintain the instrument in the operate mode.

Before suspecting a leak, the following points should be noted:

- The turbomolecular pumps will not operate if the rotary pump has failed.
- If the rotary pump is not maintained, the oil may become so contaminated that optimum pumping speed is no longer possible. Initially, gas ballasting may clean the oil. If the oil in the rotary pump has become discolored, then it should be changed according to the pump manufacturer's maintenance manual.
- The turbomolecular pumps switch off if an over temperature is detected. This could be due to poor backing vacuum, failure of the water supply or a leak in the source or analyzer.

Vacuum Leaks

If a leak is suspected, the following basic points may help to locate it:

- Leaks very rarely develop on an instrument that has been fully operational. Suspect components that have recently been disturbed.

Leaks on flanges can usually be cured by further tightening of the flange bolts or by replacing the seal.

- All seals are made using O rings. When refitting flanges pay attention to the condition of O rings. Any that are cut or marked may cause a leak. The O rings should be clean and free from foreign matter.

A hair across an O ring is sufficient to prevent the instrument pumping down.

- Source components that operate at, or slightly above, atmospheric pressure are not susceptible to vacuum leaks.

In the unlikely event of a leak on a feedthrough, then the unit should be replaced or returned to Micromass for repair.

Pirani Gauge

The Pirani gauge head does not require routine maintenance.

Active Inverted Magnetron Gauge

In particular, the quadrupole analyzer gauge requires regular maintenance. For information on cleaning the active inverted magnetron (Penning) gauge, refer to the Edwards literature supplied with the instrument.

Gas Ballasting

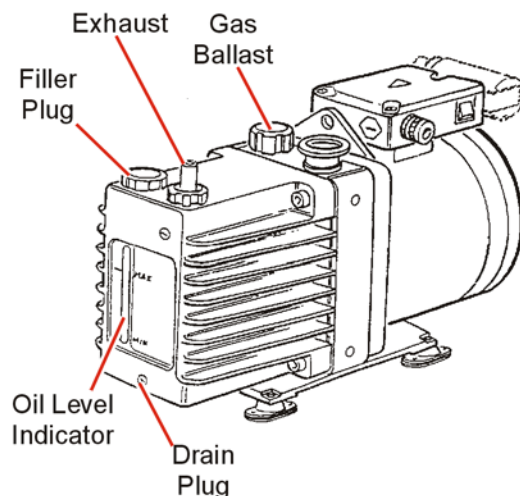


Figure 10-1 Rotary Pump

Gas ballasting serves two important purposes:

- When rotary pumps are used to pump away solvent vapors, the solvent vapor can become dissolved in the pump oil causing an increase in backing line pressure. Gas ballasting is a method of purging the oil to remove dissolved contaminants.
- Oil mist expelled from the rotary pump exhaust is trapped in the oil mist filter. This oil is returned to the rotary pump during gas ballasting.

Gas ballasting should be performed routinely on a weekly basis for 30 minutes. If the source is used in the APcI or megafLOW electrospray modes, more frequent gas ballasting is recommended.

Gas ballasting is performed on the E2M28 pump by rotating the gas ballast valve 5 to 6 turns in a counterclockwise direction.

It is normal for the rotary pump to make more noise when the gas ballast valve is open.



Caution: Failure to gas ballast the rotary pump frequently leads to shortened oil lifetime which in turn may shorten rotary pump lifetime.

Oil Mist Filter

The E2M28 rotary pump is fitted with an Edwards EMF20 oil mist filter which traps oil vapor from the rotary pump exhaust. The trapped oil is then returned to the rotary pump during routine gas ballasting. The oil mist filter contains two elements; the odor element need not be changed, but the mist element must be changed every time the rotary pump oil is changed.

- To change the element follow the instructions in the Edwards manual.

Rotary Pump Oil

The oil in the rotary pump should be maintained at the correct level at all times. Check the oil level at weekly intervals, topping up if necessary.

It is important to monitor the condition of the oil regularly. Replace the oil when it has changed to a noticeable reddish brown color, or routinely at 4 month intervals (3000 hours operation). At the same time, replace the oil mist filter's mist element (see above).

Change the oil in the rotary pump as follows:

1. Gas ballast lightly for 30 to 60 minutes.
2. Vent and shut down the instrument as described in Routine Procedures.
3. It will be found easier to drain the oil while the pump is still warm.
4. Drain the oil through the drain hole situated near the oil level sight glass.
5. Flush the pump, then replace the drain plug and refill the pump with the correct grade oil to the correct level.
6. Gas ballast lightly for 30 to 60 minutes.

For further servicing information refer to the manufacturer's manual.

Foreline Trap

This is used to protect against the chance of pump oil backstreaming into the collision gas solenoid and/or the turbo pumps. The activated alumina should be changed according to the manufacturer's instructions.

10.2 The Source

The Z-spray source is a robust assembly requiring little maintenance. The source consists of three basic parts:

- The probe adjustment flange.
- The glass tube.
- The source flange assembly.

The probe adjustment flange and the glass tube can be readily removed, without venting the instrument, to gain access to the source block and sample cone. This allows the following operations to be performed:

- Wiping the sample cone.
- Removing the sample cone.
- Fitting or removing the APcI corona discharge pin.
- Fitting or removing the exhaust liner and cleanable baffle.
- Fitting or removing the nanoflow electrospray interface.
- Enabling or disabling the purge gas.

The sample cone may be cleaned *in situ*, by gentle wiping with a cotton swab or lint tissue soaked with 50:50 acetonitrile:water. More thorough cleaning of the sample cone may be achieved by removing it from the source. This may also be done without venting the instrument, by closing the isolation valve located on the ion block. Less frequently it may be necessary to clean the ion block, the extraction cone and the hexapole lens, in which case the instrument must be vented. This should only be done when the problem is not rectified by cleaning the sample cone or when charging effects are apparent.

Charging is evidenced by a noticeable progressive drop in signal intensity, often resulting in a complete loss of signal. Switching the instrument out of and back into operate causes the beam momentarily to return.

The hexapole transfer lens should not require frequent cleaning. If it is suspected that the lens does need cleaning it may be withdrawn from the front of the instrument after removing the ion block support.



Caution: *Cleaning the various parts of the source requires the use of solvents and chemicals which may be flammable and hazardous to health. The user should take all necessary precautions.*

10.2.1 Cleaning the Sample Cone in Situ

This may be necessary due to lack of sensitivity or fluctuating peak intensity, or if deposited material is visible on the outside of the sample cone. Proceed as follows:

1. On the MassLynx Shortcut bar launch the tune page.

2. Deselect Operate to put the instrument in standby mode
3. Switch off the LC pumps.
4. Disconnect the liquid flow at the rear of the probe.
5. Set **Source Block Temp** and either **APCI Probe Temp** or **Desolvation Temp** to 20°C to switch off the heaters.



Caution: Removal of the APCI probe or desolvation nozzle when hot may cause burns.



Attention: Removal of the APCI probe when hot will shorten the probe heater's life.

The cooling time will be significantly shortened if the API gases are left flowing.

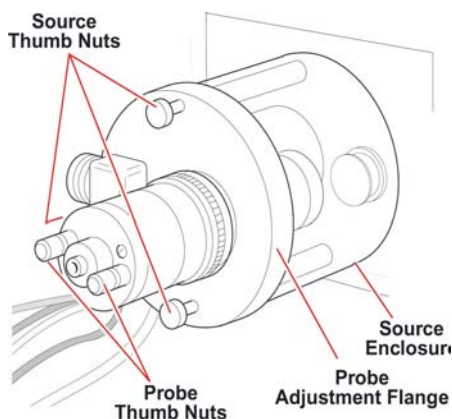


Figure 10-2 Removing The Source Enclosure

When **APCI Probe Temp** or **Desolvation Temp** has cooled below 100°C:

6. Deselect API Gas to switch off the nitrogen supply.
7. Disconnect both gas lines from the front panel by undoing the knurled nuts.
8. Disconnect both electrical connections by pulling back on the plug sleeves to release the plugs from the sockets on the front panel.
9. Undo the two knurled thumb nuts that retain the probe and withdraw it from the source. Place it carefully to one side.

10. Undo the three thumb screws and withdraw the probe adjustment flange and glass tube. Place the glass tube, end on, on a flat surface and place the probe adjustment flange on top of the glass tube.



Caution: When the source enclosure has been removed the source block is exposed. Ensure that the source block heater has cooled before proceeding.

11. If fitted, remove the APcI corona discharge pin.

The sample cone is now accessible.

12. Using a suitable flat blade screwdriver rotate the isolation valve by 90° into its fully anticlockwise position.

A small improvement in the analyzer vacuum may be observed as a result of this operation.

The isolation valve (Figure 10-3) is closed when the slot is perpendicular to the direction of flow.

13. Carefully wipe the sample cone with a cotton swab or lint free tissue soaked in 50:50 acetonitrile:water or 50:50 methanol:water.

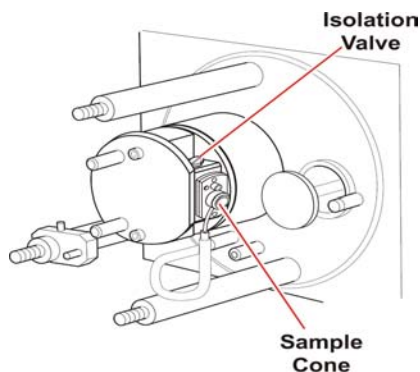


Figure 10-3 The Source, With Housing Removed



Attention: Do not attempt to remove any obstruction by poking. This may result in damage

14. Dry the cone using nitrogen.

If the sample cone is still not clean, or if the aperture is partially blocked, proceed to the following section. Otherwise, when the cone is clean and dry:

15. Open the isolation valve.
16. Replace all removed components, following in reverse order the removal procedures.

10.2.2 Removing and Cleaning the Sample Cone



Attention: *The sample cone is a delicate and expensive component and should be handled with extreme care.*

It is not necessary to vent the instrument to remove the sample cone. The source block incorporates an isolation valve for this purpose. To remove the sample cone proceed as follows:

1. Follow the procedure in the previous section, to gain access to the sample cone.
2. Using a suitable flat blade screwdriver rotate the valve by 90° into its fully anticlockwise position.

A small improvement in the analyzer vacuum may be observed as a result of this operation.

The isolation valve is in the closed position when the slot is perpendicular to the direction of flow.

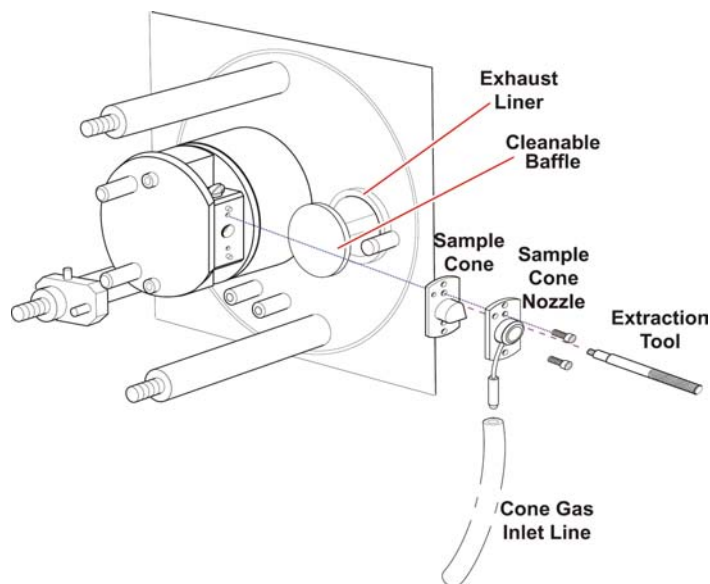


Figure 10-4 Removing the Cone

3. Disconnect the cone gas inlet line (if fitted).
4. Take the sample cone extraction tool supplied in the source spares kit and screw it to the flange of the sample cone.
5. Remove the two sample cone retaining screws using a 1.5mm Allen key and withdraw the sample cone and cone gas nozzle (if fitted) from the ion block.
6. Remove the extraction tool, and separate the sample cone from the cone gas nozzle. Place both components in an ultrasonic bath containing 40:40:10 acetonitrile:water:formic acid or 40:40:10 methanol:water:formic acid. Rinse and sonicate with 50:50 acetonitrile:water or methanol:water.
7. Dry the cone and nozzle using nitrogen.

To minimize down time fit a spare sample cone, obtainable from Micromass, at this stage.

If material has built up on the exhaust liner and cleanable baffle:

1. Remove the cleanable baffle and the exhaust liner.



Attention: Do not attempt to remove the baffle without first removing the sample cone,

2. Clean these components, or obtain replacements.
3. Fit the cleaned (or the replacement) exhaust liner and cleanable baffle to the ion block.

Note: Refitting the sample cone is a reversal of the removal procedure.

10.2.3 Removing and Cleaning the Source Block and Extraction Cone

1. From Tune Page menu bar select **Vacuum > Vent**.

The rotary pump and the turbomolecular pumps switch off. The turbomolecular pumps are allowed to run down to 50% speed after which a vent valve automatically admits dry nitrogen.

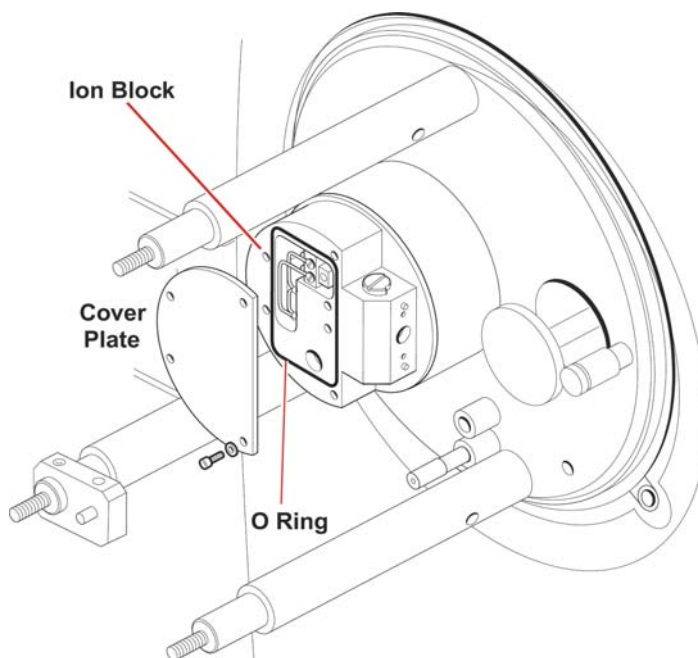


Figure 10-5 Source Block and Extraction Cone Removed

2. Remove the source enclosure and the sample cone as described in the previous section.

When the Instrument has Vented:

3. Remove the two screws which secure the ion block and remove the ion block heater and the ion block.

4. Separate the extraction cone and the PTFE insulating ring from the ion block.
5. Remove the plug and the PTFE sealing washer.
6. Remove the sample cone as described above.
7. Leaving the valve stem in place, immerse the ion block in an ultrasonic bath containing 50:50 acetonitrile:water or 50:50 methanol:water, followed by 100% methanol.
8. Clean the sample cone and the extraction cone using in turn:
 - concentrated formic acid.
 - 50:50 acetonitrile:water or 50:50 methanol:water.
 - 100% methanol.



Caution: Strong acid causes burns. Carry out this procedure in a fume cupboard using protective equipment.

9. Dry all components using a flow of nitrogen, or place them in a warm oven.

10.2.4 Removing and Cleaning the RF Lens Assembly

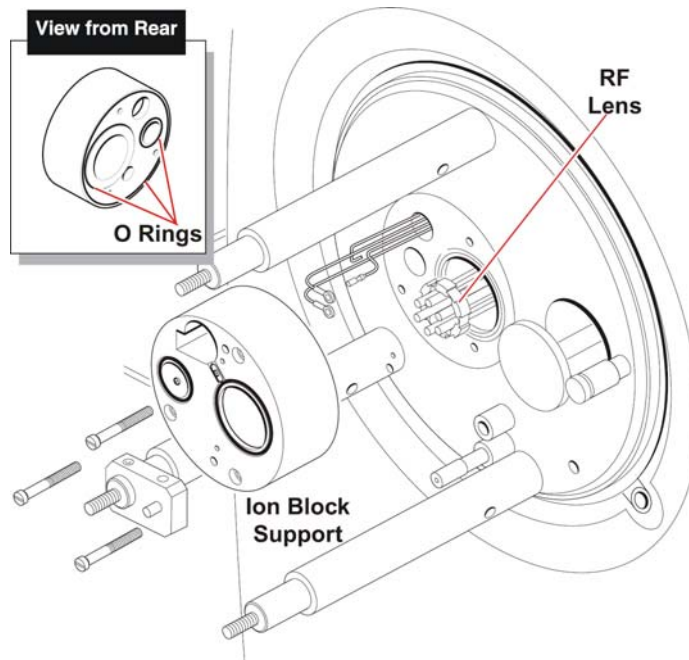


Figure 10-6 Removing The RF Lens Assembly

To remove the RF hexapole transfer lens assembly, proceed as follows:

1. Remove the ion block, as described above.
2. Remove the three screws retaining the ion block support and carefully withdraw it, together with the support liner and O rings, from the pumping block.
3. Using a lint free tissue to gently grasp the hexapole, carefully withdraw it.



Attention: Take care not to scratch the internal bore of the pumping block as the hexapole lens assembly is withdrawn.

To clean the hexapole transfer lens proceed as follows:

4. Immerse the complete assembly in a suitable solvent (100% methanol) and sonicate in an ultrasonic bath.

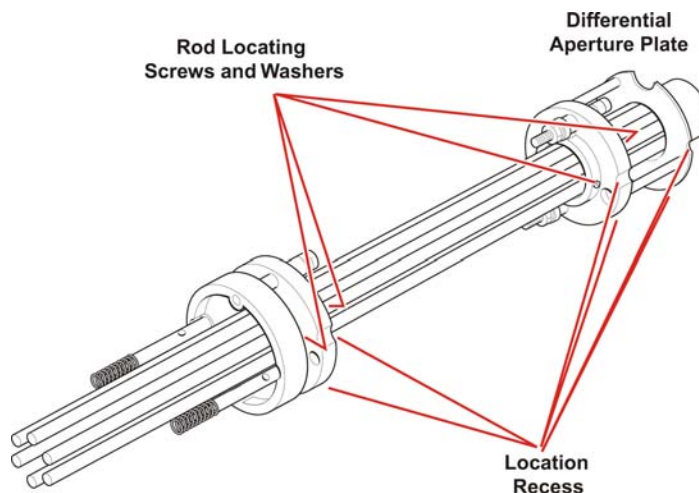


Figure 10-7 The RF Lens Assembly

5. Thoroughly dry the assembly using a flow of nitrogen.

In severe cases:

1. Remove, clean, dry and replace each rod separately (one at a time).
2. Reassemble the assembly with extreme care, checking the assembly against the diagram.

10.2.5 Reassembling and Checking the Source

1. Feed the hexapole transfer lens into the instrument, allowing the recesses in the differential aperture plate to locate onto the two support rails within the analyzer assembly. Ensure that the assembly is pushed fully in.
2. Check the condition of the O rings on the ion block support. Replace them if necessary.
3. Replace the ion block support, pushing it in against the springs of the hexapole assembly.
4. Replace the three retaining screws.
5. Fit the plug and sealing ring to the ion block.
6. Fit the insulating ring and extraction cone.
7. Offer the ion block up to the peek ion block support, locate the two dowels and push firmly.

8. Replace the ion block heater.
9. Replace and firmly tighten the two retaining screws taking care not to over-tighten the screws.
10. On the tune page select **Vacuum > Pump**.
11. Replace the PTFE exhaust liner and cleanable baffle, if removed.
12. Replace the sample cone and, if the nanoflow option is to be used, the cone gas nozzle on the ion block.
13. Reconnect the cone gas supply (nanoflow operation only).
14. Plug the purge and cone gas outlets and fit the APcI corona discharge pin.
15. Fit the source enclosure and the probe adjustment flange.
16. Insert the APcI probe and connect the **Nebuliser Gas** line.
17. Select **Inlet > Gas** and turn on **Nitrogen**. Fully open the *Nebuliser Gas* valve.
18. Set **Desolvation Gas** to read back 400 l/h (monitored on the tune page).
19. Check for gas leaks using soap solution.
20. Reduce **Desolvation Gas** to 150 l/h.
21. Set **Source Block Temp** to 150°C, and **APcI Probe Temp** to 20°C



Attention: The maximum operating temperature for the source heater is **150°C**. Do not set **Source Block Temp** higher than **150°C**.

22. Select **Operate** on the tune page.
23. With Corona set to zero, check that the Cone readback is reading the correct set value.
24. Set Corona to 4.0kV.
25. Check that the Corona readback is 4.0 kV and that the Cone readback is still reading the same set value.
26. Check that all other readbacks on the tune page agree with the set values.

10.2.6 The Corona Discharge Pin

If the corona discharge pin becomes dirty or blunt:

1. Remove it from the source.
2. Clean and sharpen it using 600 grade emery paper.

If the needle becomes bent or otherwise damaged it should be replaced.

10.3 The Electrospray Probe



Attention: *The probe tip is sharp, and may be contaminated with harmful and toxic substances. Always take great care when handling the electrospray probe.*

Indications that maintenance is required to the electrospray probe include:

- An unstable ion beam.

Nebulizing gas may be escaping from the sides of the probe tip.

Ensure that the probe tip O ring is sealing correctly.

The probe tip setting may be incorrect.

Adjust the probe tip setting as described in Electrospray.

The probe tip may be damaged.

Replace the probe tip.

There may be a partial blockage of the sample capillary or the tubing in the solvent flow system.

Clear the blockage or replace the tubing.

- Excessive broadening of chromatogram peaks.

This may be due either to inappropriate chromatography conditions, or to large dead volumes in the transfer capillaries between the LC column or probe connection.

Ensure that all connections at the injector, the column, the splitting device (if used) and the probe are made correctly.

- High LC pump back pressure.

With no column in line and the liquid flow set to 300 μ l/min the back pressure should not exceed 7 bar (100 psi). Pressures in excess of this indicate a blockage in the solvent flow system.

Samples containing particulate matter, or those of high concentrations, are most likely to cause blockages.

Check for blockages at the tube connections and couplings to the injector, the column and, if used, the flow splitter.

Concentrated formic acid can be injected to clear blockages. Rinse thoroughly afterwards.

Blockage of the stainless steel sample capillary may occur if the desolvation heater is left on without liquid flow. This is particularly relevant for samples contained in involatile solvents or high analyte concentrations. To avoid this problem it is good practice to switch off the heater before stopping the liquid flow, and flush the capillary with solvent.

A blocked stainless steel sample capillary can often be cleared by removing it and reconnecting it in the reverse direction, thus flushing out the blockage.

- Gas flow problems

Check all gas connections for leaks using soap solution, or a suitable leak searching agent such as Snoop.

10.3.1 Replacement of the Stainless Steel Sample Capillary

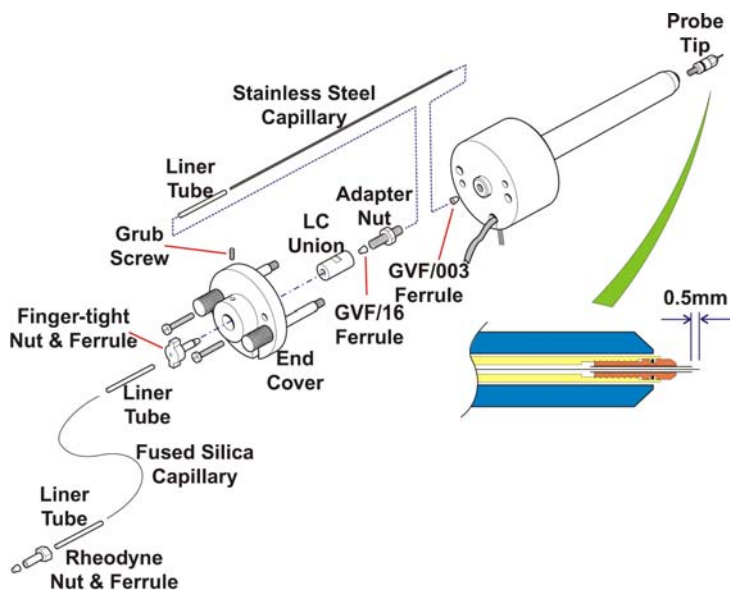


Figure 10-8 Removing the Stainless Steel Capillary from the Source

If the stainless steel sample capillary cannot be cleared, or if it is contaminated or damaged, replace it as follows:

1. Remove the probe from the source.
2. Disconnect the LC line from the probe and remove the finger-tight nut.

3. Loosen the grub screw retaining the LC union.
4. Remove the two probe end cover retaining screws, and remove the probe end cover.
5. Unscrew and remove the probe tip.
6. Remove the LC union and adapter nut. Withdraw and discard the stainless steel sample capillary.
7. Remake the LC connection to the LC union.
8. Sleeve one end of new sample capillary with the PTFE liner tube.
9. Using a GVF/16 ferrule and the adapter nut, connect the sample capillary to the LC union, ensuring that both the liner tube and sample capillary are fully butted into the LC union.
10. Disconnect the LC connection and feed the sample capillary through the probe, ensuring that a 0.3mm graphitized vespel ferrule (GVF/003) is fitted.
11. Using a Rheodyne spanner, gently tighten the adapter nut onto the probe.
12. Replace the probe tip and adjust so that 0.5mm of sample capillary protrudes from the probe tip.
13. Replace the probe end cover and tighten the grub screw to clamp the LC union.

10.4 The APcI Probe

Indications that maintenance to the APcI probe is required include:

- The probe tip assembly becomes contaminated, for example by involatile samples if the probe temperature is too low during operation (300°C).
- The appearance of chromatogram peak broadening or tailing.

Samples that give rise to a good chromatogram peak shape in APcI (for example reserpine and common pesticides) should display peak half widths of the order 0.1 minutes for 10µl loop injections at a flow rate of 1 ml/min. The appearance of significant peak broadening or tailing with these compounds is most likely to be due to a broken fused silica capillary or probe tip heater assembly.

- Low LC pump back pressure.

For 50:50 acetonitrile:water at a flow rate of 1 ml/min, a LC pump back pressure less than 14 bar (200 psi) is indicative of a broken fused silica capillary or a leaking connector.

- High LC pump back pressure.

For 50:50 acetonitrile:water at a flow rate of 1 ml/min, a LC pump back pressure above 35 bar (500 psi) is indicative of a blockage or partial blockage in the fused silica capillary, in a LC connector or in the filter. It is advisable to change the inner filter pad on a regular basis (see “Replacing the Fused Silica Capillary” in the following pages).

- Gas flow problems.

Check all gas connections for leaks using soap solution, or a suitable leak searching agent such as Snoop.

10.4.1 Cleaning the Probe Tip

Remove any visible deposits on the inner wall of the probe heater with a micro-interdental brush (supplied in the spares kit) soaked in methanol:water.

Before starting an analysis:

1. With the probe out of the instrument, connect the nebulizing gas supply line.
2. Select **API Gas** and turn on **Nitrogen**.
3. Allow the gas to flow for several seconds to clear any debris from the heater.
4. Turn off **Nitrogen**.
5. Insert the probe into the source.
6. Select **API Gas** and turn on **Nitrogen**.
7. Raise **APCI Heater** gradually, starting at 100°C and increasing in 50°C intervals to 650°C over a period of 10 minutes.



Attention: Do not set APCI Heater to 650°C immediately as this may damage the probe heater.

This procedure should remove any chemical contamination from the probe tip.

10.4.2 Replacing the Probe Tip Heater

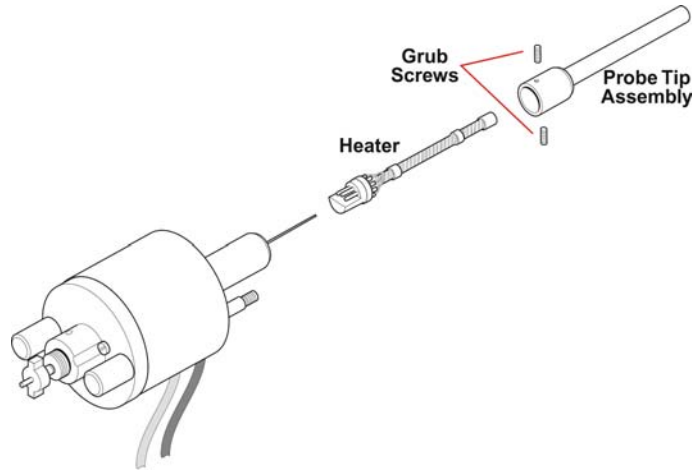


Figure 10-9 The APcl Probe Tip Heater

1. Remove the probe tip assembly by carefully loosening the two grub screws.
2. Disconnect the heater from the probe body by pulling parallel to the axis of the probe.
3. Fit a new heater assembly.
4. Reconnect the probe tip assembly.

10.4.3 Replacing the Fused Silica Capillary

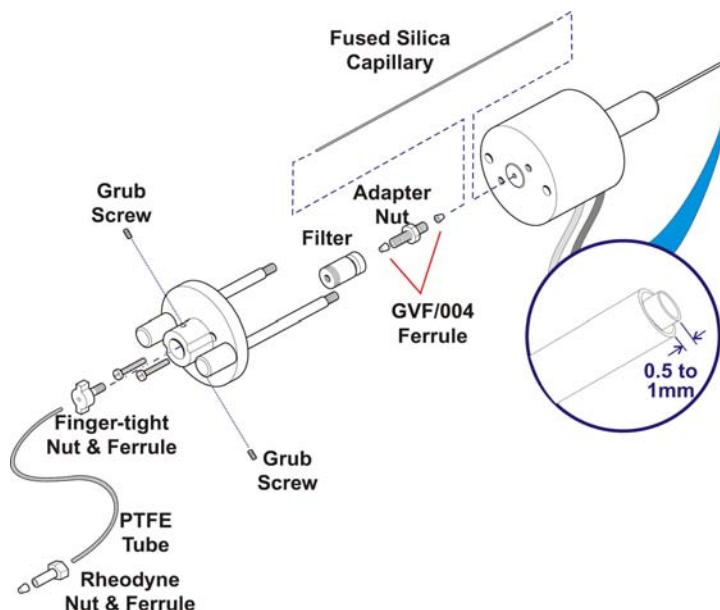


Figure 10-10 Replacing the Fused Silica Capillary on the APcI probe

With the probe removed from the source proceed as follows:

1. Remove the probe tip assembly and the heater, as described in the preceding section.
2. Remove the probe end cover by removing the two screws and the grub screws that retain the LC filter.
3. Loosen the filter from the adapter nut.
4. Unscrew the adapter nut from the probe.
5. Remove and discard the fused silica capillary.
6. Using a ceramic capillary cutter, cut a new length of $300\mu\text{m o.d.} \times 100\mu\text{m i.d.}$ fused silica capillary, about 1 centimetre excess in length.
7. Using a GVF/004 ferrule and the adapter nut, connect the capillary to the filter ensuring that the capillary passes through the ferrule but stops short of the filter.
8. Feed the sample capillary through the probe, ensuring that a 0.4mm graphitized vespel ferrule (GVF/004) is fitted.
9. Using a ceramic capillary cutter, cut the capillary at the nebuliser so that between 0.5 and 1.0mm of capillary is protruding from the nebuliser.

It is important to cut the capillary square. This should be examined using a suitable magnifying glass.

10. Undo the adapter nut from the probe and withdraw the capillary from the probe.
11. Remove 20mm of polyamide coating from the end of the capillary using a flame and clean with a tissue saturated with methanol.
12. Carefully re-feed the sample capillary through the probe ensuring that the graphitized vespel ferrule is still fitted.
13. Using a Rheodyne spanner, gently tighten the adapter nut to the probe.
14. Replace the probe end cover and retaining screws.
15. Using a 1.5mm Allen key, tighten the grub screw in the probe end cover to clamp the filter.
16. Replace the heater and probe tip assembly.

10.5 The Analyzer

The analyzer element of any high performance quadrupole mass spectrometer is, of necessity, a precisely machined and aligned assembly.

Q-ToF micro is fitted with a prefilter assembly designed to protect the main analyzer by absorbing the majority of any contamination. The prefilter is not as mechanically critical as the main rods, as it has only RF applied and is, therefore, not resolving. It does, however, act as a high pass filter, and will reject ions of low mass before they enter the main rods.

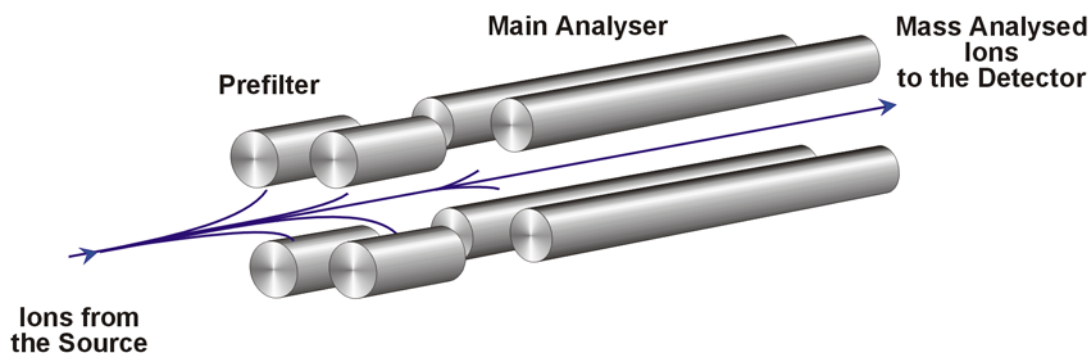


Figure 10-11 The Prefilter and Main Analyzer of the Q-ToF *micro*



Attention: Cleaning the analyzer on site is not a task to be undertaken lightly. It should be done only when there is clear evidence that the analyzer is the source of lack of performance. The main rods of the analyzer are precision ground, as are the ceramics in which they are mounted. Under no circumstances should any of the analyzer rods be removed from the ceramic mountings.

It is unlikely that there will be any means on site for measuring the assembly to the micron level required. If analyzer charging effects cannot be resolved by the techniques below, it is almost certain that the analyzer will have to be returned to Micromass for refurbishment.

10.5.1 Removing the MS1 and Collision Cell Optical Bench Assembly

1. Vent the system as described at the beginning of this chapter.
2. Remove the analyzer housing top plate.

Working inside the MS1 analyzer housing proceed as follows:

3. Disconnect the two push-on connectors which take the analyzer RF through to the RF lens at the source end of the analyzer.
4. Disconnect the three push-on connections bringing RF and DC voltages to the gas cell (analyzer housing right hand side).
5. Disconnect the three push-on connections bringing RF and DC voltages to the hexapole transfer lens following the gas cell (analyzer housing right hand side).
6. Remove the gas line connection to the gas cell where it enters the gas cell.
7. Disconnect the push-on DC (pole bias) connection to the main filter.
8. Disconnect the heavy gauge copper RF/DC main filter supply from the feedthroughs on the left side of the analyzer housing.

All connections to the ion optical bench (main filter assembly/gas cell/RF lens assembly) should now have been removed.

9. At the front and rear of the optical bench, remove the two 4mm Allen socket screws that secure the bench to the mounting flanges.



Attention: Before removing the optics from the housing prepare a clean area where it may be placed in preparation for removing and cleaning the prefilter rods.

10. Holding the optical bench assembly, using two of the clamps that secure the optical elements to the aluminium base plate, lift the rear end upwards to clear the vacuum housing and then withdraw the optical bench from the housing.
11. Once removed replace the vacuum housing lid to prevent the ingress of dust into the housing.

10.5.2 Dismantling and Cleaning the Entrance Prefilter

Dismantling and cleaning of the entrance prefilter is a skilled procedure which should be entrusted only to a Micromass engineer.

10.5.3 Cleaning the MS1 Analyzer Assembly



Attention: Under no circumstances should any of the analyzer rods be removed from the ceramic mountings.

1. Having removed the optical bench assembly as described above, remove the two top clamps that secure the main filter (including pre and post-filter assemblies) to the optical bench.
2. Roll up a narrow strip of absorbent lint-free paper. Pass one end through the gap between two adjacent main filter rods and back through the next gap so that the paper is wrapped one half turn around a rod.
3. Wet the paper with a solvent such as methanol, and move the strip up and down along the analyzer rod.
4. Remove the strip and use dry nitrogen gas from a cylinder (not compressed air) to blow out any dust or particles.



Attention: Tools, carborundum paper or micromesh must not be used to remove contamination from the main filter rods.

It is not anticipated that the gas cell or hexapole transfer lens will require cleaning.

10.5.4 Replacing the MS1 and Gas Cell Optical Bench Assembly

Reassembly is the reverse of the appropriate dismantling procedure. Take extra care to ensure that all electrical connections are correctly and securely made, and that the various mechanical assemblies are accurately aligned within the housing on their locating dowels.

10.5.5 The MS2 Analyzer and Detector Assembly

The orthogonal acceleration cell, reflectron and microchannel plate assemblies are separated from the MS1 and gas cell region by a differential pumping orifice. It is anticipated that this region of the instrument will not require routine maintenance under normal operating conditions.

The microchannel plate system is at risk from failure properly to condition the detector following venting of the system to atmosphere or vacuum failure.



Attention: *It is strongly recommended that assistance is sought from Micromass if maintenance to any of the components within the TOF analyzer housing are thought to be necessary (e.g. excessive noise, spikes, loss of detector gain or abnormal TOF peak shapes).*

10.6 Fault Finding

The majority of faults that occur can be traced to a malfunction of the ion source or inlet system. On systems equipped with more than one source, this can often be confirmed by changing sources to see if the fault “moves” with the source.

Should a fault occur soon after a part of the system has been repaired or otherwise disturbed, it is advisable first of all to ensure that this part has been correctly refitted and adjusted, and that adjacent components have not been inadvertently disturbed.

10.6.1 No Beam

Refer to the relevant chapters of this manual and check the following:

- The tune page real time display is activated by pressing the appropriate button on the tool bar of the tune page.
- Normal tuning parameters are set and, where appropriate, readback values are acceptable.
- All necessary cables have been correctly attached to the source and probe.
- Solvent is reaching the probe tip and the solvent flow rate is as required.

For solvent flow rates below 100 $\mu\text{l}/\text{min}$ it may be necessary temporarily to turn off the nebulizing gas and remove the probe from the source to allow the solvent to be seen at the probe tip.

- The flows of desolvation gas and nebuliser gas are on and are set to the correct flow rates.
- The source has been assembled correctly and is clean.
- The source isolation valve is open.

If, after performing the above checks, the beam is still absent:

1. Acquire TOF data with a mass range extending down to m 20.
2. Check that there is an interference ‘peak’ at approximately m 55 due to the pusher pulse being switched off.

If this interference peak is not present, either the pusher is not pulsing or the output from the detector is not reaching the TDC (time to digital converter).

The most likely cause of an absent pusher interference pulse is a faulty attenuator.

If the pusher interference peak is not present no data will be acquired.

10.6.2 Unsteady Beam

Refer to the relevant chapters of this manual and check that:

- Capillary (electrospray) and Sample Cone are tuned correctly.
- The capillary is not protruding too far from the end of the probe.
- The probe is not too far into the source.
- The flow of solvent from the HPLC pump is correct and steady.

To do this, remove the probe, degas the solvent, increase the flow rate for several minutes to purge any trapped air then reset and re-measure the flow rate.

- Solvents have been adequately degassed.
- The nitrogen flow of desolvation gas and nebuliser gas is steady. The nitrogen supply pressure should be 7 bar (100 psi) $\pm 10\%$.
- Desolvation Temp is not set too high for the liquid flow rate used.

High temperatures can vaporize solvent within the electrospray probe.

Should the preceding checks fail to reveal the cause of the problem, proceed to the following section.

10.6.3 High Back Pressure

For electrospray, a higher than normal back pressure readout on the HPLC pump, together with a slowing of the actual solvent flow at the probe tip, can imply that there is a blockage in the capillary transfer line or injection loop due to particulate matter from the sample.

To clear the blockage:

Remove the probe from the source and increase the solvent flow to 50 $\mu\text{l}/\text{min}$ to remove the blockage.

Often, injections of neat formic acid help to redissolve any solute which has precipitated out of solution.

If the blockage cannot be cleared in this fashion:

Remove the finger-tight nut and tubing from the back of the probe.

If the back pressure remains high:

Replace the tubing with new tube (or first try removing both ends of the tube).

If the back pressure falls:

Replace the stainless steel sample tube inside the probe (or try reversing the tube to blow out any blockage).

Reconnect the tubing to the probe.

The solvent flow can now be readjusted and the probe replaced into the source.

To check the flow rate from the solvent delivery system, fill a syringe barrel or a graduated glass capillary with the liquid emerging from the probe tip, and time a known volume, say 10 μl .

Once the rate has been measured and set, a note should be made of the back pressure readout on the pump, as fluctuation of this reading can indicate problems with the solvent flow.

For APcI a higher than normal back pressure readout on the HPLC pump can imply that, after a long period of use, the filter pad requires replacement.

10.6.4 Loss of Sensitivity

As the ion source becomes dirty after prolonged use, the performance will degrade.

Unstable or reduced ion currents are indicators that the source needs cleaning. The usual remedy is to clean the source as described earlier in this chapter.

An increase in the analyzer pressure above 4×10^{-6} mbar can also cause loss of sensitivity, although the pressure at which this occurs will be sample dependent.

10.6.5 Incorrect Isotope Distributions

Incorrect isotope distributions can be caused by:

- The TDC Stop (mV) threshold being set too high.

Refer to the tune page settings section of Routine Procedures for information regarding the setting of this parameter.

- A faulty attenuator.

Attenuators can fail so that they are open circuit (no beam or pusher interference 'peak' present), or they can fail such that they stop attenuating. The latter failure mode gives rise to incorrect isotope distributions.

When the attenuator fails in this way the TDC **Stop (mV)** threshold can be increased to a significantly higher value than that used previously without reducing the beam intensity.

In normal operation setting the TDC threshold above 200 or 250mV will start to reduce the beam intensity. If the attenuator has failed the TDC threshold can be increased to 500mV or higher before the beam intensity is reduced.

10.6.6 High Noise Levels

High noise levels can either be chemical or electronic in nature.

Chemical Noise

Chemical noise usually originates from contaminated samples, solvents or source gases.

Chemical noise can be distinguished from electronic noise simply by stopping source ionization. If no liquid or gases are entering the source and all the source voltages are set to zero then the remaining noise will be electronic in nature.

Electronic Noise

Electronic noise can be caused by setting the TDC **Stop (mV)** threshold too low. Refer to the Tuning Chapter for information regarding the setting of this parameter.

The microchannel plate detector can be damaged by failure to properly condition the detector following venting of the system to atmosphere. If the detector is producing microdischarges, excessive noise will be apparent on the baseline of mass spectra in the absence of any ion beam. Reducing the detector voltage will reduce the number of discharges and reduce the noise.



Attention: *It is strongly recommended that assistance is sought from Micromass if maintenance to the detector system is thought necessary.*



Attention: *Assistance from Micromass should be sought if, due to symptoms such as excessive noise, spikes, loss of detector gain or abnormal peak shapes, maintenance to any of the components within the TOF analyzer housing is thought to be necessary.*

10.6.7 Poor Analyzer Vacuum

Before suspecting a pump fault or vacuum leak (see Vacuum System earlier in this chapter) it is worth checking the inverted magnetron (Penning) gauge. If this gauge has become dirty it will indicate a poor vacuum, or even fail to register at all.

For information on cleaning the gauge, refer to the Edwards literature supplied with the instrument.



Caution: *The instrument must be vented and electrically isolated at the supply outlet before removing the instrument's covers to gain access to the active inverted magnetron gauge.*

Note: *If the instrument has been vented to atmosphere (instead of dry nitrogen) it may take one to two days before reaching the vacuum levels obtained prior to venting.*

10.7 Cleaning Materials

It is important when cleaning internal components to maintain the quality of the surface finish. Deep scratches or pits can cause loss of performance. Where no specific cleaning procedure is given, fine abrasives should be used to remove dirt from metal components. Recommended abrasives are:

- 600 and 1200 grade emery paper.
- Lapping paper (produced by 3M).

After cleaning with abrasives it is necessary to wash all metal components in suitable solvents to remove all traces of grease and oil. The recommended procedure is to sonicate the components in a clean beaker of solvent and subsequently to blot them dry with lint-free tissue. Recommended solvents are:

- Isopropyl Alcohol (IPA)
- Methanol
- Acetone

Following re-assembly, components should be blown with oil-free nitrogen to remove dust particles.



Attention: Many of the procedures described in this chapter involve the removal of possibly toxic contaminating deposits using flammable or caustic agents. Personnel performing these operations should be aware of the inherent risks, and should take the necessary precautions.

10.8 Preventive Maintenance Check List

- Avoid venting the instrument when the rotary pump is gas ballasting.
- Do not gas ballast the rotary pump for more than 2 hours under any circumstances.

For full details of the following procedures, consult the relevant sections of this chapter and / or refer to the manufacturer's literature.

10.8.1 Daily

- Gas ballast the rotary pump for 30 minutes at the end of a day's megafLOW or APcI operation.

It is normal for the rotary pump noise level to increase during gas ballasting.

10.8.2 Weekly

- Gas ballast for at least 30 minutes by rotating the gas ballast knob anticlockwise by 5 to 6 turns.

When gas ballast is complete, check the rotary pump oil level and color.

Oil that has become noticeably red in color should be replaced.

- Check the water chiller level and temperature (if fitted).

10.8.3 Monthly

Check all cooling fans and filters.

10.8.4 Four-Monthly

- Change the mist element in the oil mist filter.
- Change the oil in the rotary pump.

Gas ballast lightly for 30 to 60 minutes both before and after changing oil.

Appendix A

Reference Information

A

The reference files listed in this chapter have all ion intensities set to 100%. Actual ion intensities are not, of course, all 100%, but the calibration software does not take account of the ion intensities and this is a convenient way to store the reference files in the required format.

Most samples can be purchased from the Sigma chemical company. To order, contact Sigma via the internet, or by toll-free (or collect) telephone or fax:

Internet:

<http://www.sigma.sial.com>

This site contains a list of worldwide Sigma offices, many with local toll-free numbers.

Toll-free telephone:

USA & Canada 800-325-3010

Outside USA & Canada ++1 314-771-5750 (call collect)

Toll-free fax:

USA & Canada 800-325-5052

Outside USA & Canada ++44 314-771-5750 call collect and ask for the fax machine)

Direct fax:

Outside USA & Canada ++1 314-771-5757 (this is a toll call)

A.1 Positive Ion

Ref. File Name	Chemical Name [Sigma Code #]	Molecular Mass	/	Uses
UBQ	Bovine Ubiquitin [U6253]	8564.85	650-1500	General
HBA	Human α globin [H753]	15126.36	700-1500	Hb analysis
SOD	Superoxide dismutase [S2515]	15591.35	900-1500	Hb (internal cal.)
HBB	Human β globin [H7379]	15867.22	800-1500	Hb analysis
MYO	Horse heart myoglobin [M1882]	16951.48	700-1600	General
PEGH1000	Polyethylene glycol + ammonium acetate mixture PEG 200+400+600+1000		80-1000	ES+ and APcI+ calibration
PEGH2000	Polyethylene glycol + ammonium acetate mixture PEG 200+400+600+1000+1450		80-2000	ES+ calibration
NAICS	Sodium Iodide / Caesium Iodide mixture		20-4000	General, ES+ calibration
NAIRB	Sodium iodide / Rubidium Iodide mixture		20-4000	ES+ calibration

A.1.1 Horse Heart Myoglobin

Reference File: MYO.REF - Molecular Weight: 16951.48

Charge State	Calculated m/z Value	Charge State	Calculated m/z Value	Charge State	Calculated m/z Value
28 ⁺	606.419	21 ⁺	808.222	13 ⁺	1304.969
	616.177	20 ⁺	848.583	12 ⁺	1413.633
27 ⁺	628.841	19 ⁺	893.192	11 ⁺	1542.053
26 ⁺	652.989	18 ⁺	942.758	10 ⁺	1696.158
25 ⁺	679.068	17 ⁺	998.155	9 ⁺	1884.508
24 ⁺	707.320	16 ⁺	1060.477	8 ⁺	2119.945
23 ⁺	738.030	15 ⁺	1131.108	7 ⁺	2422.651
22 ⁺	771.531	14 ⁺	1211.829		

A.1.2 Polyethylene Glycol

PEG + NH₄⁺

Reference Files: PEGH1000.REF, **PEGH2000.REF**

Calculated m/z Value					
63.04	459.28	855.52	1251.75	1647.99	
107.07	503.31	899.54	1295.78	1692.01	
151.10	547.33	943.57	1339.80	1736.04	
195.12	591.36	987.60	1383.83	1780.07	
239.15	635.39	1031.62	1427.86	1824.09	
283.18	679.41	1075.65	1471.88	1868.12	
327.20	723.44	1119.67	1515.91	1912.15	
371.23	767.46	1163.70	1559.94	1956.17	
415.25	811.49	1207.73	1603.96	2000.20	

A.1.3 Sodium Iodide and Caesium Iodide Mixture

Reference File: NAICS.REF

Calculated m/z Value					
22.9898		772.4610		1671.8264	2571.1918
132.9054		922.3552		1821.7206	2721.0861
172.8840		1072.2494		1971.6149	2870.9803
322.7782		1222.1437		2121.5091	3020.8745
472.6725		1372.0379		2271.4033	3170.7688
622.5667		1521.9321		2421.2976	3320.6630

A.1.4 Sodium Iodide and Rubidium Iodide

Reference File: NAIRB.REF

Calculated m/z Value					
22.9898		772.4610		1671.8264	2571.1918
84.9118		922.3552		1821.7206	2721.0861
172.8840		1072.2494		1971.6149	2870.9803
322.7782		1222.1437		2121.5091	3020.8745
472.6725		1372.0379		2271.4033	3170.7688
622.5667		1521.9321		2421.2976	3320.6630

A.2 Negative Ion

Ref. File Name	Chemical Name [Sigma Code #]	Molecular Mass	/	Uses
MYONEG	Horse heart myoglobin [M1882]	16951.48	700-2400	General

SUGNEG	Sugar mixture of: maltose [M5885] raffinose [R0250] maltotetraose [M8253] corn syrup [M3639]		100-1500	Low mass range
NAINEG	Sodium Iodide / Caesium Iodide (or Rubidium Iodide) mixture		200-3900	ES- calibration

A.2.1 Horse Heart Myoglobin

Reference File: MYONEG.REF

Calculated m/z Value				
891.175		1209.812		1882.490
940.741		1302.952		2117.927
996.138		1411.615		2420.632
1058.460		1540.036		
1129.091		1694.140		

A.2.2 Mixture of Sugars mixture

Reference File: SUGNEG.REF

Calculated m/z Value				
179.06		665.21		1151.37
341.11		827.27		1313.42
503.16		989.32		1475.48

A.2.3 Sodium Iodide and Caesium Iodide (or Rubidium Iodide Mixture

Reference File: NAINEG.REF

Calculated m/z Value					
126.9045		1026.2699		1925.6353	
				2825.0008	
					3724.3662

276.7987		1176.1641		2075.5296		2974.8950		3874.2604
426.6929		1326.0584		2225.4238		3124.7892		
576.5872		1475.9526		2375.3180		3274.6835		
726.4814		1625.8469		2525.2123		3424.5777		
876.3757		1775.7411		2675.1065		3574.4719		

A.3 Preparation of Calibration Solutions

A.3.1 PEG + Ammonium Acetate for Positive Ion Electrospray and APcI

Prepare a solution of polyethylene glycols at the following concentrations:

PEG 200	25 ng/ml
PEG 400	50 ng/ml
PEG 600	75 ng/ml
PEG 1000	250 ng/ml

Use 50% acetonitrile and 50% water containing 2 mmol ammonium acetate.

Use reference file PEGH1000.REF.

A.3.2 PEG + Ammonium Acetate for Positive Ion Electrospray (Extended Mass Range)

Prepare a solution of polyethylene glycols at the following concentrations:

PEG 200	25 ng/μl
PEG 400	50 ng/ml
PEG 600	75 ng/ml
PEG 1000	250 ng/ml
PEG 1450	250 ng/ml

Use 50% acetonitrile and 50% water containing 2 mmol ammonium acetate.

Use reference file PEGH2000.REF.

A.3.3 Sodium Iodide Solution for Positive Ion Electrospray

Method 1

Prepare a solution of Sodium Iodide at a concentration of 2 µg/µl in 50:50 propan-2-ol (IPA):water with no additional acid or buffer.

Add Caesium iodide to a concentration of 0.05 µg/µl.

The purpose of the caesium iodide is to obtain a peak at m/z 133 (Cs^+) to fill the gap in the calibration file between m/z 23 (Na^+) and the first cluster at m/z 173, which would lead to poor mass calibration in this mass range.

Do not add more CsI than suggested as this may result in a more complex spectrum due to the formation of NaCsI clusters.

Use reference file NAICS.REF.

Method 2

Prepare a solution of Sodium Iodide at a concentration of 2 µg/µl in 50:50 propan-2-ol (IPA):water with no additional acid or buffer.

Add Rubidium iodide to a concentration of 0.05 µg/µl.

The purpose of the rubidium iodide is to obtain a peak at m/z 85 ($^{85}\text{Rb}^+$) with an intensity of about 10% of the base peak at m/z 173. Rubidium iodide has the advantage that no rubidium clusters are formed which may complicate the spectrum. Note that rubidium has two isotopes (^{85}Rb and ^{87}Rb) in the ratio 2.59:1, giving peaks at m/z 85 and 87.

Use reference file NAIRB.REF.

A.3.4 Sodium Iodide Solution for Negative Ion Electrospray

Either of the above solutions is suitable for calibration in negative ion mode. In both cases the first negative reference peak appears at m/z 127 (I^-) and the remaining peaks are due to NaI clusters.

Use reference file NAINEG.REF.

A

Appendix B

Instrument Specifications

To enable the instrument to give its best performance the recommended environmental conditions and power and water supplies are outlined below. Preparation of the laboratory in advance will assist us in efficiently installing the instrument.

B

B.1 Room Location

Dimensions

The Micromass Q-ToF *micro* MS is **180 mm** wide by **635 mm** deep by **1163 mm** long and weighs 200 kg. It is mounted on 6 supporting feet. An external Edwards EM28 requires an additional **650 mm** by **200 mm** of floor space.

In the event that the unit needs to be lifted, the side panels should be removed and suitable lifting bands passed beneath the bottom of the frame. These should be attached to a suitable device (hoist etc.) to raise the analyser unit in a safe and controlled way. Only trained personnel with the correct equipment should carry this out.

A separate table **1200 mm** by **730 mm** is supplied for the computer terminal.

Doorways through which the instrument is to be transported should be a minimum of **820 mm** wide.

In the laboratory a minimum clearance of **500 mm** should be allowed all round the bench for service access, apart from at the rear of the instrument where **900 mm** is required.

Note: *The instrument should not be placed close to heavy machinery (compressors, generators etc.) which give excessive floor vibration.*

B.2 Environment

It is recommended that the instrument be sited in an air conditioned laboratory, in a draught free position and away from excessive amounts of dust.

B.2.1 General

Altitude: upto 2000m

Pollution degree 1 in accordance with IEC 664

Rotary Pump

15⁰ to 40⁰C

Instrument

The maximum ambient laboratory temperature should not exceed 30⁰C, optimum temperature lies in the range **19⁰ - 22⁰C**. Short term (1.5 hour) variations should be no more than **2⁰C**.

The relative humidity should not exceed 70%.

Heat dissipated into the laboratory from the instrument is about 1.2 kW.

The instrument conforms to IEC 1010 - 1, Pollution Degree 1, and Installation Category II.

B.2.2 Magnetic Fields

OA-Tof instruments are relatively tolerant to stray magnetic fields. We would advise an upper limit of 10 Gauss for both AC and DC components of magnetic field measured at the mass spectrometer.

B.2.3 Radio Emissions

The instrument should not be placed within a RF field greater than 0.2 V/metre. This approximates to a 1W hand held unity gain transmitter at a distance of 10 m.

Possible sources of RF emission include RF linked alarm systems or LANs, portable telephones and hand held transmitters.

B.3 Water Supplies

The heat dissipated into the cooling water is about 400 Watts. The water flow required to dissipate the heat generated by the turbo pumps is 35 L/hour for an inlet temperature of about 200 C or 23 L/hour at an inlet temperature of about 150 C assuming an outlet temperature of about 300 C.

The water may be supplied by a recirculating chiller with the following characteristics:-

- Heat dissipation into system: 400 W
- Temperature stability: $\pm 2^{\circ}\text{C}$
- Minimum reservoir volume: 5 L
- Minimum supply pressure: 10 psi (outlet at atmospheric pressure).
- Maximum supply pressure: 60 psi
- Minimum flow rate at 15°C : 0.4 L/min

The above assume that the outlet water temperature will be no more than 30°C . Inlet water temperatures below 15°C are not recommended since excessive condensation may form on exposed pipework.

Alternatively, when there is a cooled water supply available it may be used either directly through the instrument or indirectly via a water-to-water heat exchanger. In this latter case the chilled water supply to the heat exchanger must be at least 10°C below the required inlet temperature for the instrument.

One inlet and one outlet are required for the instrument. Reinforced 10mm (3/8 inch) flexible hose is preferred.

To prevent blockage of the water pipes suitable in-line filters will be required to remove particulate matter from town water supplies when these are used.

Operation above 2000 metres altitude may adversely affect the cooling of the system.

B.4 Power Requirements

The instrument requires a single phase **50 - 60 Hz, 230 V** nominal power supply rated at **13 A** (UK) or **15 A** (Europe).

An additional single phase **50 – 60 Hz, 230 V** or **115 V** nominal power supply rated at 5A (UK) is required to run the embedded PC.

In the USA and Canada a single-phase **50 – 60 Hz** supply at **230 V** phase to neutral fused and rated at **15 A** is required. Alternatively two phases of a **50 – 60 Hz 208 V** phase to phase, **3** phase supply, rated and fused at **15 A** may be used. It is mandatory that no other apparatus is connected to this supply.

Circuit breakers are an acceptable alternative to fuses. The supply should be terminated in the laboratory no more than **2 m** from the instrument with either a wall mounted isolator or socket and plug to be fitted to the instrument.

Other supply voltages can be accommodated using a transformer to change the primary supply voltage to **230 V**. Advance notice is required and Micromass should be contacted.

On single-phase supplies the power supply should ensure that the line and neutral wires cannot be transposed.

On pump start-up currents of up to **30 A** may be drawn for several seconds. Time delay fuses and breakers are recommended to prevent nuisance tripping.

A safety earth (ground) correctly rated must be provided in all cases.

Data system components, chromatographs, syringe pumps etc. should be connected directly to laboratory power outlets (no ancillary outlets are provided on the instrument).

A residual current device (RCD) is recommended for additional protection. In the case of instruments fitted with a transformer the RCD should be fitted in the supply side of the transformer.

Supply brownout should not fall to less than half main voltage for greater than **20 msec** duration.

B.5 Gases and Regulators

Nitrogen

The instrument requires oil free dry nitrogen regulated at 7 bar (100 psi) minimum outlet pressure to provide nebulising and drying gas to the instrument.

During API operation typical usage of nitrogen is about 400 L/hour, but under high flow rate conditions (Megaflow/APcI) this may increase to 650 L/hour. This equates, approximately, to the consumption of a large cylinder of compressed nitrogen each day and it may be preferred to use a liquid nitrogen dewar which may last several weeks.

Collision Gases

Typically **Argon** is used as the collision gas for CID experiments. This should be 99.9% pure, regulated at no more than **50 psi**. Connection is via **1/8 inch** OD stainless steel or copper tubing (NOT SUPPLIED).

B.6 Exhaust Outlets

Rotary Pump Outlet

The rotary pump exhaust outlet must be vented to the atmosphere external to the laboratory clear from any air intakes for air conditioning systems. A **12 mm (1/2 inch)** hose connection is required. If the length of exhaust exceeds **4 m** then the internal diameter of the pipe should be increased to **48 mm (2 inch)** for the excess distance.

Nitrogen Outlet

Severe damage to the instrument will result if the electrospray/APCI exhaust is connected to the rotary pump exhaust line. This will occur when the nitrogen supply is off and rotary pump oil vapour will migrate via the source exhaust to the ion source and then through the sampling orifice into the quadrupole and gas cell assembly.

A separate exhaust for the ion source gas (nitrogen) must be provided to the atmosphere external to the laboratory clear from any air intakes for air conditioning systems. A (6mm OD) hose connection is required. If the length of exhaust exceeds **3 m** then the internal diameter of the pipe should be increased to **12 mm (1/2 inch)** for the excess distance.

B

B.7 Performance Specifications

Time of Flight Mass Resolution, Positive Ion

5000 (FWHM) on $(M+H)^+$ ion from Leucine Enkephalin.

Time of Flight Mass Resolution, Negative Ion

5000 (FWHM) on $(M-H)^-$ ion from Raffinose.

Full Scan MS Sensitivity, Positive Ion

The signal height obtained from a sample consumption of 200 fmol of horse heart Myoglobin (16952 Da) will be greater than 166 ion counts on the most intense peak in the charge state envelope.

A solution of 200 fmol/ μ L horse heart Myoglobin (in 50/50 acetonitrile/water + 0.2% formic acid) will be introduced at a flow rate of 5 μ L/min.

Full Scan MS Sensitivity, Negative Ion

The signal height obtained from the sample consumption of 1 ng of raffinose will be greater than 200 counts on the $(M-H)^-$ peak at m/z 503. This will correspond to a signal to noise ratio of greater than 200:1 (after a 1x3 smooth).

The instrument will be tuned at 5000 resolution (as demonstrated in specification 2a) and a solution of 5 ng/L in 50/50 acetonitrile/water (no additives) will be introduced at 10 μ L/min.

Full Scan MS/MS Sensitivity, Positive Ion

The signal height obtained from a consumption of 20 fmol of $[Glu^1]$ -Fibrinopeptide B (1569 Da) will be greater than 6 counts on the most intense y_1 sequence ion from the MS/MS spectrum of the doubly charged precursor ion.

This will correspond to a signal to noise ratio of greater than 30:1 (*after a 3x9 smooth*) on the most intense y_1 sequence ion.

A solution of 100 fmol/ μ L will be introduced at a flow rate of 5 μ L/min.

Mass Measurement Accuracy (with internal reference)

The RMS error between the measured and the accepted masses of peaks which have sufficient intensity and are free from interference from other masses, over the range from 150-900 Daltons, will be less than 5 ppm.

One suitable peak of known mass will be used as an internal reference.

The instrument will be tuned at 5000 resolution as demonstrated in specification 1a.

B.7.1 Electrospray Option - Nanoflow

Full Scan MS/MS Sensitivity, Positive Ion

The signal to noise from a consumption of 2 fmol of $[Glu^1]$ -Fibrinopeptide B (1569 Da) will be greater than 30:1 (*after a 3x9 smooth*) on the most intense y_1 sequence ion from the MS/MS spectrum of the doubly charged precursor ion.

A solution of 500 fmol/L concentration in MeOH/H₂O + 0.2% formic acid solution and with glass micropipettes with 1 or 2 μ m tips will be used.

The integration period per spectrum will be about 5 sec and data will be summed over a period appropriate for the required consumption of sample.

B.7.2 Electrospray Option – Transform Software

(Software for the determination of molecular weight from a spectrum containing a series of multiply charged ions on a m/z scale by a transform of the data to a true mass scale.)

Mass Measurement Accuracy (no internal reference)

The mean measured mass of transformed data shall be 15867.2 ± 0.5 Da and the standard deviation (σ) of the mean < 0.5 Da.

The transform data will be created from five repeat analyses of the globin from normal human haemoglobin. Mass calibration to be performed using the multiply charged α globin peaks from a separate analysis. The raw data should be transformed over the range 15,000 - 16,000 Da and smoothed appropriately.

The instrument will be tuned at 5000 resolution. It is recommended that a solution containing 10 pmol/ μ L of each globin in 50/50 acetonitrile/water + 0.2% formic acid is used.

B

Index

A

- Accept Calibration 105
- Accurate Mass 93
- Acquisition Status Window 114
- Add Task 34
- Ammonium Acetate 192
- Analog Channels 120
- Analyzer
 - Maintenance 176
- APCI 67
 - Obtaining an Ion Beam 69
- APCI Probe
 - Maintenance 172
- Argon 198
- Atmospheric Pressure Chemical Ionization
 - see APCI 67
- Automatic Startup / Shutdown 28

B

- Borosilicate Capillary 55

C

- Caesium Iodide 190, 191
- Calibration 93, 102, 113
- Calibration Parameters 104
- Calibration Solutions 192
- Calibration Window 102
- Camera System
 - Nanoflow 60
- Capillary 82
- Capillary Voltage
 - ElectroSpray 48
- Centroid 123

- Charge State 129, 146
- Charge State Recognition 154
- Checking the ESI Probe 44
- Checking the Probe
 - APCI 68
- Checking the Resolution 90
- Chromatogram
 - Real-Time Update 115
- Chromatographic Interfacing 52
- Cleaning Materials 184
- Cleaning Sample Cone 160
- Collision Energy
 - Parent 152
 - Survey 135
- Collision Energy Profile 153
- Collision Gases 198
- Combine 91
- Components 10
 - Electrical 10
 - Mechanical 10
- Cone Gas
 - ElectroSpray 48
- Configure Contact Closures 31
- Configure Shutdown on Error 30
- Continuum 123
- Conventions, documentation 3
- Copy Function 118

D

- Daily Maintenance 185
- Data Acquisition 111
- Data File Name 112
- Data System 18
- Deadtime Correction 106
- Delete Task 34

- Desolvation Gas
 - APCI 72
 - ElectroSpray 47
- Detected Precursor Inclusion 146
- Dialogs
 - Accept Calibration 105
 - Analog Data 120
 - Calibration Parameters 105
 - Display Calibration Graphs 103
 - Gas Threshold 33
 - MCP Conditioning 23
 - Scan Report 114
 - Shutdown on Error Configuration 31
 - Solvent Delay 119
 - TDC Settings 88
 - TOFMS Scan 122
 - Tune Page Acquisition 91, 112
- Dimensions 195
- Documentation
 - conventions 3
 - related 1

E

- Electrical Components 10
- Electronics 17
- Electronics Maintenance 156
- Electrospray 39
- ElectroSpray Operation 43
- Electrospray Probe
 - Maintenance 170
- Electrospray Tip 45
- Environment 195
- Error
 - Shutdown 30
- ESI interface 39
- ESI Source 43
- Event Out 17
- Exclude Masses

- Parent 151
 - Survey 137
- Exhaust Outlets 199
- Experiment Setup Editor 115, 116
 - Toolbar 116
- Experimental Record 78
- Extraction Cone Voltage 49

F

- Fault Finding 155, 179
- Four-Monthly Maintenance 185
- Front Panel Connections 14
- Front Panel Controls 15
- Function List 116

G

- Gas Ballasting 158
- Gas Thresholds 32
- Gases and Regulators 198
- Glass Capillary Option 61

H

- Horse Heart Myoglobin 189, 191

I

- Include Masses
 - Parent 149
 - Survey 133
- Infusion Pump 40
- Injection Valve 40
- Insert Task 34
- Installing Nanoflow 57
- Instrument Specifications 195
- Instrument Warm-up 24

Intensity Threshold 104
Internal Layout 9
Ion Optical System 76
Ion Optics 8
Ionisation Techniques
 APci 6
Ionization Techniques
 Electrospray 6
 Nanoflow Electrospray 7

L

LC column 39
LC-MS Sensitivity 53
Lock Mass 89
LockMass Correction 106
Log Parameters 37
Lteff 89, 94

M

Maintenance
 Electrospray Probe 170
 Source 159
Maintenance 155
 Analyzer 176
 APcI Probe 172
Manual Pusher 87
MCP 87
MCP Conditioning 23
Mechanical Components 10
MegafLOW 43
Microscope
 Nanoflow 60
Mobile phase 39
Modify Function 117
Modify Task 34
Monthly Maintenance 185
MS to MSMS

 Parent 145
MS to MSMS Switch Criteria 129, 145
MSMS Template
 Parent 144
 Survey 132
MSMS to MS 131
 Parent 148
 Survey 131
Multi Channel Analysis (MCA). 123
Multiple Sample Acquisition 114
Myoglobin 189, 191

N

NanoFlow Electrospray 55
Nanoflow Interface 55
Nano-LC 55
Nano-LC Assembly 64
Nano-LC Option 63
Nanovial Tip Position 61
Nebuliser Gas
 ElectroSpray 47
Negative Ion 190
Neutral Loss
 Parent 143
New Function 117
Nitrogen 198
Nitrogen Gas
 In 16
Nitrogen Outlet 199
Nitrogen Supply 28
Nominal Mass Accuracy 94
NP Multiplier 88

O

Obtaining an Ion Beam 89
 APcI 69
 ElectroSpray 45

Oil Mist Filter 159
 Operation
 ElectroSpray 43
 Nanoflow 65
 Origin 113
 Outlets 199

P

Parent
 Collision Energy 152
 Exclude Masses 151
 Include Masses 149
 MS to MSMS 145
 MSMS Template 144
 MSMS to MS 148
 Neutral Loss 143
 Product Ions 142
 Parent Ion Discovery 139
 Parent Survey Window 140
 Peak Display
 Tune Page 79
 Performance Specifications 199
 Pirani Gauge 157
 Polyethylene Glycol 189
 Polyethylene glycol 189
 Positive Ion 188
 Post-column Splitting 41
 Power Failure 27
 Power Requirements 197
 Precursor Selection 146
 Preparation
 APCI Operation 26
 Electrospray Operation 25
 For Startup 21
 Preparation of Calibration Solutions 192
 Preventive Maintenance 184
 Probe Position
 ElectroSpray 46

Product Ions
 Parent 142
 Purge Gas
 ElectroSpray 48

Q

Quadrupole Tuning Page 83

R

Rear Panel Connections 16
 Reciprocating Pump 40
 Reference compound 187
 Related documentation 1
 Remove Function 118
 Resolution 90
 Restore Function List 120
 Room Location 195
 Rotary Pump 13, 17
 Rotary pump 13
 rotary pump 13
 Rotary Pump Oil 159
 Rotary Pump Outlet 199
 Routine Procedures 21
 Rubidium Iodide 190

S

Sample Cone 82
 Sample Cone Voltage 49
 Shutdown Editor 28
 Toolbar 29
 Shutdown Log 36
 Shutdown Page 29
 Sodium Iodide 190, 193
 Sodium iodide 191
 Solvent Delay 119
 Source

- Maintenance 159
- Source Temp 82
- Source temperature 48
- Source Tuning Page 81
- Spectrum
 - Real-Time Update 115
- Standard Tune Parameters 90
- Start Up
 - After Complete Shutdown 21
 - After Overnight Shutdown 24
- Starting an Acquisition
 - Multiple Samples 114
 - Tune Page 111
- Status Display 15
- Sugar mixture 191
- Sugars 191
- Supply Inlet 17
- Survey
 - Collision Energy 135
 - Exclude Masses 137
 - Include Masses 133
 - MS to MSMS 128
 - MSMS Template 132
 - MSMS to MS 131
- Switch Criteria
 - Parent 142
- Syringe Pump 40
- System Manager 115

T

- TDC Settings 88
- The Auto Control Tasks Page 33
- The Source 159
- Time Of Flight Page 85
- TOFMS Scan 122
- Transient Pressure Trip 27
- Tune Page 77
- Tune Page Acquisition 111

- Tuning 75
- Tuning and Optimization
 - ElectroSpray 46
- Turbomolecular pumps 12

V

- Vacuum Leaks 156
- Vacuum Protection 13
- Vacuum System 12
- Veff 89

W

- Water 16
- Water Supplies 196
- Weekly Maintenance 185

